

MAXI BATHROOM RANGE

QUICK REFERENCE GUIDE


BTS118020 Aspire MAXI Shower Chair BTS119030 Aspire MAXI Shower Stool BTT145020 Aspire MAXI Over Toilet Aid


Carton Contents:

- 1 x Aspire MAXI Shower Chair or Stool or Overtoilet Aid
- 1x Seat
- 1x Backrest
- 1 x User Guide

1. Installing Seat (Over Toilet Aid Only):

- a. Line up the 2×10^{-2} x holes in the seat base with the inner two holes on the mounting plates and fasten with the bolts provided (shown with red arrows).
- b. Secure in place with the two **self-taping screws** provided (holes shown in red circles).

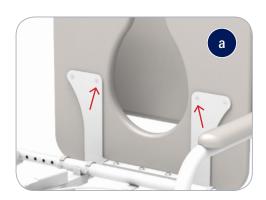

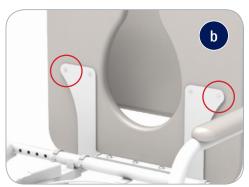





WARNING! Ensure that the seat is secured before use.

2. Installing Backrest:

- a. Insert the backrest tube into the holes at the rear of the seat. Once the desired depth is reached, insert the safety pins.
- b. Ensure the fastening hand wheels are tightened and the backrest is rigid before sitting.

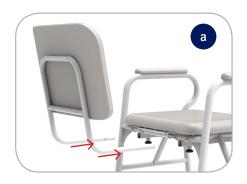



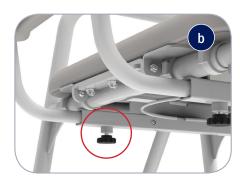



WARNING! Ensure that the fastening hand wheels are always tight when the chair is in use to prevent injury.







MAXI BATHROOM RANGE

3. Adjusting Height:

- Remove C-clip (x4) and adjust legs to desired height, ensuring that all legs are located in the same adjustment points to prevent an unstable frame.
- b. Reattach C-clip (x4) to fasten and secure in place.

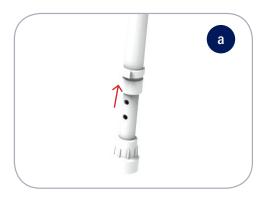

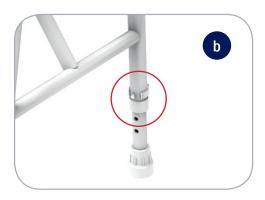





WARNING! Ensure that the C-clip is secure at all times.

4. Adjusting Width:

- Remove C-clip (x4) and adjust frame to desired width.
- Ensure both sides of the frame are uniformly adjusted to maintain stability.
- Reattach C-clip (x4) to fasten and secure in place.

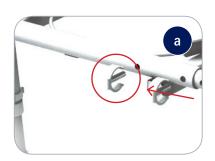



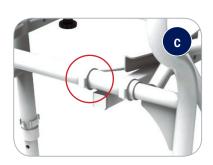



WARNING! Ensure that the C-clips are secure at all times.





MAXI BATHROOM RANGE

5. Service & Maintenance

To prolong the working life of MAXI Whitegoods and for trouble-free operation, it is recommended that routine inspections be performed by a person who is suitably qualified and well acquainted with the design, use and care of whitegoods.

If inspection reveals that the safety of the product is jeopardised in any way, the whitegoods frame should be IMMEDIATELY WITHDRAWN FROM SERVICE.

Maintenance Guide & Checklist

Clean surfaces with a pH neutral detergent - chlorine and phenol based detergents should not be used as these may damage commode components

The Aspire MAXI Bathroom range should be wiped clean and kept dry as much as possible after use.

AS A MINIMUM A FULL SERVICE SHOULD BE COMPLETED EVERY 12 MONTHS TO MAINTAIN WARRANTY.

6. Warranty

The Aspire MAXI Bathroom range carries a comprehensive warranty as detailed below. Please keep proof of purchase for proof of warranty commencement.

Consumer Warranty Against Defects:

Aidacare is bound by Australian Consumer Law in respect to Warranty Against Defects. A copy of the Aidacare Warranty Against Defects statement is available on the Aidacare homepage at www.aidacare.com.au. All warranty matters should be directed to Aidacare Customer Service on 1300 133 120.

Aidacare does not warrant against excessive or incorrect use, modification or any situation that could not be deemed as fair wear and tear.

This is a back to base warranty and does not cover freight costs pertaining to the return of any items under warranty. Aidacare acknowledges and conforms to all ACCC guidelines around consumer rights.

For further information contact Aidacare on 1300 133 120 or email product@aidacare.com.au.

Aidacare will not warrant the safety and or correct functioning of products where any original components have been changed or modified by non-Aidacare approved and trained service & maintenance staff or external providers. Furthermore, safety is not guaranteed where components have been replaced with non-original Aidacare approved parts. If any faults are detected upon receipt of this product please phone Aidacare. Any faults that are detected during normal use should be reported to Aidacare immediately to determine if warranty conditions apply and if so, the necessary repair or replacement work to be completed.

Spare parts list is available upon request.

For servicing, preventative maintenance and any other questions regarding this or any other Aidacare product please contact:

Aidacare Pty Ltd

4 Noonan Rd, Ingleburn NSW 2565

T - 1300 133 120

W - www.aidacare.com.au

E - product@aidacare.com.au



ALL WARRANTIES MADE ARE BACK TO BASE WARRANTIES MEANING THAT ANY PARTS OR FRAMES REQUIRING REPAIR NEED TO BE RETURNED TO THE NEAREST AIDACARE SERVICE CENTRE AT THE CUSTOMER'S EXPENSE - PLEASE CONTACT AIDACARE FOR A CURRENT LIST OF SERVICE CENTRE LOCATIONS.

