

Bed Alertamat

This mat is antibacterial and latex free

For safety, replace this mat when mat life expires

CAUTION: This product is designed as a monitor alerting the staff when a patient has risen from their bed. This product is not a restraint and therefore will not prevent a patient from getting out of their bed. This product is to be used in conjunction with the total fall prevention program.

DO NOT ASSUME: This product must be tested each time it is used with a patient to make sure it functions properly.

NOTE: For patient safety, fill in the start date below when mat is first put into use. It is recommended to discontinue use of mat and replace one year from start date.

Start Date End Date


Use start date to track mat use.



Tested and approved to
IEC 60601-1-2
Medical Devices

Manufactured to
ISO 9001

Alertamat and Monitor Instructions

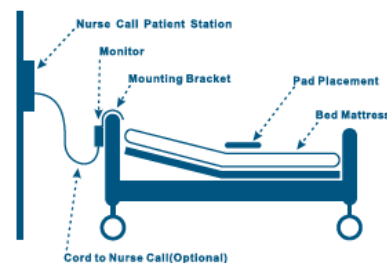

SET-UP INSTRUCTIONS

1. Position the sensor mat under the user. To be alerted when the user sits up in bed, place the **bed Alertamat** under the user's shoulder. To be alerted when the user leaves the bed, place the **bed Alertamat** under the user's hips.
2. Install batteries in the alarm monitor at back cover using screwdriver. Put **bed Alertamat** connecting pin into opening of the alarm monitor. **CAUTION: ALWAYS** remove **bed Alertamat** pin from alarm monitor by pushing down the locking clip of the connecting pin, before pulling it out. DO NOT attempt to remove connecting pin any other way or you may damage pin and void warranty.
3. Secure Alarm Monitor to bottom or back of bed, ensuring it is out of sight or out of reach of the user. Secure excess wire cord to avoid tripping and damage.
4. Carefully assist the user into the bed.
5. Push "Test/Reset" on alarm monitor to test the set-up. The alarm will sound and shut off to indicate it is functioning properly. If the user gets off the **bed Alertamat**, the alarm will sound to alert the staff. When the user gets back on the mat, the alarm will re-set automatically.
6. When staff remove the user from the mat, the alarm will sound. To silence alarm, press "Test/Reset" button.



MAINTENANCE

- When using a **bed Alertamat** in conjunction with foam or gel cushions, test the system to make sure it functions properly.
- To minimize spreading infection, it is recommended to use a new **bed Alertamat** for each patient.
- To clean **bed Alertamat**, use disinfectant wipes or anti-bacterial cleaners.
- **DO NOT** bend, fold or roll **bed Alertamat**.
- **DO NOT** immerse **bed Alertamat** pad in any liquid cleaning solution.
- **DO NOT** expose the **bed Alertamat** to high temperature sanitizers.

TROUBLESHOOTING

1. Make sure that the **bed Alertamat** is plugged into the fall monitor.
2. Check the power supply. Ensure correct batteries are used in fall monitor and batteries are in good condition and replace if necessary. Make sure the power adaptor is plugged into power outlet and switched on, if present.
3. Check monitor settings behind the cover of the battery. To hear monitor ringing tone, select Mute 'OFF'.
4. If the steps above do not resolve your problem, call the supplier to troubleshoot. Make sure both the **bed Alertamat** and the fall monitor are with you at the time of the call.


www.alertamedical.com

WARNING

FAILURE TO COMPLY WITH WARNING MAY RESULT IN INJURY OR DEATH.

The effectiveness of the **bed Alertamat** depends on proper equipment installation and function. Always test the system prior to use. This mat is designed to help alert staff of situations that lead to or result in falls, and is not a substitute for proper caregiver supervision. Please ensure that patients at high risk of falls are sufficiently monitored and help limit environmental hazards that can add to fall risk. Please read set-up instructions and Legal Disclaimer before each use.

LEGAL DISCLAIMER

The **Alerta bed Alertamat System** is designed to be installed by the end-user. As such, it is the entire responsibility of the buyer to ensure that the system is properly installed and tested, further, the system is not designed to replace good caregiving practices including, but not limited to:

- a) Direct patient supervision
- b) Adequate training for staff personnel for resident fall prevention and elopement
- c) Testing of the system before each use

WHERE PERMISSIBLE, IT IS THE SELLER'S INTENT TO LIMIT ANY LIABILITY FOR SPECIAL, CONSEQUENTIAL OR INCIDENTAL ECONOMIC DAMAGES TO REFUND OF THE PURCHASE PRICE OR REPLACEMENT OF THE GOODS. SELLER DISCLAIMS ANY LIABILITY FOR COMPENSATORY OR OTHER DAMAGES ARISING FROM THE USE OF OUR GOODS.

LIMITED WARRANTY

All products supplied are guaranteed to be free from defects in material & workmanship for a period of 1 year from time of delivery. Our liability shall be limited to replacement of defective material ONLY. Neither the seller nor manufacturer shall be liable for any injury, loss or damage, direct or consequential, arising out of the use or the inability to use these products. Before using, the buyer shall determine the suitability of the product for the intended use. Buyer shall assume all risk and liability whatsoever in connection herewith.

Except for this limited warranty, NO EXPRESS OR IMPLIED WARRANTIES ARE GIVEN AND NO AFFIRMATION OF SELLER BY WORDS OR ACTION WILL CONSTITUTE A WARRANTY. THE ENTIRE RISK AS TO THE QUALITY AND PERFORMANCE OF THE GOODS ARE WITH THE BUYER. THE BUYER, AND NOT THE MANUFACTURER, DISTRIBUTOR OR RETAILER, ASSUMES THE ENTIRE COST OF ALL NECESSARY INSTALLATION, SERVICING AND REPAIRS.

Distributed by

aidacare™
HEALTHCARE EQUIPMENT

1300 133 120 | aidacare.com.au