

Wireless Bed Alertamat

This mat is antibacterial and latex free

For safety, replace this mat when mat life expires

CAUTION: This product is designed as a monitor alerting the staff when a patient has risen from their bed. This product is not a restraint and therefore will not prevent a patient from getting out of their bed. This product is to be used in conjunction with the total fall prevention program.

DO NOT ASSUME: This product must be tested each time it is used with a patient to make sure it functions properly.

NOTE: For patient safety, fill in the start date below when mat is first put into use. It is recommended to discontinue use of mat and replace one year from start date.

Start Date	End Date


Use start date to track mat use.



Tested and approved
to EN 60601-1-2
Medical Devices

Manufactured to
ISO 9001

Alertamat and Monitor Instructions

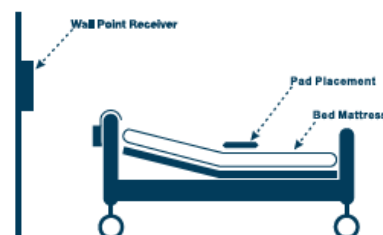

SET-UP INSTRUCTIONS

1. Position the sensor mat under the user. To be alerted when the user sits up in bed, place the **bed Alertamat** under the user's shoulder. To be alerted when the user leaves the bed, place the **bed Alertamat** under the user's hips.
2. Plug bed mat cable connection in to **Alerta Wireless Transmitter**, install battery into **Alerta Wireless Transmitter** and turn transmitter to 'ON' position. The indicator blinks once to show that the power is normal. Then close the flap over the transmitter and stud shut. Please also refer to separate **Alerta Wireless Transmitter** user manual.
3. **Alerta Wireless Transmitter** now must be 'PAIRED' with **Alerta Wall Point Receiver**. Please refer to separate **Alerta Wall Point Receiver** user manual.
4. Once successful pairing, please test by applying pressure to the **bed Alertamat**, when pressure is lifted alarm will sound. Once tested working OK, carefully assist the user into the bed.
5. If the user gets off the **bed Alertamat**, the alarm will sound to alert the staff. When the user gets back on the mat, the alarm will re-set automatically.
6. When leaving the mat, the alarm will sound. To silence alarm, press "Reset" button.



MAINTENANCE

- * When using a **bed Alertamat** in conjunction with foam or gel cushions, test the system to make sure it functions properly.
- * To minimize spreading infection, it is recommended to use a new **bed Alertamat** for each patient.
- * To clean **bed Alertamat**, use disinfectant wipes or anti-bacterial cleaners.
- * **DO NOT** bend, fold or roll **bed Alertamat**.
- * **DO NOT** immerse **bed Alertamat** pad in any liquid cleaning solution.
- * **DO NOT** expose the **bed Alertamat** to high temperature sanitizers.

TROUBLESHOOTING

1. Make sure that the **bed Alertamat** is plugged into the **Alerta Wireless Transmitter**.
2. Check the power supply. Ensure correct batteries are used in **Alerta Wall Point Receiver** and **Alerta Wireless Transmitter** and batteries are in good condition and replace if necessary. Make sure the power adaptor is plugged into power outlet and switched on, if present.
3. Check **Alerta Wall Point Receiver** settings: behind the cover of the battery. To hear monitor ringing tone, select volume low or high.
4. Reprogram the **Alerta Wireless Transmitter** with the **Alerta Wall Point Receiver** unit.
5. If the steps above do not resolve your problem, call the supplier to troubleshoot. Make sure both the **bed Alertamat** and the **Alerta Wall Point Receiver** are with you at the time of the call.

Product Code: W-BAM



www.alertamedical.com

Alerta Medical, 4 Symington Place, Riverside Business Park, Irvine KA11 5DE

WARNING

FAILURE TO COMPLY WITH WARNING MAY RESULT IN INJURY OR DEATH.

The effectiveness of the **bed Alertamat** depends on proper equipment installation and function. Always test the system prior to use. This mat is designed to help alert staff of situations that lead to or result in falls, and is not a substitute for proper caregiver supervision. Please ensure that patients at high risk of falls are sufficiently monitored and help limit environmental hazards that can add to fall risk. Please read set-up instructions and Legal Disclaimer before each use.

LEGAL DISCLAIMER

The **Alerta bed Alertamat System** is designed to be installed by the end-user. As such, it is the entire responsibility of the buyer to ensure that the system is properly installed and tested, further, the system is not designed to replace good caregiving practices including, but not limited to:

- a) Direct patient supervision
- b) Adequate training for staff personnel for resident fall prevention and elopement
- c) Testing of the system before each use

WHERE PERMISSIBLE, IT IS THE SELLER'S INTENT TO LIMIT ANY LIABILITY FOR SPECIAL, CONSEQUENTIAL OR INCIDENTAL ECONOMIC DAMAGES TO REFUND OF THE PURCHASE PRICE OR REPLACEMENT OF THE GOODS. SELLER DISCLAIMS ANY LIABILITY FOR COMPENSATORY OR OTHER DAMAGES ARISING FROM THE USE OF OUR GOODS.

LIMITED WARRANTY

All products supplied are guaranteed to be free from defects in material & workmanship for a period of 1 year from time of delivery. Our liability shall be limited to replacement of defective material ONLY. Neither the seller nor manufacturer shall be liable for any injury, loss or damage, direct or consequential, arising out of the use or the inability to use these products. Before using, the buyer shall determine the suitability of the product for the intended use. Buyer shall assume all risk and liability whatsoever in connection herewith.

Except for this limited warranty, NO EXPRESS OR IMPLIED WARRANTIES ARE GIVEN AND NO AFFIRMATION OF SELLER BY WORDS OR ACTION WILL CONSTITUTE A WARRANTY. THE ENTIRE RISK AS TO THE QUALITY AND PERFORMANCE OF THE GOODS ARE WITH THE BUYER, THE BUYER, AND NOT THE MANUFACTURER, DISTRIBUTOR OR RETAILER. ASSUMES THE ENTIRE COST OF ALL NECESSARY INSTALLATION, SERVICING AND REPAIRS.

Distributed by

aidacareTM
HEALTHCARE EQUIPMENT

1300 133 120 | aidacare.com.au