


JAY - J3

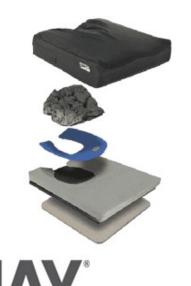

- The Jay-J3 series of cushions are highly customisable to suit your client's needs
- Customisations cover 6 spectrums: (1) Cushion height; (2) Cushion width & depth; (3) PLA size (Pelvic Loading Area); (4) PLA insert choice (4 options); (5) Cover Options (3 options); (6) Positioning components
- We strongly recommend that you consult with your Aidacare C.A.S.E team member for any inquiries relating to these cushions
- The J3 features two models (1) the Gel and (2) the Air (with ROHO DRY FLOATATION[®] insert (see below for more detail)
- The J3 range features outstanding pressure distribution, comfort, heat/moisture dissipation, and optimal stability
- With 230 kg weight capacity the J3 is also a bariatric option
- Pelvic loading area (PLA): Optimally shaped and sized to evenly distribute pressure under the pelvis. Available in 3 sizes, the PLA can be selected and sized to fit the user's pelvis regardless of the overall cushion width
- The PLA insert is available in four options: (1) Factory Filled JAY Flow[™] fluid;
 (2) Field Variable JAY Flow fluid; (3) Air Single ROHO DRY FLOATATION[®] Technology; and, (4) Air Dual ROHO DRY FLOATATION[®] Technology. These four inserts provide two important benefits: hydrostatic loading (even distribution of load) and adjustability (in-field adjustment or custom factory adjustment) to ensure the best possible immersion, pressure distribution, and stability
- The J3 cushion features three cover options: (1) Stretch; (2) Microclimatic; or
 (3) Incontinent-resistant. All 3 covers are designed to minimise shear risks;
 provide comfort; and wick away heat and moisture
- The J3 is available with a wide range of positioning options to help the user maintain appropriate postural stability. Optional lateral and medial thigh supports as well as lateral pelvic supports are affixed securely to the cushion

Contracted DVA product

\oslash

ndis Item Code:

05_181006190_0105_1_2 Pressure Care Cushion – Composite:05_429_0103_1_2 Click here for more information

Compliance

No applicable compliance standard



base

Warranty Information

- 2 Years Cushion
- 6 Months Cover

Code	Display Name	Safe Working Load	Pressure Care Rating	Material
PTC536100	JAY J3 - As Per Script	225 kg	Very High Risk	Foam