

Softform Premier Active Mattress






- The Softform Premier Active features an alternating air insert beneath the castellated foam insert pad
- The mattress retains its properties as a high risk static mattress, but should clinical judgement determine that the patient requires stepping up to a dynamic surface, then a discreet pump can be fitted to the air insert, transforming the mattress into an alternating surface delivering additional levels of pressure relief
- There is no requirement to move the patient, thereby minimising any discomfort
- The pump will assess individual patients' weights, and adjust levels of inflation accordingly
- The air cells alternate on a ten minute cycle beneath the foam insert, delivering a very low impact at the point of change, thereby delivering high levels of comfort in addition to improved pressure reduction
- Two-way stretch PU waterproof cover
- High frequency welded seams
- Full length flap conceals the zip
- Toughened PU coated base
- Double coated PU sides

Warranty Information

- Cover 4 Years
- Foam 8 Years
- Pump 2 Years

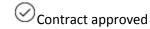

Support @ Home

Function: Manage tissue integrity

Category: Prescribed

Click **here** for more information

NDIS

Item Code:

05 043306002 0103 1 2

Power Bed – Composite:

05_423_0103_1_2

Click here for more information



Code	Display Name	Total Length	Total Width	Total Height	Safe Working Load
PTM636395	Softform Premier Active Mattress	1970 mm	880 mm	152 mm	248 kg
PTM636396	Alternating Air Pump	-	-	_	-