



Fixed Hoists up to 450kg


Aidacare is an experienced supplier and installer of ceiling hoist systems. We understand the risks involved in the installation and use of ceiling hoists and can help during early-stage specification and design to reduce risk and improve user-acceptance.

Our systems are world-class and are accompanied by accredited installation guidelines which allow us to work with a range of installers. We are an ISO 9001 accredited company and consistently provide high-quality outcomes for Australia's largest Aged Care and Hospital providers.

A wide range of motor & track options are available to suit client specifications. Some common examples are provided below.

Ceiling Hoist Considerations

- A. Structural Attachment Point
- B. Threaded Rod
- C. Ceiling Bracket
- D. Ceiling Track
- E. Traverse Boom (Track)
- F. Charging Station
- G. X-Y Gantry
- H. Ceiling Lift Unit
- I. Carry Bar
- J. Lift Strap
- K. Sling
- L. Gate Assembly
- M. Linear Track
- N. Modified Door Header

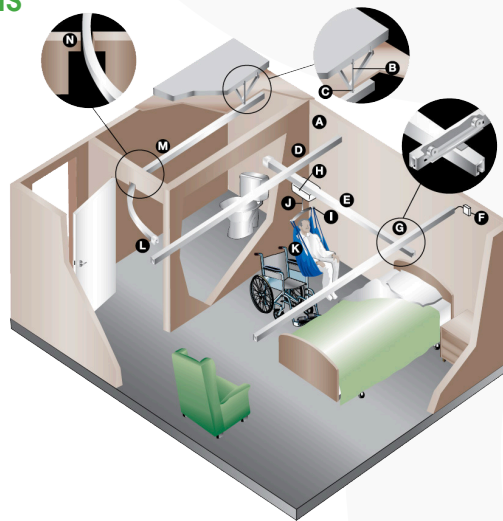


Portable Hoists up to 270kg


Consulting Solutions

Aidacare has extensive experience with installation and user-acceptance of ceiling systems. Simple things can be identified during planning to reduce risk in the installation and improve user and patient outcomes. We can assist in developing the best solution. Some examples of installations are provided below.


Linear Track with Turn Table

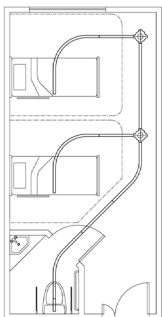


Straight Monorail

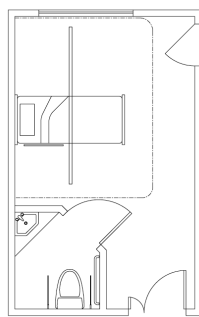

J Track




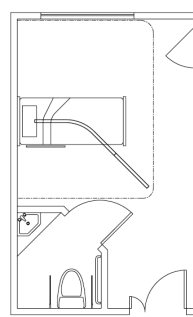


Turntable

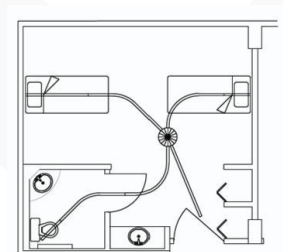


Installation Solutions


No two installations are alike. Aidacare has extensive experience dealing with many different installation substrates and scenarios. Aidacare utilises manufacturer-certified installation guidelines. Some of the more complex examples are shown below. This leads to a consistent high-quality outcome for our clients.

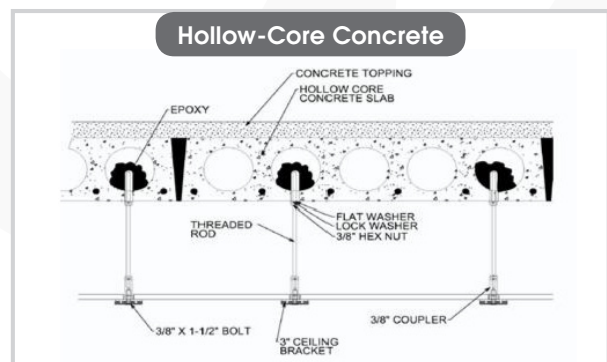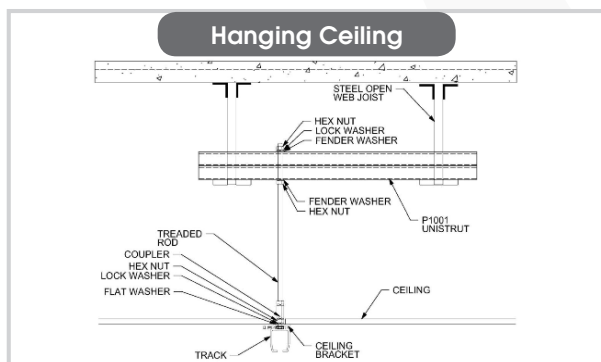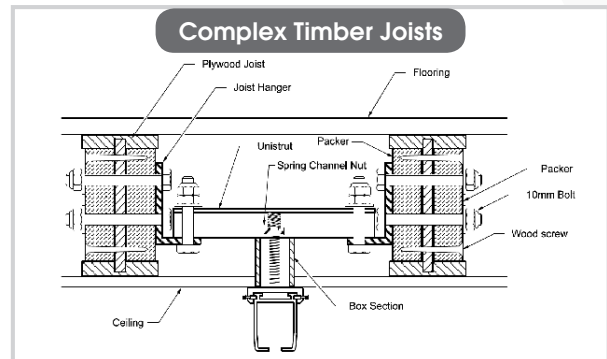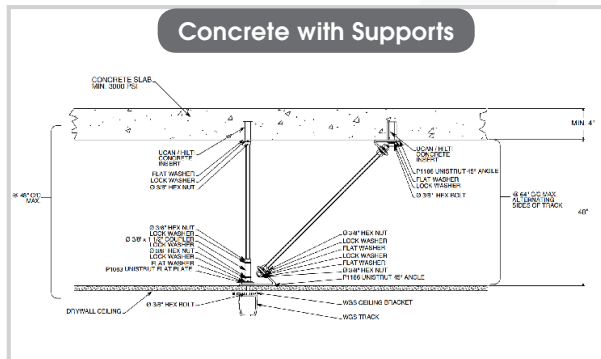


Track Options

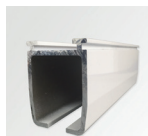

Exposed Track

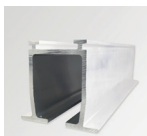

Concealed Track

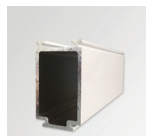


Rebar



45° Curve Track



90° Curve Track

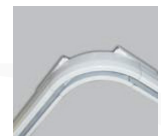

Quick Curve

Maintenance Safety

Aidacare fully-supports its installations with ongoing maintenance programs. Aidacare also assists our clients to develop their own internal maintenance regime.

Training & User Acceptance

Aidacare Equipment Specialists offer user training and ongoing support. We also provide WHS signage with all equipment.

Quality & Compliance

- All systems are tested and installed in compliance with Australian Standard AS/NZ ISO 10535-2011 Hoists for the transfer of disabled persons.
- All systems are maintained in accordance with manufacturer's guidelines and in accordance with AS/NZ 3551-2012 Management Programs for Medical Equipment.
- Aidacare is an ISO9001 compliant company.
- Aidacare uses an accredited network of trained installers, engineers, and load testers.

