

200S STANDING LIFTER

The 200S Standing Lifter empowers patient mobility and rehabilitation through its specialised bullhorn lifting boom with six attachment points. This steel-constructed mobility solution features interchangeable clip sling compatibility, enabling optimal standing transfers.



Interchangeable Clip Sling Bar

One-for-one swap of existing sling bar converts stander to accept clip slings only (retro-fit).

Product code: LSS390655



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset, and Cables





Intelligent Diagnostics

State of the art Linak Jumbo Care control system. Tracks usage and service requirements for safety and compliance.



Build Quality

Strong powder-coated aluminium frame with industry leading Linak electronics, Tente castors and warranty provision.



Bullhorn Lifting Boom

Offset patient handles provide stability and 6 sling attachment points give the carer flexibility to perform the optimal lift in any scenario.




Adjustable Padded Legrests

Provide comfort and support during transfers and can be adjusted to suit the patient.




Ergonomic Handles

Provide multiple grip points to suit carers of different heights and strengths.

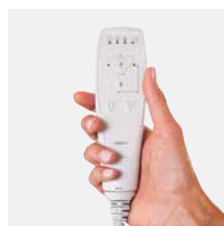

Smart Handset

The latest slimline technology from Linak fits neatly in the hand with an intuitive control panel. Battery and service indicator lights assist with service scheduling.

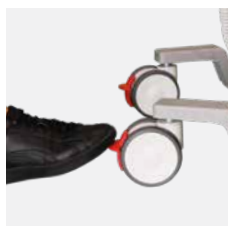

TENTE

2 x 75mm front castors for maximum under bed clearance and accessibility.

2 x 75mm rear-locking castors for safety and manoeuvrability.



Safety


Emergency stop button, emergency lowering, and anti-entrapment technology, keeps both patient and caregiver safety at the forefront.



Powered Leg Function

Opens and closes the legs for lifting stability or accessibility.

ACCESSORY

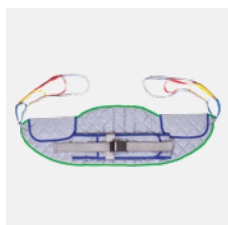

Sling Compatibility

Compatible with a range of Aspire slings.



See our Sling Fitting and Safe Use Guide for full details

Code	Frame Material	Lifting Range	Under Bed Clearance	Turning Diameter	Base Length	External Base Width	Internal Base Width	Overall Height	Lifter Weight	SWL
LSS390650	Steel	875 - 1777	110	1361	1277	656 - 1203	561 - 1103	1832	40kg	200kg

All measurements are in millilitres unless otherwise specified.