

A205 FLOOR LIFTER

The A205 Floor Lifter combines portability with professional performance through its lightweight aluminium construction and tool-free disassembly design. Featuring three adjustable height settings, intelligent Linak diagnostics, and powered leg functionality, this versatile lifter seamlessly adapts to home, community, and facility environments.



Manual and Powered Pivot Frame Upgrades

Integrated hoist system with rotating frame and clip sling help to reduce carer effort.

Manual: LSS390670

Powered: LSS390665



Transport or Storage

Is quick and convenient with tool-free disassembly.



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset, and Cables





Adjustable Height

Three adjustable height settings provides the caregiver with the flexibility to adjust the height range to suit the specific lift and patient scenario.



Build Quality

Strong powder-coated aluminium frame with industry leading Linak electronics, Tente castors and warranty provision.



Intelligent Diagnostics

State of the art Linak Jumbo Care control system. Tracks usage and service requirements for safety and compliance.




Quick-release Yoke Attachment

One-touch release and attach mechanism, designed for the busy caregiver. Padded yoke for added patient safety.




Powered Leg Function

Opens and closes the legs for lifting stability or accessibility.

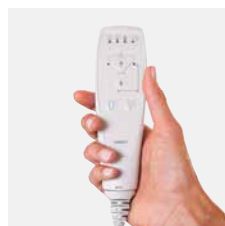

Smart Handset

The latest slimline technology from Linak fits neatly in the hand with an intuitive control panel. Battery and service indicator lights assist with service scheduling.

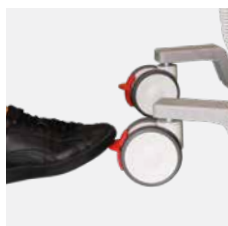

TENTE Castors

2 x 75mm front castors for maximum under bed clearance and accessibility.
2 x 100mm rear-locking castors for safety and manoeuvrability.



Safety

Emergency stop button, emergency lowering, and anti-entrapment technology, with both patient and caregiver safety at heart.

ACCESSORIES



In-line Weigh Scale

Clear LCD display with a 320kg SWL.


Product Code: LSS390680



Integrated Weigh Scale

Simple LCD screen paired with easy-to-clean padding, and a 200kg SWL.

Product Code: LSS390675

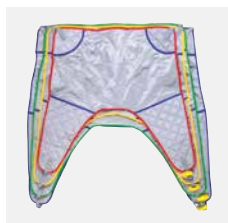

Sling Compatibility

Compatible with a range of Aspire slings.



See our Sling Fitting and Safe Use Guide for full details

Code	Frame Material	Lifting Range*	Under Bed Clearance	Turning Diameter	Base Length	External Base Width	Internal Base Width	Overall Height	Lifter Weight	SWL
LSS390610	Aluminium	415 - 515 1700 - 1800	110	1332	1238	651 - 1173	547 - 1073	2040	38.1kg	205kg

All measurements are in millimetres unless otherwise specified | *Three mast height settings are possible hence the range of minimum and maximum heights.