

A320 BARIATRIC FLOOR LIFTER

The A320 Bariatric Floor Lifter excels with its high 320kg safe working load through its reinforced aluminium construction and broad-base design for enhanced patient clearance. Equipped with folding support arms, adjustable height settings, and 4-point yoke system, it delivers uncompromising stability for complex bariatric transfers.



Quick-release 4-point Yoke

Padded 4-point yoke for added patient safety. One-touch release and attach mechanism assists the busy caregiver.



Manual and Powered Pivot Frame Upgrades

Integrated hoist system with rotating frame and clip sling help to reduce carer effort.

Manual: LSS390670
Powered: LSS390665



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset, and Cables





Intelligent Diagnostics

State-of-the-art Linak Jumbo Care control system. Tracks usage and service requirements for safety and compliance.



Build Quality

Strong powder-coated aluminium frame with industry leading Linak electronics, Tente castors and warranty provision.



Adjustable Height

Three adjustable height settings provides the caregiver with the flexibility to adjust the height range to suit the specific lift and patient scenario.



Smart Handset

The latest slimline technology from Linak fits neatly in the hand with an intuitive control panel. Battery and service indicator lights assist with service scheduling.



Ergonomic Handles and Support Arms

Provide multiple grip points to suit carers of different heights and strengths. Folding support arms assist with stability or mobility but can also be folded out of the way.



Safety

Emergency stop button, emergency lowering, and anti-entrapment technology, with both patient and caregiver safety at heart.



TENTE Castors

2 x 100mm front castors for reduced rolling resistance under high load.
2 x 125mm rear-locking castors for safety and manoeuvrability.



Powered Leg Function

Opens and closes the legs for lifting stability or accessibility.

ACCESSORIES



In-line Weigh Scale

Clear LCD display with a 320kg SWL.


Product Code: LSS390680



Integrated Weigh Scale

Simple LCD screen paired with easy-to-clean padding, and a 200kg SWL.

Product Code: LSS390675

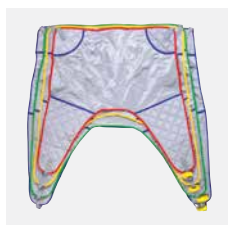

Sling Compatibility

Compatible with a range of Aspire slings.



See our Sling Fitting and Safe Use Guide for full details

Code	Frame Material	Lifting Range*	Under Bed Clearance	Turning Diameter	Base Length	External Base Width	Internal Base Width	Overall Height	Lifter Weight	SWL
LSS390630	Aluminium	540 - 640 1840 - 1940	115	1595	1500	765 - 1190	640 - 1065	2155	63kg	320kg

All measurements are in millimetres unless otherwise specified | *Three mast height settings are possible hence the range of minimum and maximum heights.