

AC4 BED

The AC4 Bed offers a perfect blend of modern design, clinical-grade functionality, and everyday practicality. It's a dependable choice for aged care and residential environments. With a sleek gunmetal frame and optional teak surrounds or upholstered panels, the AC4 delivers a warm, homely feel that fits seamlessly into any interior setting.

Built with trusted TiMotion electronics, the bed provides smooth, reliable adjustments and is designed for long-term performance with minimal maintenance required.




10 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics


1 YEAR WARRANTY
Handset



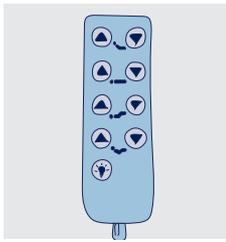

Designed for Carers and Residents

TiMotion electronics are up to 20% faster than other brands, saving carers time. Simple and intuitive handset is easy to operate.

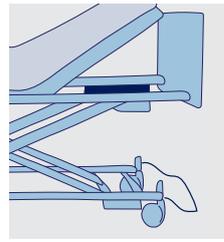


Considered Cable Management

Ensures cables are kept safely out of the way, minimising damage and downtime.




Maintains Infection Control

Wipe-clean mattress platform allows easy cleaning for infection control compliance.

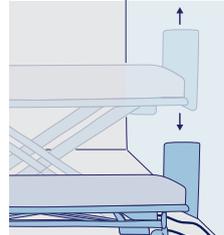


Ensures Wall Protection

Wall locator bar included. Bed rises vertically, thereby minimising costly wall collisions while maximising functional bedroom space at the foot-end.

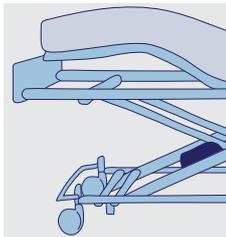

Easy Access to the Control Box


The control box is located conveniently within reach underneath the foot section, allowing easy access for maintenance staff.

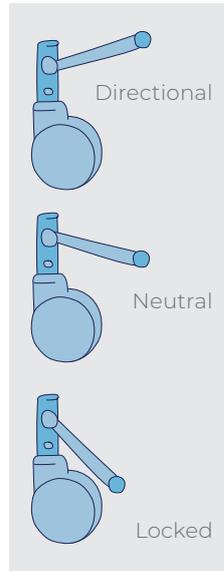

Enhanced Manual Handling

125mm Twin Castors
Offer lower rolling resistance for easier bed repositioning.

Central Locking - Total and Directional
Activated centrally via an easily accessible full-width foot pedal, push down for total lock, or lift up for directional lock to facilitate safe and easy tracking through corridors and doorways.

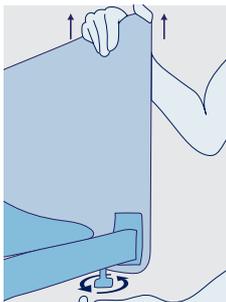

Personalised Sleeping Solution

Quick-attaching and easily removable headboards and footboards.

200mm mattress platform extension caters for taller residents.

Optional Melamine ends or upholstered Allure surround available for an alternative homely style.

ACCESSORIES



Upholstered Allure Surround

Single BEA829150
King Single BEA829160




Universal 3/4 Folding Bed Side Rails

BEA829100



Clamp-on Bed Rail

BEA829105

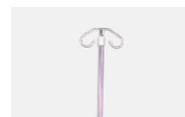

IV Pole

BEA029710



Self-help Pole

BEA829140



Battery Back-up

GSP405130

Code	Mattress Dimensions	Overall Length	Overall Width	Height Range	Backrest Angle	Knee Raise Angle	Product Weight
Single BEB050400	2000 x 900	2100*	920	300 - 730	70°	40°	112kg
King Single BEB050410	2000 x 1060	2100*	1070	300 - 730	70°	40°	116kg

All measurements are in millilitres unless otherwise specified | *External length of 2100mm + 200mm extension