
# FLX Floorline Bed

Refined Floorline Care





### FLX Floorline Bed

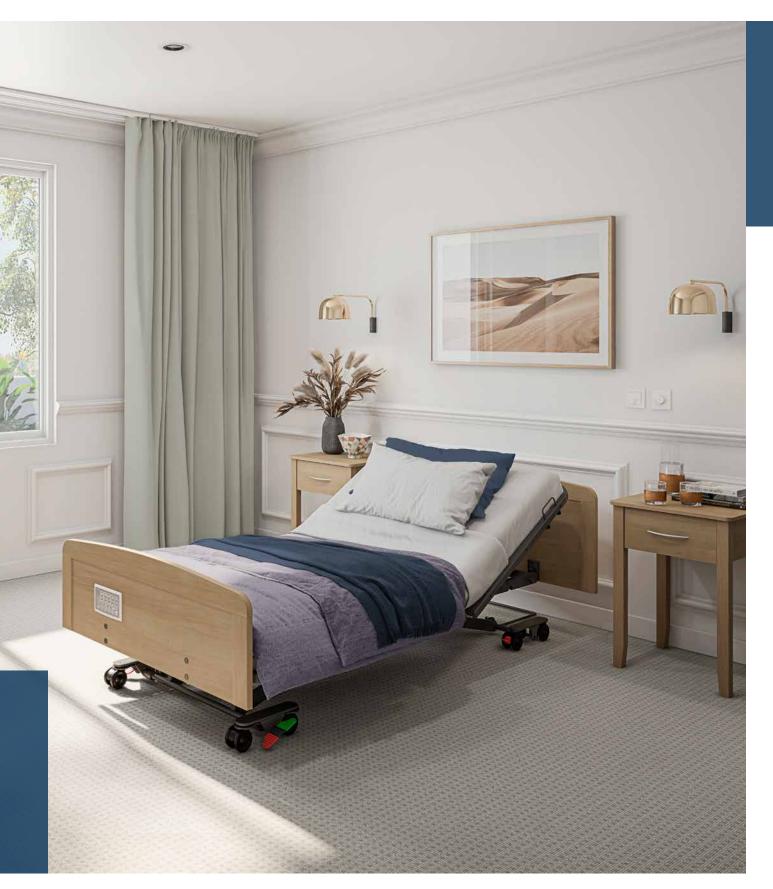

#### Included Features

- Respite Ready 275kg Safe Working Load
- 2 UltraLow 100mm to 750mm Height Range
- 3 Single and King Single Sizes
- 4 Four Section Wipe Clean Mattress Platform
- 5 Adjustable Powered Backrest

- 6 Adjustable Powered Legrest
- 7 Anti-Shear Auto-Regression Backrest
- 8 Trendelenburg / Reverse Trendelenburg
- 9 Backlit Handset
- 10 Lockout Controls on Nurse Control Panel

- Gentle Halo Underbed Lighting
- 12 Wall Buffer Bar
- 13 Integrated Bed Extension from 100-200mm
- 14 Side and Foot-End Mattress Retainers
- One-Touch Foot-End Brakes with On / Off Highlights

- 16 Stylish, Custom and Bespoke Bed Ends and Side Panels
- One Button Chair Positioning
- 18 One Button Entry / Exit Height
- 19 Adjustable Height Settings

### Best Practice Care

#### Safety First

- · Ultra Low Height: Market leading 100mm low height whilst maintaining full ground clearance.
- · Attendant Control Panel: Allows caregivers functional control and the ability to lock out certain bed functions from being operated via the handset.
- · Underbed Clearance: Full clearance under the bed minimises the risk of cord impingement and protects flooring from damage.
- · Custom Entry/Exit Height Function: Set the exact safe transfer height for each resident. Effective across a large range of user heights.
- · SmartDescent: Bed lowering speed is reduced once the bed reaches 300mm maximizing carer potential to remove obstacles as the bed lowers to the floor.
- · Underbed Lighting: Creates a gentle halo around the bed to reduce falls risk and improve carer safety.

#### **Clinical Considerations**

- · Backrest Auto-Regression: Reduces shear forces on the spine and sacral regions. Minimizes abdominal compression when back is raised supporting cardio pulmonary and respiratory function.
- · One Touch Chair Positioning: Improves social interaction, assists respiratory and cardio pulmonary function and creates a comfortable position for sitting, dining and transfers.
- · Integrated Bed Extension: 100-200mm length extension is easily operated and ensures a large range of user heights will fit comfortably without compromise.
- · Trendelenburg (and Reverse): Augments patient positoning, comfort and assists with bed mobility\*.

\*Carers should seek advice from a trained Allied Health professional for correct manual handling techniques.

#### Infinite Positioning









LEG RAISE



SEATED

TRENDELENBURG







## Interior Styling

#### **Contemporary Designs**

Choose from our contemporary standard bed ends and sides; made to order styles; or completely bespoke designs of your choosing to meet any interior design requirements.

The FLX bed offers a suite of aesthetic options from our subtle Como Radius bed ends (standard) and side panels, to custom shapes and finishes.

#### Standard Bed Ends - Como Radius





#### Standard Side Panels





#### **Custom Designs**



Generous outer frame and recessed inner panel with void for a luxurious aesthetic and a softer silhouette. Timber or laminate.



Eyre Elegant frame and recessed inner panel with curved void detail creates an aperture to the surrounding room. Timber or laminate.



Como Curved Sweeping curves provide a gentle aesthetic to compliment a robust outer frame design. Timber or laminate.



Sophia Luxurious upholsterv and subtle lines provide a gentle aesthetic to compliment interior design accents.

#### Bespoke Finishes

Don't see a design that suits your needs? We can customise as required to incorporate your choice of material, colour, finish, size or shape and match any interior design specifications.

FLX Floorline Bed - Refined Floorline Care

### Performance and Durability

- Quick Swap Actuators: Allow for quick and easy inspection and maintenance, minimising downtime.
- TiMotion Electronics: Medical grade electronics from a leading global manufacturer.
- Control Box: Easily accessible beneath the back rest panel for efficient servicing.
- Quick Connect Bed Ends: Simple and quick to set-up upon bed installation. Easily removed for access or maintenance.
- Wall Buffer: Length adjustable to suit room configuration.

- · Cable Conduit: Allows other cables like air mattress power cords to be safely tucked away, preventing trip hazards and cord damage.
- · EasyMove Transport System: Facilitates safe no-lift transportation and movement of beds, as well as safe and easy access underneath the bed platform for routine inspection and maintenance.
- · Width Extension Kit: Ability to increase a single to a king single, ideal for resident changes in a facility.
- · High Performance Twin Castors x 8: Head end total locking (x2), foot end central locking (x4)
- · Respite Ready: Generous 275kg SWL facilitates care of bariatric residents and augments standard bed assets to facilitate respite care.

#### **Technical Specifications**

Patient Surface – Single	2000 x 900mm
Patient Surface – King Single	2000 x 1060mm
Overall Length	2130 – 2330mm
Length Extension	100/200mm
Overall Width – Single	964mm
Overall Width – King Single	1124mm
Patient Surface Height Range	100 – 750mm
Safe Working Load	275kg
Backrest Range	0 – 62°
Knee Bend Range	0 – 40°
Underbed Clearance	Variable. Min 20mm

Mattress Platform	4 Sections
Product Weight – Single	105kg
Product Weight – King Single	115kg
Castors	8 x Twin 75mm
Electrical Parameters	
Power Input	240V AC 50Hz, Max current 4000mA
Rate Load	Interruption 10%, Max 2 minutes / 18 minutes
Leakage Protection	Class II Type B
Waterproof Grade	IP54

Comprehensive user manual, troubleshooting guide, safety checklist, assembly instructions available.

#### Beds

BEB048600 FLX Floorline Bed – Single BEB048610 FLX Floorline Bed – King Single


#### Standard Bed Ends (Sold Separately)

BEB048700	Bed End Set – Como Radius – Sublime Teak – Single
BEB048701	Bed End Set – Como Radius – Milano Walnut – Single
BEB048710	Bed End Set – Como Radius – Sublime Teak – King Single
BEB048711	Bed End Set – Como Radius – Milano Walnut – King Single

#### Standard Side Panels (Sold Separately)

BEB046265	Side Panel Set – Sublime Teak
BEB046275	Side Panel Set – Milano Walnut



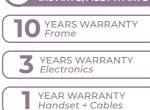

### Accessories



















- 1. **EasyMove Transport System** BEA030055 Is an ideal transport solution designed to minimise manual handling. It allows effortless lifting and tilting of the bed to its transport position, while the large 100mm castors offer smooth and easy manoeuvring for one person. The EasyMove is also a perfect maintenance platform facilitating easy access to mechanical parts under the bed platform.
- 2. Split Folding Grab Rails (Pair) BEB046525 Left/Right pair of moulded rails provides manoeuvring and repositioning assistance to a resident. Can be positioned at the head end or foot end. Folds down to lowest position beside the bed platform.
- 3. **Self Help Pole** BEA829140 Self Help Pole Mounting Bracket BEA829145 Provides a firm support handle for improved bed mobility

- 4. **IV Pole** BEA029710 IV Pole Mounting Bracket BEA829145
- 5. Battery Backup BEB046592 Prolong bed operation by up to 15 lifts during a power outage.
- 6. Width Extension Kit BEB046500 Allows conversion of a Single FLX into a King Single.
- 7. **USB Port Upgrade Kit** BEB046616 Add USB charging functionality to the side of the bed, perfect for charging mobile phones and tablets.
- 8. Handset Holder BEA030030 Mounts discreetly behind the head board for easy access and improved room aesthetics.
- 9. Directional Castor BEB046630 Replaces total locking castor at the head end for improved bed tracking during bed movement.

06 FLX Floorline Bed - Refined Floorline Care 07



ABBR13.17 010823

