


# aidacare | AC3 SINGLE PRODUCT CODE: BEB045900

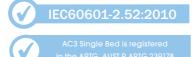

The Aidacare AC3 bed is a three function bed tailored to the Aged Care environment. Designed in Australia to meet International Standards, the premium frame and Linak electronics will ensure long term reliability and safety.

Multiple functions and accessories combine to exceed the comfort, support and safety expectations for a diverse range of users.

### **Bed Functions** BACKREST INCLINE/DECLINE HIGH SAFE WORKING LOAD **BED HEIGHT** INCLINE/DECLINE RANGE 220kg 285 - 685 mm 70° 35°

#### **LEADING QUALITY FRAME & ELECTRONICS**

Australian designed durable steel frame with European Linak electronics and handset allows for market leading warranty provision.



#### **REMOVABLE MATTRESS PLATFORM**

Facilitates easy cleaning and infection control compliance.



#### REMOVABLE LAMINATE BED ENDS

In a contemporary beech finish. Any finish available to order.





#### **INTEGRATED BED EXTENSION**

Provides an additional 195mm length to cater for all users.



#### **CABLE MANAGEMENT**

Ensures cables are kept safely out of the way for residents and care staff. Easy access service compartment for electronics.



#### **CENTRAL LOCKING**

125mm twin castors with full width foot pedal.



A wall locator bar facilitates bed positioning and helps prevent wall damage.



#### **WARRANTY**

**ELECTRONICS 3 YEARS** 









BEA029700

Self Help Pole BEA029705

IV Pole BEA029710

BEA030030




# aidacare | AC3 KING SINGLE

The AC3 King Single bed is a three function bed with the added comfort of a larger sleeping surface.

It combines the benefits of the AC3 Single bed, including the premium frame and Linak electronics with the added safety of a wider mattress platform.

Multiple functions and accessories combine to exceed the comfort, support and safety expectations for a diverse range of users.

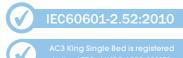



#### **Bed Functions** BED HEIGHT RANGE BACKREST INCLINE/DECLINE KNEE BEND INCLINE/DECLINE











#### **EXTENDABLE MATTRESS RETAINER**



### **STYLE**

285 - 685 mm

The AC3 King Single bed is designed to facilitate a homely environment for residents. Head and Foot boards can be removed to meet interior design requirements and the bed frame is compatible with valances, wooden surrounds, wall mounted head boards and comforter bed covers.

#### **SAFETY**

Mobility challenges, dementia, vertigo, memory loss and vision impairment are common ailments and contributors to falls risk. The wide surface of the AC3 King Single bed reduces the risk of falling out of bed. Residents feel more comfortable. whilst carers can have greater peace of mind for a resident's safety.

#### **COMFORT**

Transitioning from a home to a care environment is life changing and can be daunting. Most residents are used to the space of a double or queen size bed so moving to a single bed can be disconcerting. The AC3 King Single bed provides more room for residents to move and reposition in familiar comfort.

### **INCLUDES:**

**Integrated Bed Extension** 

**Removable Laminate Bed Ends** 

**Central Locking** 

**Removable Mattress Platform** 

**Leading Quality Frame** & Electronics

**Cable Management** 

**Wall Protection System** 

## **Optional Accessories**



BEA029700



Self Help Pole

BEA029705





Handset Holder Bracket BFA030030

#### **WARRANTY**

FRAME

**ELECTRONICS 3 YEARS** 


## aidacare | AC2 SINGLE



The Aidacare AC2 bed offers streamlined functionality without compromise, tailored for the Aged Care environment.

Designed in Australia to meet International Standards, the premium frame and Linak electronics will ensure long term reliability and safety.

Multiple features and accessories combine to meet the comfort, support and safety expectations for a diverse range of users.

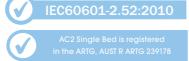

# BED HEIGHT RANGE 285 - 685 mm



BACKREST INCLINE/DECLINE



HIGH SAFE WORKING LOAD



#### **LEADING QUALITY FRAME & ELECTRONICS**

Australian designed durable steel frame with European Linak electronics and handset allows for market leading warranty provision.



#### **REMOVABLE MATTRESS PLATFORM**

Facilitates easy cleaning and infection control compliance.



#### REMOVABLE LAMINATE BED ENDS

In a contemporary beech finish. Any finish available to order.





#### **INTEGRATED BED EXTENSION**

Provides an additional 195mm length to cater for all users.





kept safely out of the way for residents and care staff. Easy access service compartment for electronics.





**TOTAL LOCK** SINGLE CASTORS with dust guards

#### WALL PROTECTION **SYSTEM**

A wall locator bar facilitates bed positioning and helps prevent wall damage.



### **Optional Accessories**



RFA029700





Self Help Pole BFA029705



BFA029710



Handset Holder Bracket BFA030030

#### **WARRANTY**

ELECTRONICS 3 YEARS

Model	AC3 Single	AC3 King Single	AC2 Single	
Image				
Code	BEB045900	BEB045905	BEB045901	
Height Range	285 – 685 mm	285 – 685 mm	285 – 685 mm	
Overall Length	2120 mm	2140 mm	2120 mm	
Overall Length including Wall Protection System	2230 mm	2250 mm	2230 mm	
Overall Width	920 mm	1090 mm	920 mm	
Overall Width (Rails Fitted)	1030 mm	1200 mm	1030 mm	
Backrest Incline	0 – 70°	0 – 70°	0 – 70°	
Knee Bend	0 – 35°	0 – 35°	N/A	
S.W.L	220 kg	220 kg	220 kg	
Recommended Mattress Size	2000 x 900 mm	2030 x 1050 mm 'King Single'	2000 x 900 mm	
Warranty	Frame – 10 years Electronics – 3 years Handset – 1 year	Frame – 10 years Electronics – 3 years Handset – 1 year	Frame – 10 years Electronics – 3 years Handset – 1 year	

Features	AC3 Single	AC3 King	AC2 Single
Central Locking Twin Castors	✓	✓	-
'Total Lock' Single Castors (125mm)	-	-	$\checkmark$
Full Width Central Locking Pedal	✓	✓	-
Integrated Bed Extension	✓	✓	✓
Wall Protection System	✓	✓	✓
Mattress Retainers	✓	✓	✓
Removable Contemporary Bed Ends - Beech	✓	✓	✓
Profiling Bed Base	✓	✓	✓
Large Lifting Range	✓	✓	✓
Linak Electronics And Handset	✓	✓	✓
High S.W.L (220kg)	✓	✓	✓
Removable Mattress Platform	✓	✓	✓







### MATTRESS RANGE

#### **FEATURES**

- Australian made covers using the highest quality European two way stretch fabric.
- Australian made foam
- · Fire retardant
- Unique serial number
- · Removable Cover:
- 1. Easy inspection
- 2. Facilitates better hygiene
- 3. Fully washable

- Optimum infection control:
  - 1. Fully welded seams
  - 2. Water resistant zippers
  - · Reinforced waterfall flaps
  - Mattresses do not require rotation, minimising OH&S handling risks.



#### **LIFECOMFORT**

Dual layer support surface incorporating a full size memory foam overlay and durable high density foam base. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for maximum infection control

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046410	Lifecomfort - AC Series Single	2000mm	900mm	130mm	5 Years Foam 2 Years Cover	120kg
BEM046460	Lifecomfort - King Single	2030mm	1050mm	130mm	5 Years Foam 2 Years Cover	120kg



#### LIFECOMFORT DELUXE

Three layer PUPPS compliant support surface incorporating a full size memory foam overlay and two layers of high quality foam for maximum comfort and support. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for optimum infection control.

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046510	Lifecomfort Deluxe - AC Series Single	2000mm	900mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046455	Lifecomfort Deluxe - King Single	2030mm	1050mm	150mm	5 Years Foam 2 Years Cover	150kg



#### LIFECOMFORT BARIATRIC

Four layer foam, designed for larger users with a maximum user weight of 300kg. Full size memory foam overlay and side walls for added stability and ease of transfers. European fabric cover with fully welded seams, water resistant zippers and reinforced waterfall flaps for optimum infection control.

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046700	Lifecomfort Bariatric - King Single	2030mm	1050mm	200mm	5 Years Foam 2 Years Cover	300kg



#### LIFECOMFORT PRESSURE REDUCING

Top of the range PUPPS2 compliant mattress featuring a quad layer support surface with additional side wall supports. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for optimum infection control.

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046610	Lifecomfort Press Red – AC Series Single	2000mm	900mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046660	Lifecomfort Press Red - King Single	2030mm	1050mm	150mm	5 Years Foam 2 Years Cover	150kg

#### **ASSET MANAGEMENT**

Covers have individual serial numbers heat transferred into base material. This allows accurate age tracking on fixed asset registers.



#### **INFECTION CONTROL**

- 1. Fully welded seams
- 2. Reinforced waterfall flaps
- 3. Water resistant zippers





## ALTERNATING PRESSURE SYSTEMS



#### **FEATURES**

- Compact, powerful and ultra-quiet control units
- Two way stretch fully welded top cover – high temperature washable
- Waterfall flaps over all zippers for protection and infection control
- Medical Grade power cords run inside covers to avoid hazards
- CPR emergency release function
- Static head sections to eliminate dizziness

Air-Lo Alternating Pressure Systems offer the latest in pressure prevention and treatment. The unique patented technology is designed to provide maximum support, release and reperfusion to prevent and heal pressure ulcers in patients up to 200kg.



#### **AIR-LO 3 OVERLAY**

The Air-Lo 3 Alternating Overlay is effective in the prevention of pressure areas and skin breakdowns. Air-Lo 3 supports patient loads of up to 120kg and should be used in conjunction with a foam mattress or underlay. 21 pure TPU transverse air cells alternately allow pressure to be released across the skin surface. A great value choice for safely caring for patients assessed as at risk of developing pressure areas.

Code	Description	Length	Width	Height	Warranty	S.W.L
PTM636629	Air-Lo 3	2000mm	900mm	100mm	1 Year Cover & Softgoods, 2 Years Pump	120kg



#### **AIR-LO 5 OVERLAY**

The Air-Lo 5 Alternating Overlay leads the field in the prevention of pressure areas and skin breakdowns. Air-Lo 5 supports patient loads of up to 170kg and should be used in conjunction with a foam mattress or underlay. 18 pure TPU transverse air cells alternately allow pressure to be released across the skin surface. When only the very best prevention and support is required, choose Air-Lo 5 to safely care for patients assessed as at high risk of developing pressure areas.

Code	Description	Length	Width	Height	Warranty	S.W.L
PTM636630	Air-Lo 5	2000mm	900mm	130mm	1 Year Cover & Softgoods, 2 Years Pump	150kg
PTM636631	Air-Lo 5 King Single	2000mm	1070mm	130mm	1 Year Cover & Softgoods 2 Years Pump	170kg



#### AIR-LO 8 MATTRESS REPLACEMENT

The Air-Lo 8 Alternating Mattress Replacement system leads the field in the critical prevention and treatment of pressure areas and skin breakdowns. Air-Lo 8 will support patient loads of up to 200kg with 20 pure TPU transverse air cells which alternately allow pressure to be released across the skin surface. The 8 inch depth allows the Air-Lo 8 to be used without an underlying foam mattress or underlay. When only the very best prevention, therapy and support is required, choose an Air-Lo 8 Mattress Replacement to safely care for patients assessed as high risk or with existing pressure ulcers. Also available in King Single size.

Code	Description	Length	Width	Height	Warranty	S.W.L
PTM636635	Air-Lo 8	2000mm	900mm	210mm	1 Year Cover & Softgoods, 2 Years Pump	200kg
PTM636636	Air-Lo 8 King Single	2000mm	1070mm	210mm	1 Year Cover & Softgoods, 2 Years Pump	200kg





## ASSOCIATED PRODUCTS













































Distributed by:

ABBR02.18 050417

