



DELTA SCOOTER

SMALL PULL APART

Smart Design. Strong Mobility. Travel Ready.

Engineered for convenience and performance, the Aspire Delta Small Scooter combines durable construction with practical features that make everyday mobility easier and more enjoyable both at home and on the go.





Class-leading 150kg SWL

Offers enhanced support and stability for a wider range of users.



Adjustable Tiller for Effortless Control

Easily adjusts to your preferred angle for a more comfortable and personalised ride.



Compact Pull-apart Design

Easily disassembles into five parts for effortless lifting and transport.



Additional Battery Options Available

Standard range battery up to 24km. Extended and airline-approved batteries also available (up to 35 km).




4-wheel Independent Suspension

Enhanced suspension absorbs bumps and uneven surfaces, for a smoother ride.



Tight Turning Radius

Experience seamless navigation in tight spaces, perfect for confident indoor mobility.

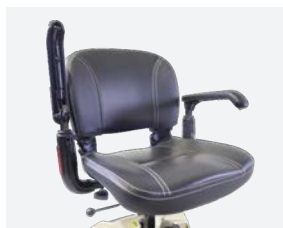

Smooth Transfers, Ultimate Comfort

Featuring a height-adjustable, swivelling padded seat, this scooter ensures easy side transfers every time.



Intuitive Controls

Navigate with confidence using an intuitive dashboard featuring an LED battery display, speed controls, and a horn.

Overall L x W x H	1048 x 550 x 900mm
Seat Size D x W	400 x 450 - 550mm
Seat to Footrest Height	400 - 475mm
Armrest Height	180mm
Min Turning Radius	1200mm
Ground Clearance	70mm
Max Climbing Angle	6°
Max Driving Speed	6.4km/h
Driving Range	24km (35km optional)
Front/Rear Wheels (puncture proof)	190mm
Battery	1 x 24V 20Ah (Lithium)
Product Weight (excl. battery)	35.6kg
Heaviest Part	15.7kg
Battery Weight	3.2kg
Motor	270W, 2 Pole, Brushed
Controller	PG 45A
Brake System	Electromagnetic
Free Wheel Mode	Yes
Suspension	Yes
Charger Size	24V 2A



SCT671600
Red



SCT671610
Champagne

Optional Batteries

SCA669400
Airline-approved battery (11.2km)

SCA669410
Extended range battery (35km)

*Scooter Frame/Motors/Axles/Brakes 2 Years | Battery 6 Months | Electronics 1 Year

2 YEAR WARRANTY*



1300 133 120 | aspirecare.com.au

ASFY74;16.300725