




# LIFECOMFORT HYBRID COMFORT MATTRESS

The Hybrid Comfort is an Acute mattress system offering foam-in-cell technology enabling flexible operation as a high comfort foam static support surface or a powered-up alternating pressure support surface. Specifically re-engineered for increased patient comfort, foam layering now effectively supports patient weight ranging from 30 to 380kg.

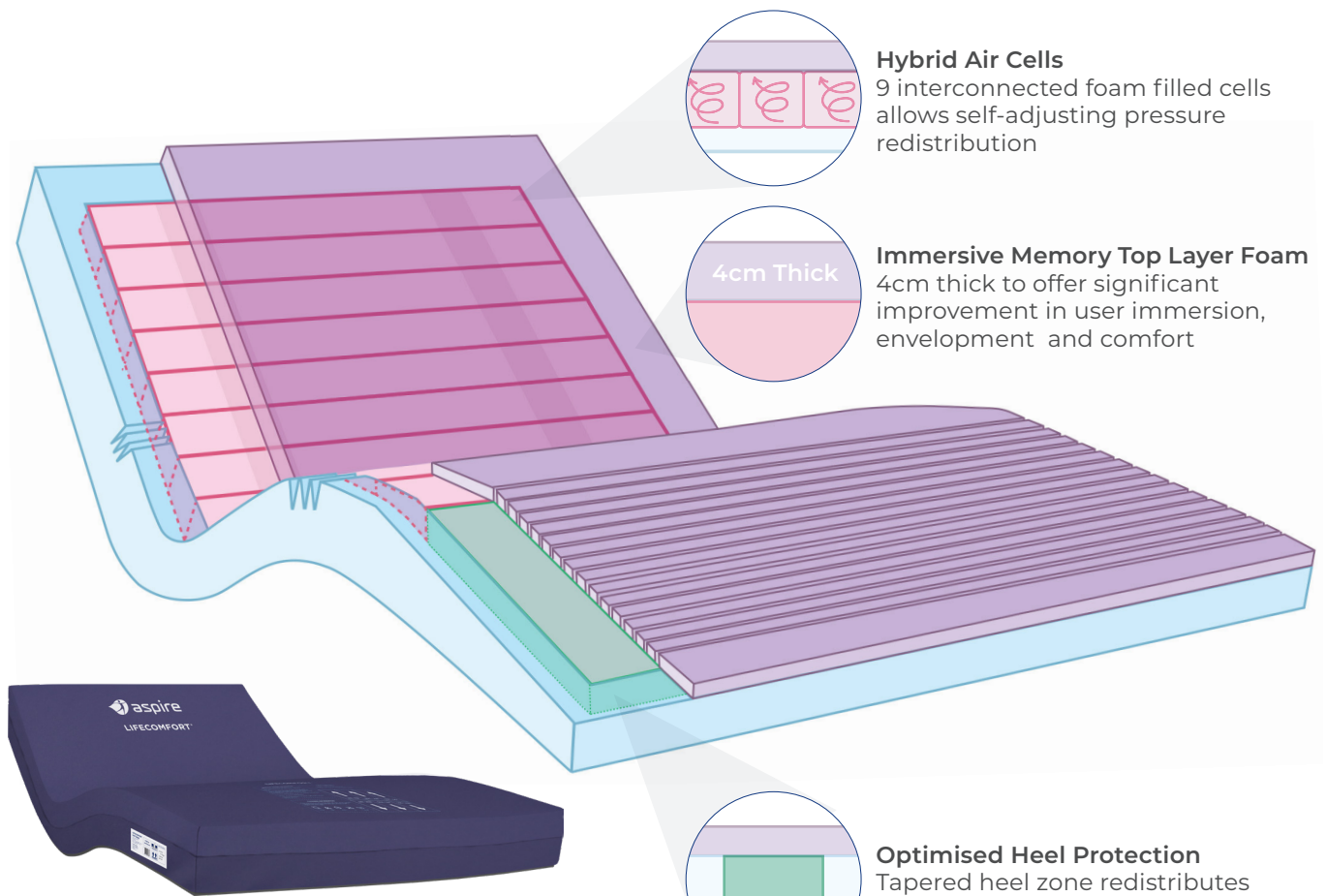

**4-way Stretch Waterproof Cover** For enhanced shear mitigation and friction reduction

**Intelligent Hose Management:** When not in use the air-hose can be stored internally, reducing infection control and OH&S risks

**Integrated Shear Reduction:** Memory foam layer moves independently of the cell surface reducing shear when bed is adjusted

**Innovative Pump Technology:** Reduces noise and vibration for undisturbed sleep



Digital Pump sold separately  
Single / King Single PTM638030\*  
Bariatric - PTM638032\*\*

## APPLICATIONS

- Aged Care
- Home Care
- Ambulant patients
- Palliative Care




**5 YEAR WARRANTY**  
Cover and Foam





**Asset Management** Enabled through unique batch serial numbering

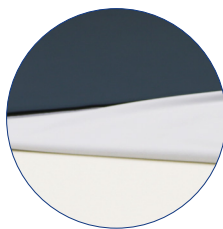

**Hospital-grade Antimicrobial Foam** Ensures superior infection prevention



**Infection Control** Fully welded seams with waterfall zip covers, protect against fluid ingress

## CLINICAL BENEFITS

### Reduced Pressure Injury Risk:

- Medium-High mattress risk rating with 4-way stretch cover providing high MVTR for optimal skin microclimate
- Memory foam layer with integrated shear reduction moves independently during bed adjustment
- Tapered heel zone redistributes pressure; castellated surface reduces shear through wave-like movement
- Enhanced visco-elastic memory foam surface (doubled thickness) improves immersion and envelopment
- Self-adjusting pressure redistribution via 9 interconnected foam cells allowing free air movement
- Hybrid system enables seamless conversion from static to alternating pressure support

### Reduced Clinical Impact:

- Prevents multidisciplinary intervention requirements (medical, nursing, dietetics, AHP)
- Eliminates extended offloading periods (weeks to months for healing)
- Reduces increased support needs: personal care, domestic assistance, meal preparation

### Enhanced Comfort:

- Thickened visco-elastic memory foam improves user immersion and comfort
- Integrated memory foam overlay with self-adjusting pressure redistribution
- Innovative pump technology minimizes noise/vibration for uninterrupted sleep

### Improved Mobility:

- Firm foam bed edges enhance transfer safety and independence



Clinical Assistance Guide

## LIFE CYCLE & COMPLIANCE

**Warranty:** 5 years Foam | 5 years Cover

**Machine Washable Covers:** Up to 95°

**Durability + Performance:** ISO 3385

**Safe Working Load Verified:** AS/NZ 2282

**Flame Retardant Foam + Covers:** AS/NZ 4088.1 Cribs 5

**TGA:** Lifecomfort Mattresses are registered in the ARTG, AUST R ARTG 288486

Code	Description	Length	Width	Height	Effective Weight
PTM638005	Lifecomfort Hybrid Comfort - Single	1980mm	895mm	170mm	30 - 230kg
PTM638006	Lifecomfort Hybrid Comfort - King Single	1980mm	1050mm	170mm	30 - 230kg
PTM638007**	Lifecomfort Hybrid Comfort - Bariatric	1980mm	1050mm	170mm	30 - 380kg

DVA AB14 \*AB19 NDIS 05\_043306211\_0103\_1\_2 \*05\_043306211\_0103\_1\_2 \*\*05\_043306804\_0103\_1\_2