




LIFECOMFORT® ACUTE CARE

The **LIFECOMFORT® ACUTE CARE** mattress has three foam layers including a castellated patient surface providing users with immersive, pressure relieving support and air flow for enhanced skin microclimate. The U-shaped base and sidewall structure provides support to assist with patient transfers as well as creates roll resistance at the sleep surface edges.

- 1. Ultra-soft Heel Protection Zone** improves pressure relief for vulnerable points. By gently tapering the heel area, pressure is redistributed effectively onto the less vulnerable calf area.
- 2. Castellated Surface** collapses in a wave-like movement reducing shear.

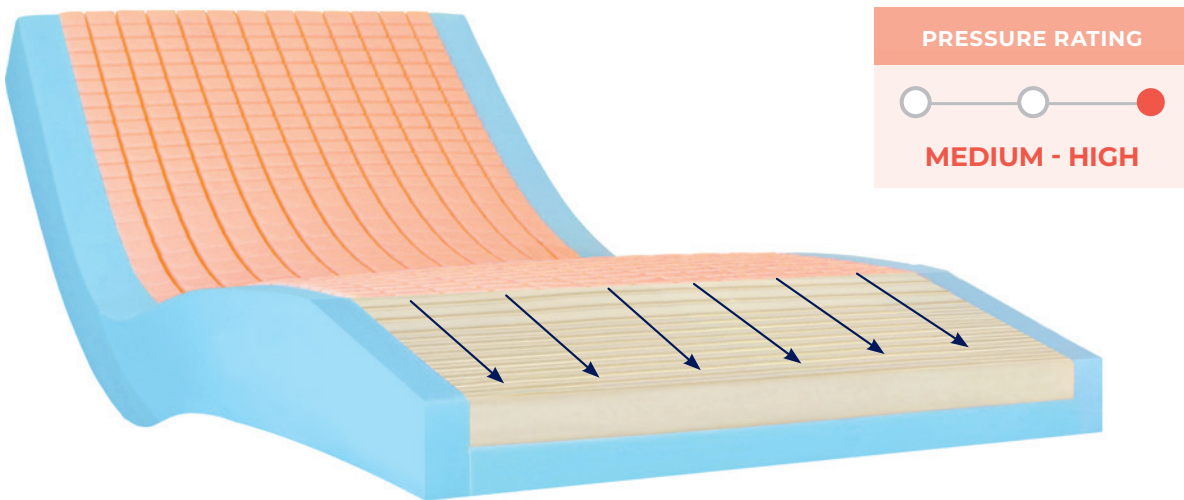

FOAM LAYERS

- IMMERSIVE CASTELLATED SURFACE**
Castellated cells offer greater airflow and transitional support
- COMFORT LAYER**
Visco-elastic comfort and envelopment
- SAFETY & SUPPORT BASE**
Strong U-shaped base for maximum durability and assistive support
- INTEGRATED SIDE WALLS**
Strong rigid side layers for transfer efficiency and safety

APPLICATIONS:

- Hospital
- Aged Care
- Home Care
- Palliative Care



EFFECTIVE WEIGHT RANGE: 30-250KG

BEM046615	Single 1980 x 880 x 150mm
BEM046616	King Single 2030 x 1050 x 150mm



ASFY113.2 051224