

The configurable lift recline chair for unique and changing needs



Personalise your comfort, **optimise your posture!**

PostureFit Configurable Lift Recline Chairs offer adjustable and adaptable seating options designed to provide unparalleled postural support and comfort.





PRESSURE RELIEF SEATING


Hybrid Cushion Technology: Unique foam-in-cell air cushion automatically redistributes pressure and maximises user comfort.

After-Market Cushion Integration: Pre-cut seat structure accommodates third-party cushions such as ROHO without compromising performance or infection control.

Medical Grade Foam: High density foam allows for prolonged periods of sitting.

PRESSURE MAPPING

| 0 | 0 | 0 | 0 | 0 | 0 | 0 | 0 | 0 | 0 | 0 | 0 | 0 | 0 | 0 | |
|-----|----|-----|---|-----|---|-----|----|------|-----|----|----|-----|-----|------|--|
| 5 | 0 | 0 | 0 | ø | 0 | 0 | 0 | 0 | ą | 0 | e. | \$. | 0. | 0 | |
| 2 | 1 | 1 | * | | | . 4 | 2 | 3 | E. | 1 | 10 | | 1 | | |
| 4 | 1 | 5 | | | 5 | 4 | ż. | 3 | - 5 | | | | | | |
| 4 | | 2 | | | | 4 | 2 | 1 | | 1 | | 7 | | 7 | |
| 4 | .6 | | | | 4 | 3 | 12 | | | | | | | 3 | |
| ٩., | 4 | | | | 4 | 4 | 1 | 1 | | | | | | | |
| 6 | 7 | 1 | | 3 | | ñ. | -3 | -15- | | * | | | | 2 | |
| 4 | 4 | * | | | | 14 | -1 | 7 | | | | 1. | . 6 | 7 | |
| έ., | | 8. | 7 | | | | 4 | | | | | | 1 | 7 | |
| 4 | 1 | .t. | | | | | | | | | t | | 4 | 6 | |
| 4 | 1 | | | | | | 7 | | 7 | 18 | 7 | | Ť | 5 | |
| έ., | | | | | | ×. | 5 | | - 2 | | 3. | | | | |
| 5 | | | | -14 | | | 4 | | | | 5. | | | . 61 | |
| 9 | 0 | 0 | œ | 0 | 0 | 0 | 0 | | 2 | a | c | | 0 | 0 | |
| à | à | 6 | 0 | 0 | 0 | ġ. | 0 | ò | 8 | 0 | 0 | ¢ | 0 | 0 | |

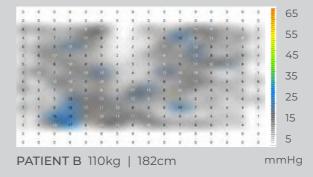

*The pressure mapping scale is set to accurately perform with respect to the clinically recognised 32mmHg average capillary closing pressure. Clinicians should always ensure pressure mapping scale is set to accurately represent true performance.



VERTICAL LIFTING MECHANISM

Direct rise motion offers sitting height adjustment for taller users and creates a stable and reassuring sit to stand experience.



INFINITE POSITIONING

Adjustable Headrest: Electrically adjustable via the handset, the headrest cushion can be tilted forward and backward for optimal support.



Zero Gravity Positioning: Creates a feeling of weightlessness and redistributes pressure off the lower back for pain management.



Oedema Management: Leg elevation above the heart promotes drainage of lower limb swelling, aids circulation and reduces fatigue.



NO MORE KNEELING DOWN OR MANUAL ADJUSTMENTS!

Integrated Footrest Extension: Suits most user heights without the need to alter the legrest or height of the chair.



TOOL-FREE SEATING ADJUSTMENT

Therapists or carers can quickly and easily make adjustments without the need for any tools.



Seat Depth: Backrest can be quickly and easily adjusted to achieve two seat depths with a 50mm range.

Seat Width: Armrests are interchangeable for left/right side, simply swap sides to achieve two seat widths with a 50mm range.

MOVE WITH EASE

Effortless Mobility and Safety: Large directional castors are engaged by activating a rear carer foot pedal. These castors are not in contact with the floor when in a seated position, ensuring stability and preventing accidental tipping or sliding.



 \geq

Ó

Transport mode: For safety, all electrical functions are inactivated while castors are engaged to prevent the chair being activated while mobile.

CLINICAL CONSIDERATIONS

- Medical grade foam for all day sitting
- · No lateral movement in front castors for safe sit-to-stand transfer
- Tool free adjustment
- High level pressure relieving cushion as standard
- · Auto lockout of electronics when transport wheels are engaged



Scan for rescources **Clinical Assistance Guide** and Assessment Form

PERFORMANCE ENHANCING ACCESSORIES



Classic Backrest CHP198290 - Medium CHP198295 - Medium CHP198310 - Medium CHP198291 - Large



Advance Immersion Backrest

CHP198296 - Large



Lateral

Support Blocks

CHP198315

Profiled Headrest CHP198311 - Large



Channel Legrest CHP198320

*All accessories are black, blue highlights are used for identification only

TECHNICAL SPECIFICATIONS

| | SEAT HEIGHT* | SEAT DEPTH | SEAT WIDTH | SEAT TO TOP OF BACK | OVERALL HEIGHT | OVERALL DEPTH | OVERALL WIDTH | DISTANCE FROM WALL | SWL |
|---------------------------|-----------------|---------------|---|------------------------|-------------------|------------------|------------------|-----------------------|-------|
| | | | | | | | | | |
| Small/Medium | 490mm | 475– 525mm | 485– 535mm | 740mm | 1140mm | 900mm | 760– 810mm | 580mm | 160kg |
| CHP198250 CHP198250-UK | 19.3" | 18.7"–20.7" | 19"–21" | 29" | 44.9" | 35.5" | 30"-32" | 22.8" | 25st |
| Tall Medium/ Large | 490mm | 525– 575mm | 535– 585mm | 790mm | 1180mm | 950mm | 820– 870mm | 580mm | 160kg |
| СНР198260 СНР198260-UK | 19.3" | 20.7"-22.6" | 21"-23" | 31" | 46.5" | 37.4" | 32.3"–34.3" | 22.8" | 25st |
| *60mm electric vert | ical heigh | t adjustmer | 2 YEARS WARRANTY Frame & Electronics 1 YEAR WARRANTY Handset & Cables | | | | | | |