

# FL250 FLOORLINE BED

Experience a premium sleeping experience with our FL250 Floorline Bed. Pair it with a selection of soft furnishings to suit the modern home and facility environments.



Mattress sold separately Bed ends sold separately



Ultra Low  
110mm



**10** YEAR WARRANTY  
Frame

**3** YEAR WARRANTY  
Electronics

**1** YEAR WARRANTY  
Handset and Cables



SWL  
250KG

**Low Height:** With its impressive low height of 110mm, injury to the resident can be significantly reduced in the event of a fall from bed, therefore reducing the need for side-rails

**Attendant Control Panel:** Allows caregivers functional control and the ability to lock out certain bed functions from being operated via the handset

**Underbed Clearance:** Market leading 110mm low height whilst maintaining full ground clearance

**ExitSafe:** By pushing the ExitSafe button the FL250 moves to a flat and ergonomic position at the optimum height for resident ingress / exit

**Rapid Hi-lo:** Takes 26 seconds

**SmartDescent:** Bed lowering speed is reduced once the bed reaches 300mm maximising carer potential to remove obstacles as the bed lowers to the floor

**BrakeFirst:** Brakes automatically lock when bed reaches its lowest height. Brakes are located within the footprint of the bed, minimising any trip hazard for the caregiver and resident alike. Brakes are easily accessible at all heights

**Underbed Lighting:** Creates a gentle halo around the bed to reduce falls risk and improve carer safety

**Auto-Regression:** AutoRegression in the backrest and leg section helps to reduce sacral pressure, shear forces and abdominal compression during bed profiling

Infinite Positioning



HI LO



BACK RAISE



LEG RAISE



SEATED



TRENDELENBURG

### Standard Bed Ends Como Radius



Sublime Teak



Milano Walnut

### Side Panels (sold separately)



Sublime Teak




Milano Walnut

**Accessories** - See brochure to view full range of accessories and custom head and footboard options



Full Length Rails

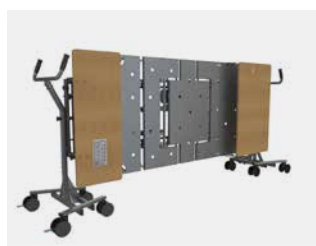

Easy Move Transport System




Split Folding Grab Rails



Self Help Pole & IV Pole

| FL250 Floorline Bed | Single      | King Single |
|---------------------|-------------|-------------|
| Code                | BEB046100   | BEB046150   |
| Mattress Dimensions | 2000 x 900  | 2000 x 1060 |
| Overall Width       | 900         | 1060        |
| Overall Length      | 2088 - 2288 | 2088 - 2288 |
| Length Extension    | 100 - 200   | 100 - 200   |
| Height Range        | 110 - 740   | 110 - 740   |
| Product Weight      | 105kg       | 115kg       |

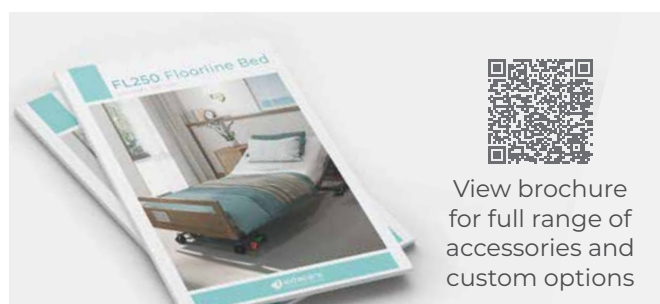

View brochure for full range of accessories and custom options