



PATIENT LIFTERS & SLINGS

Mobile patient lifters and slings designed for a range of environments, from community to facility.



1300 133 120 | aspirecare.com.au

A150F FOLDING FLOOR LIFTER ALUMINIUM

The Aidacare Aspire A150F is a lightweight and portable folding patient lifter engineered to support both healthcare professionals and care recipients, it features a rapid-folding frame for streamlined transport and compact storage, alongside a robust powder-coated aluminium structure ensuring long-term durability.



See page 22 for detailed technical specifications and lifter comparisons



FOLDS FOR TRANSPORT & STORAGE
Folds quickly and conveniently for storage or transport in the boot of the car.



BUILD QUALITY
Strong powder coated aluminium frame with industry leading Linak electronics, Tente castors and warranty provision.



OFFSET LEG DESIGN
Provides additional accessibility around chairs and furniture in the home.



QUICK RELEASE YOKE ATTACHMENT
One touch release and attach mechanism, designed for the busy caregiver. Padded yoke for added patient safety.



INTELLIGENT DIAGNOSTICS
State of the art Linak Jumbo Care control system. Tracks usage and service requirements for safety and compliance.



SMART HANDSET
The latest slimline technology from Linak fits neatly in the hand with an intuitive control panel. Battery and service indicator lights assist with service scheduling.




LEG FUNCTION
Simple foot operation to open or close legs for lifting stability or accessibility.



SAFETY
Emergency stop button, emergency lowering and anti-entrapment technology, with both patient and caregiver safety at heart.



TENTE CASTORS
2 x 75mm front castors for maximum under bed clearance and accessibility.
2 x 100mm rear locking castors for safety and manoeuvrability.

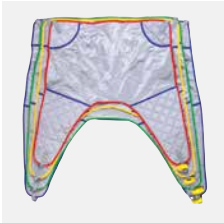

SLING COMPATIBILITY
Compatible with a large range of Aspire Slings, as well as other brands.
See page 24 for the full range of Aspire Slings



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset & Cables

Registered in the
ARTG, AUST R ARTG 224749

200 FLOOR LIFTER STEEL

The Aidacare Aspire 200 Floor Lifter delivers uncompromising reliability for demanding healthcare environments. This steel-constructed workhorse combines intelligent monitoring capabilities with ergonomic design and quick-release mechanisms, ensuring optimal patient care outcomes while maximising caregiver efficiency and operational compliance.



INTELLIGENT DIAGNOSTICS
State of the art Linak Jumbo Care control system. Tracks usage and service requirements for safety and compliance.



QUICK RELEASE YOKE ATTACHMENT
One touch release and attach mechanism, designed for the busy caregiver. Padded yoke for added patient safety.



ERGONOMIC HANDLES
Provide multiple grip points to suit carers of different heights and strengths.



SMART HANDSET
The latest slimline technology from Linak fits neatly in the hand with an intuitive control panel. Battery and service indicator lights assist with service scheduling.



POWERED LEG FUNCTION
Opens and closes the legs for lifting stability or accessibility.



SAFETY
Critical features include emergency stop button, emergency lowering and anti-entrapment technology, with both patient and caregiver safety at heart.



TENTE CASTORS
2 x 75mm front castors for maximum under bed clearance and accessibility.
2 x 100mm rear locking castors for safety and manoeuvrability.



BROAD SLING COMPATIBILITY
Works with Aspire and most third-party slings.
See page 24 for the full range of Aspire Slings



BUILD QUALITY
Strong powder coated aluminium frame with industry leading Linak electronics, Tente castors and warranty provision.



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset & Cables

Registered in the
ARTG, AUST R ARTG 224749

A205 FLOOR LIFTER ALUMINIUM

The Aidacare Aspire A205 Floor Lifter combines portability with professional performance through its lightweight aluminium construction and tool-free disassembly design. Featuring three adjustable height settings, intelligent Linak diagnostics, and powered leg functionality, this versatile lifter seamlessly adapts to home, community, and facility environments.



TRANSPORT OR STORAGE
Is quick and convenient with tool free disassembly.



BUILD QUALITY
Strong powder coated aluminium frame with industry leading Linak electronics, Tente castors and warranty provision.



ADJUSTABLE HEIGHT
Three adjustable height settings provides the caregiver with the flexibility to adjust the height range to suit the specific lift and patient scenario.



QUICK RELEASE YOKE ATTACHMENT
One touch release and attach mechanism, designed for the busy caregiver. Padded yoke for added patient safety.



INTELLIGENT DIAGNOSTICS
State of the art Linak Jumbo Care control system. Tracks usage and service requirements for safety and compliance.




SMART HANDSET
The latest slimline technology from Linak fits neatly in the hand with an intuitive control panel. Battery and service indicator lights assist with service scheduling.




POWERED LEG FUNCTION
Opens and closes the legs for lifting stability or accessibility.



SAFETY
Emergency stop button, emergency lowering and anti-entrapment technology, with both patient and caregiver safety at heart.

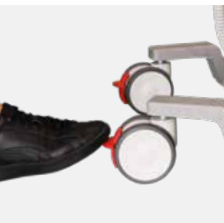

TENTE CASTORS
2 x 75mm front castors for maximum under bed clearance and accessibility.
2 x 100mm rear locking castors for safety and manoeuvrability.

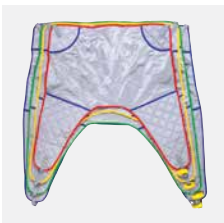

BROAD SLING COMPATIBILITY
Works with Aspire and most third-party slings.
See page 24 for the full range of Aspire Slings



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset & Cables

Registered in the
ARTG, AUST R ARTG 224749

AURA FLOOR LIFTER ALUMINIUM


Engineered for demanding commercial healthcare environments, the Aura Floor Lifter delivers unmatched quality, safety, and ease of use. The compact base design is easy to move, offering market leading read and lift ranges.



TiMOTION

Bonus Second TiMotion Battery & External Charging Cradle
Doubles the operational capacity and reduces lifter downtime

Product code: MBT401050

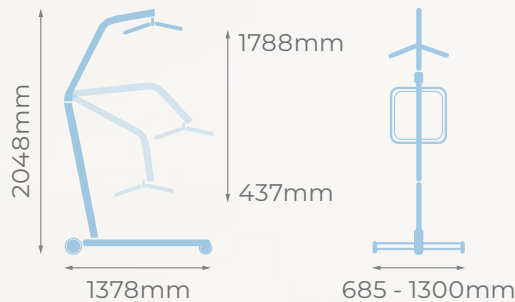

See page 22 for detailed technical specifications and lifter comparisons

AS/NZS ISO 10535:2024 COMPLIANT




38KG LIGHTWEIGHT ALUMINIUM FRAME

Supports a 250kg safe working load, even at maximum extension, making standout in its class.



EMERGENCY LOWERING & STOP FUNCTIONS

Built-in safety for power or system failures

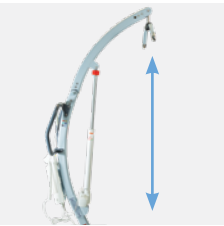

CURVED ERGONOMIC FRAME & SLING BAR

Enhances reach and user headroom, ideal for transfers from king single beds or around bulky furniture in tight spaces.



DURABLE LEG SHROUDS

Replaceable covers protect components when going under beds and furniture, maintains product appearance.



TRUE FLOOR LIFT CAPABILITY

Safe, solo-carer lifts from the floor, essential in falls recovery and emergency care.



COLOUR-CODED SLING LOCKS

Quick visual checks for safer, more confident transfers.



POWERED LEG OPENING & TWIN CASTORS

Effortlessly positions around wheelchairs and beds with smooth manoeuvrability and stability.



SMART HANDSET

With battery levels and service status indicators, supports efficient product servicing

ACCESSORIES



4 POINT SLING BAR

Optional accessory increases and precision stability for bariatric users.
Product code: LSS361000



INTEGRATED WEIGH SCALE

Optional accessory supports efficient weight assessment.
(Available Q3 2025)



ELECTRIC PIVOT FRAME

Optional accessory enables smoother transfers via natural positioning and minimises user strain on joints
(Available Q3 2025)



BROAD SLING COMPATIBILITY

Works with Aspire and most third-party slings.
See page 24 for the full range of Aspire Slings



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset & Cables

Registered in the
ARTG, AUST R ARTG 224749

A320 FLOOR LIFTER

BARIATRIC ALUMINIUM

The Aidacare Aspire A320 Bariatric Floor Lifter excels with its high 320kg safe working load through its reinforced aluminium construction and broad-base design for enhanced patient clearance. Equipped with folding support arms, adjustable height settings and four-point yoke system, it delivers uncompromising stability for complex bariatric transfers.



INTELLIGENT DIAGNOSTICS
State of the art Linak Jumbo Care control system. Tracks usage and service requirements for safety and compliance.



BUILD QUALITY
Strong powder coated aluminium frame with industry leading Linak electronics, Tente castors and warranty provision.



ADJUSTABLE HEIGHT
Three adjustable height settings provides the caregiver with the flexibility to adjust the height range to suit the specific lift and patient scenario.



QUICK RELEASE 4 POINT YOKE
Padded 4 point yoke for added patient safety. One touch release and attach mechanism, designed for the busy caregiver.



ERGONOMIC HANDLES AND SUPPORT ARMS
Provide multiple grip points to suit carers of different heights and strengths. Folding support arms assist with stability or mobility but can also be folded out of the way.



SMART HANDSET
The latest slimline technology from Linak fits neatly in the hand with an intuitive control panel. Battery and service indicator lights assist with service scheduling.




TENTE
2 x 100mm front castors for reduced rolling resistance under high load.
2 x 125mm rear locking castors for safety and manoeuvrability.



SAFETY
Critical features include emergency stop button, emergency lowering and anti-entrapment technology, with both patient and caregiver safety at heart.



POWERED LEG FUNCTION
Opens and closes the legs for lifting stability or accessibility.

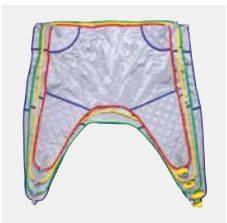

SLING COMPATIBILITY
Compatible with a large range of Aspire Slings, as well as other brands.
See page 24 for the full range of Aspire Slings



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset & Cables

Registered in the
ARTG, AUST R ARTG 224749

200S STEEL

The Aidacare Aspire 200S Standing Lifter empowers patient mobility and rehabilitation through its specialized bullhorn lifting boom with six attachment points. This steel-constructed mobility solution features interchangeable clip sling compatibility, enabling optimal standing transfers.



INTELLIGENT DIAGNOSTICS
State of the art Linak Jumbo Care control system. Tracks usage and service requirements for safety and compliance.



BUILD QUALITY
Strong powder coated aluminium frame with industry leading Linak electronics, Tente castors and warranty provision.



BULLHORN LIFTING BOOM
Offset patient handles provide stability and 6 sling attachment points give the carer flexibility to perform the optimal lift in any scenario.



ADJUSTABLE PADDED LEGRESTS
Provide comfort and support during transfers and can be adjusted to suit the patient.



ERGONOMIC HANDLES
Provide multiple grip points to suit carers of different heights and strengths.



ERGONOMIC HANDSET
The latest slimline technology from Linak fits neatly in the hand with an intuitive control panel. Battery and service indicator lights assist with service scheduling.



TENTE
2 x 75mm front castors for maximum under bed clearance and accessibility.
2 x 75mm rear locking castors for safety and manoeuvrability.



SAFETY
Emergency stop button, emergency lowering and anti-entrapment technology, with both patient and caregiver safety at heart.



POWERED LEG FUNCTION
Opens and closes the legs for lifting stability or accessibility.



SLING COMPATIBILITY
Compatible with a large range of Aspire Slings, as well as other brands.
See page 24 for the full range of Aspire Slings

ACCESSORY



INTERCHANGEABLE CLIP SLING BAR
One-for-one swap of existing sling bar converts stander to accept clip slings only (retro-fit).
Product code: LSS390655



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset & Cables

Registered in the
ARTG, AUST R ARTG 224749

A200S STANDING LIFTER

ALUMINIUM

The Aidacare Aspire A200S Standing Lifter optimises patient comfort through its depth-adjustable yoke system and lightweight aluminium construction. Designed for versatile home-to-facility use, this compact mobility solution features adjustable padded leg rests and intelligent diagnostics, enabling personalised standing assistance across diverse care environments.




INTELLIGENT DIAGNOSTICS

State of the art Linak Jumbo Care control system. Tracks usage and service requirements for safety and compliance.



BUILD QUALITY

Strong powder coated aluminium frame with industry leading Linak electronics, Tente castors and warranty provision.

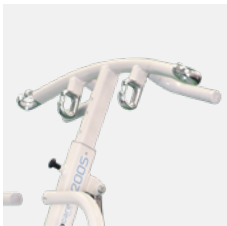

DEPTH ADJUSTABLE YOKE

Yoke is depth adjustable for optimal patient fit.



ADJUSTABLE PADDED LEGRESTS

Provide comfort and support during transfers and can be adjusted to suit the patient.




ERGONOMIC HANDLES

Provide multiple grip points to suit carers of different heights and strengths.



SMART HANDSET

The latest slimline technology from Linak fits neatly in the hand with an intuitive control panel. Battery and service indicator lights assist with service scheduling.

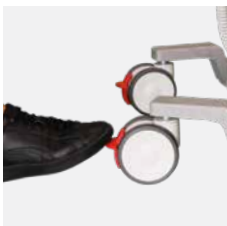

TENTE

2 x 75mm front castors for maximum under bed clearance and accessibility.
2 x 75mm rear locking castors for safety and manoeuvrability




SAFETY

Critical features include emergency stop button, emergency lowering and anti-entrapment technology, with both patient and caregiver safety at heart.



POWERED LEG FUNCTION

Opens and closes the legs for lifting stability or accessibility

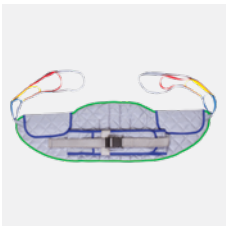

SLING COMPATIBILITY

Compatible with a large range of Aspire Slings, as well as other brands.
See page 24 for the full range of Aspire Slings



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset & Cables

Registered in the
ARTG, AUST R ARTG 224749

NOVA STANDING LIFTER ALUMINIUM


The Nova Standing Lifter combines a lightweight aluminium frame with a high Safe Working Load (SWL), delivering a leading solution for commercial healthcare environments. Its modern design is ideal for facilities seeking a safer and more efficient way to perform standing lifts, even in tight spaces and for users of varying sizes.



0% TiMOTION

Bonus Second TiMotion Battery & External Charging Cradle
Doubles the operational capacity and reduces lifter downtime

Product code: MBT401050

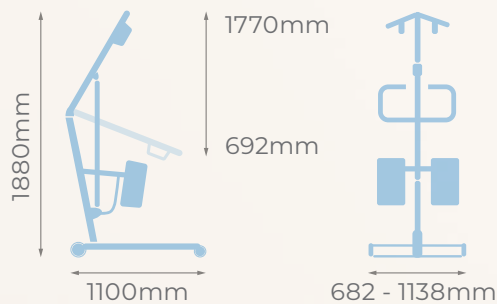

See page 22 for detailed technical specifications and lifter comparisons

AS/NZS ISO 10535:2024 COMPLIANT




39.6KG ALUMINIUM FRAME WITH COMPACT BASE

Supports a 250kg safe working load, reduces carer effort in facilities with crowded space.



MOTORBIKE-STYLE HANDLES

Enhances user confidence and stability during transfers.

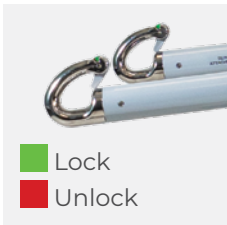

COLOUR-CODED SLING LOCKS

Quick visual checks for safer, more confident transfers.



EMERGENCY LOWERING & STOP FUNCTIONS

Built-in safety for power or system failures.



OVER ONE METRE LIFTING RANGE

Supports a broad spectrum of user mobility needs, ensures safe, efficient clinical transfers.



DURABLE LEG SHROUDS

Replaceable covers protect components when going under beds and furniture, maintains product appearance.



ADJUSTABLE KNEE PADS & LEG STRAP

Comfort-focused design with infection control surfaces.



SMART HANDSET

With battery levels and service status indicators, supports efficient product servicing.



POWERED LEG OPENING & LOW-MOUNTED FOOTPLATE

Allows smooth positioning around wheelchairs and chairs, improving both access and safety.



TWIN FRONT AND REAR LOCKING CASTORS

Smaller front castors enhance underbed clearance and manoeuvrability, with rear locks for added control.

ACCESSORIES



INTERCHANGEABLE CLIP SLING BAR

One-for-one swap allows compatibility with clip slings for flexible lift handling.

Product Code: LSS390660



BROAD SLING COMPATIBILITY

Works with Aspire and most third-party slings.

See page 24 for the full range of Aspire Slings



5 YEAR WARRANTY
Frame

3 YEAR WARRANTY
Electronics

1 YEAR WARRANTY
Battery, Handset & Cables

Registered in the
ARTG, AUST R ARTG 224749

STAR MOVER


ALUMINIUM

The Aspire Star Mover (Patient Mover), provides a compact solution for transferring users from a bed or chair in comfort and safety through its unique four-point contact system. Featuring adjustable knee pads, dual grab bars and wide seat entry it enables safe lateral transfers without lifting across homecare, acute care and aged care environments.



CLIP SLING ATTACHMENT POINTS

Allows a standing sling (optional accessory) with clips to be used for additional support where required.

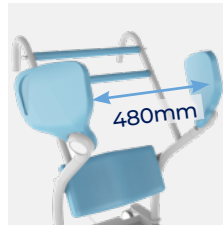


WIDE SEAT ENTRY SPACE

Wide seat entry suits a larger range of body shapes.



FOUR USER CONTACT POINTS

Enhance user confidence during transfers via foot platform, knee pad, seat pads and grip handles, delivers superior stability and comfort for safer transfers.

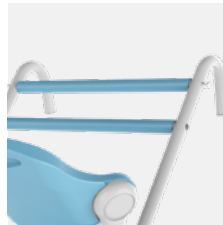

SOFT TOUCH GRAB BARS

Two bars to suit multiple user heights provide a comfortable patient support for improved stability and confidence during transfers.



HEIGHT ADJUSTABLE SOFT KNEE PAD

Provides comfortable counter balance stability during transfers.



LOW FOOTPLATE HEIGHT

Low footplate height makes patient entry and exit easier.




LOCKABLE 75MM CASTORS

The platform remains stationary while a user is transferring, while minimising the leg height for sufficient clearance beneath beds and furniture.



LARGE CASTORS

2 x large 100mm locking tente castors.

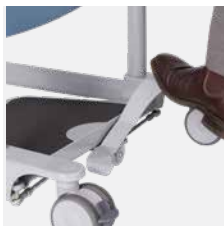

WIDTH ADJUSTABLE LEGS

Activated via 2 foot pedals, provide the ability to manoeuvre closer to furniture, toilets or beds for safer access.



EASY TO CLEAN SURFACES

For effective and hassle-free infection control.



SWL
185kg

2 YEARS
WARRANTY

AS/NZS Approved

Registered in the ARTG, AUST R ARTG 265455

WEIGH SCALE UPGRADES

Our range of weigh scale models are designed to support precise and efficient patient care. Each unit enables caregivers to obtain accurate weight measurements while the patient is safely suspended in a sling, eliminating the need for additional transfers. Featuring a clear, easy-to-read LCD display, these scales streamline clinical workflows and support better-informed care decisions across a variety of care environments.



In-Line Weigh Scale [LSS390680](#)

DIVERSE COMPATIBILITY
Designed to be fitted on A150F, A205, 200 and A205

HIGH WEIGHT CAPACITY
Large weight capacity of up to 320kg.




Integrated Weigh Scale [LSS390675](#)

DIVERSE COMPATIBILITY
Designed to be fitted on A150F, 200 and A205, in use with general-purpose slings featuring loop attachments.


OPTIMISED LIFTING PERFORMANCE
By reducing height loss compared to conventional hoist scale systems those unable to use conventional scales.

INFECTION CONTROL
With non-absorbent polyurethane padding designed for easy cleaning.

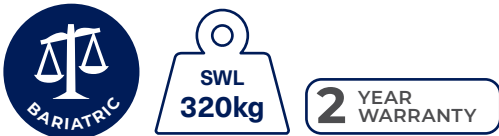

PIVOT FRAMES UPGRADES

Engineered for effortless use with clip slings, this system simplifies the attachment process to enhance operational efficiency. The rotating frame makes both sling application and patient positioning easier and more precise. They are compatible with a wide range of lifters, supporting seamless integration into diverse care environments.

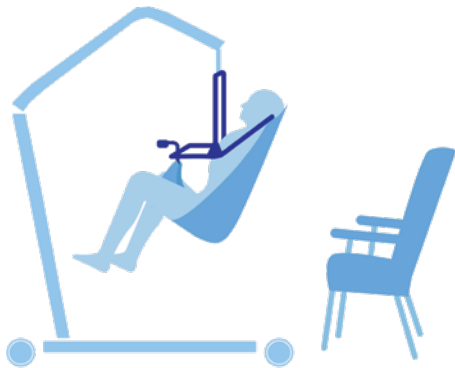

Manual Pivot Frame [LSS390670](#)

DIVERSE COMPATIBILITY
Designed to be fitted on 200, A205 and A320.

IMMEDIATE USAGE
Battery-free design allows immediate usage at all times with minimal maintenance required.



Powered Pivot Frame [LSS390665](#)

DIVERSE COMPATIBILITY
Designed to be fitted on 200, A205 and A320.

MINIMISES CARER EFFORT
Assists staff to conduct repetitive transfers in institutional care environments.

INTEGRATED ELECTRONICS
Seamlessly integrates with lifter control system for use with one handset.



COMPARE OUR LIFTER RANGE



MODEL	A150F	200	A205	AURA
Lifter Type	Folding Floor Lifter	Floor Lifter	Floor Lifter	Floor Lifter
Frame Material	Aluminium	Steel	Aluminium	Aluminium
Product Code	LSS390600	LSS390620	LSS390610	LSS390615
Safe Working Load – SWL	150kg	200kg	205kg	250kg
Product Weight	32.9kg	40kg	38.1kg	38kg
Base Length	1280mm (1220mm folded)	1320mm	1238mm	1378mm
External Base Width	660 - 1040mm (615mm folded)	665 - 1120mm	651 - 1173mm	682 - 1300mm
Internal Base Width	580 - 950mm	580 - 1040mm	547 - 1073mm	518 - 1150mm
Base Width Adjustment	Manual	Electric	Electric	Electric
Base Height / Under Bed Clearance	133mm	110mm	110mm	131mm
Minimum Height of Spreader Bar **	428mm	520mm	415 - 515mm	437mm
Maximum Height of Spreader Bar**	1640mm	1720mm	1700 - 1800mm	1788mm
Under Base Height	72mm	50mm	50mm	49mm
Maximum Overall Height	1873mm (395mm folded)	1935mm	2040mm	2048mm
Front Castor Diameter	75mm	75mm	75mm	100mm
Rear Locking Castor Diameter	100mm	100mm	100mm	125mm
Turning Radius	1365mm	1398mm	1332mm	1453mm
Heaviest Part Weight – Base Assembly	32.9kg	23.0kg	19.6kg	19.6kg
Handle Height	1170mm	1220mm	1085 - 1185mm	883 - 1305mm
Warranty – Frame	5 Years	5 Years	5 Years	5 Years
Warranty – Electronics	3 Years	3 Years	3 Years	3 Years
Warranty – Battery Handset & Cables	1 Year	1 Year	1 Year	1 Year

ACCESSORIES

In-line Weigh Scale – SWL 320kg	LSS390680	LSS390680	LSS390680	Not Available
Integrated Weigh Scale – SWL 220kg	LSS390675	LSS390675	LSS390675	Available in Q3 2025
Manual Pivot Frame – SWL 227kg	Not Available	LSS390670	LSS390670	Not Available
Powered Pivot Frame – SWL 227kg	Not Available	LSS390665	LSS390665	Available in Q3 2025
4 Point Sling Bar	Not Available	Not Available	Not Available	LSS361000
Clip Sling Bar Attachment	Not Available	Not Available	Not Available	Not Available
Leg Protection Strips / Shrouds*	Not Available	GSP390816	GSP390815	GSP242902*

**Three mast height settings are possible hence the range of minimum and maximum heights.




A320	200S	A200S	NOVA	STAR
Floor Lifter	Standing Lifter	Standing Lifter	Standing Lifter	Patient Mover
Aluminium	Steel	Aluminium	Aluminium	Aluminium
LSS390630	LSS390650	LSS390640	LSS390652	LSS390580
320kg	200kg	200kg	250kg	185kg
63kg	40kg	33kg	39.6kg	40kg
1500mm	1277mm	1070mm	1100mm	925mm
765 - 1190mm	656 - 1203mm	590 - 1040mm	682 - 1138mm	520mm
640 - 1065mm	561 - 1103mm	490 - 940mm	519 - 986mm	740mm
Electric	Electric	Electric	Electric	Manual
155mm	110mm	110mm	136mm	60mm
540 - 640mm	857mm	880mm	692mm	-
1840 - 1940mm	1777mm	1415mm	1770mm	-
155mm	50mm	50mm	53mm	-
2155mm	1832mm	1480mm	1800mm	1030mm
100mm	75mm	75mm	100mm	100mm
125mm	100mm	75mm	125mm	100mm
1595mm	1361mm	1160mm	1245mm	-
38.9kg	19.6kg	19kg	17.w6kg	-
1200 - 1300mm	1045mm	1160mm	810 - 1095mm	1030mm
5 Years	5 Years	5 Years	5 Years	2 Years
3 Years	3 Years	3 Years	3 Years	-
1 Year	1 Year	1 Year	1 Year	-

LSS390680	Not Available	Not Available	Not Available	Not Available
LSS390675	Not Available	Not Available	Not Available	Not Available
LSS390670	Not Available	Not Available	Not Available	Not Available
LSS390665	Not Available	Not Available	Not Available	Not Available
Standard Inclusion	Not Available	Not Available	Not Available	Not Available
Not Available	LSS390655	Not Available	LSS390660	Not Available
GSP390817	Not Available	GSP390815	GSP242806*	Not Available

To view our comprehensive list of spare parts visit aidacare.com.au

ASPIRE SLINGS


The Aidacare Aspire range of slings is designed to meet patient needs in the home, aged care and hospital environments. Manufactured under the highest quality control standards, the extremely durable slings come in a range of fabrics and sizes to suit the individual users and caregivers' needs.

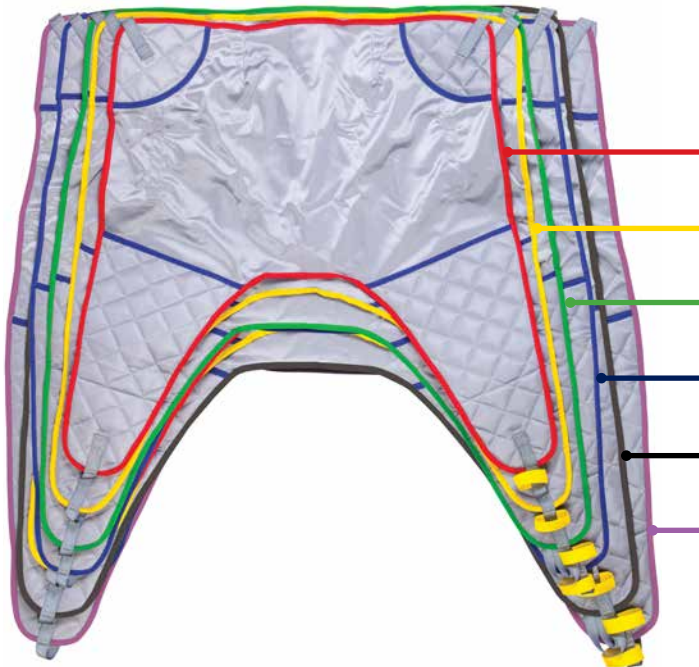


Size	User Girth	User Height	SWL
SMALL	65 – 95cm	120 – 160cm	200kg
MEDIUM	77 – 110cm	140 – 175cm	200kg
LARGE	98 – 135cm	158 – 185cm	300kg
X LARGE	118 – 157cm	165 – 193cm	300kg
XX LARGE	155 – 184cm	170 – 193cm	300kg
XXX LARGE	180 – 210cm	170 – 193cm	300kg

Sling Fabric Options


Your choice of fabric will also depend on the task at hand and the physical characteristics of the user.




POLYESTER
A smooth, knitted polyester fabric of solid consistency for multi-disciplinary use



MESH
A versatile, breathable and durable polyester mesh. Ideal for use in wet environments

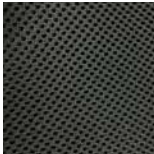

SUPER SOFT
A “super soft” breathable fabric, which is padded for comfort during longer periods of transfer



PARASILK
A smooth, lightweight, low friction polyester fabric. Ideal for easy positioning and optimum comfort

Aspire Sling Type	Polyester	Mesh	Supersoft	Parasilk
Deluxe General Purpose Sling	●	●		
Deluxe General Purpose Sling with Head Support	●	●		
Deluxe Hammock Sling with Head Support		●	●	
Deluxe Low Friction Sling				●
Deluxe Hygiene Access Sling	●			
Deluxe Hygiene Access Sling with Head Support	●			
Deluxe Standing Sling	●			
Hybrid Patient Sling	●			



DURABILITY ASSURANCE
Loop webbing delivers structural integrity with a breaking strain exceeding one tonne, anti-tangle engineering that prevents cork-screwing during cleaning procedures




COMPREHENSIVE TRACEABILITY SYSTEM
Individual embroidered identification labels with unique serial number enables complete product traceability and support regulatory compliance.



PREMIUM MANUFACTURING STANDARDS
Manufactured in the UK utilising European materials, computer guided stitching technology and CNC fabric cutting.



SERVICE & MAINTENANCE DOCUMENTATION
Service record labels facilitates accurate maintenance tracking and service history documentation

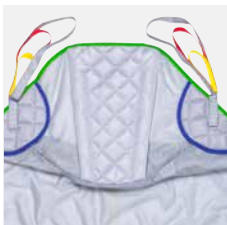

ENHANCED USER COMFORT
Strategically integrated padding in critical weight-bearing zones, optimises pressure distribution and reduces discomfort during patient handling.



COLOUR CODED SIZING
Reinforced colour coded edging for quick size identification, streamlining clinical workflow and reduces selection errors.



MANOEUVRING HANDLES
Thoughtfully positioned handles provide carers with secure grip locations to facilitate safe and controlled patient transfers. (excluding Deluxe Low Friction Slings)



PATIENT DIGNITY
Discreet positioning options allow for external leg strap placement, maintaining user modesty during transfer procedures. (applicable models)



ADVANCED PATIENT SAFETY
Innovative velcro and buckle combination systems ensure optimal positioning security and safety. (applicable models)



EASY CLEAN
Machine washable up to 85°C


Aidacare Sling Auditing

Aidacare's Sling Audit Program can assist you to ensure your sling inventory is in safe working order at all times.

Periodic sling audits are now a mandatory accreditation requirement. Ask us how we can help you manage the risk in this area of your business.

1300 133 120 | aidacare.com.au



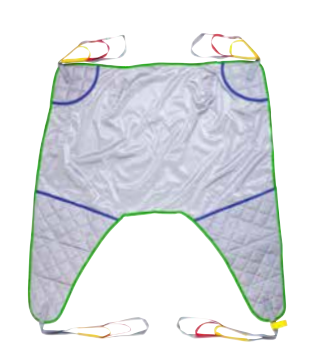

Aspire Deluxe General Purpose Sling

Code	Description	Code	Description	Edging	Size	SWL
LSS390700	Polyester – Loop	LS390716	Mesh – Loop	●	Small	200kg
LSS390705	Polyester – Loop	LS390717	Mesh – Loop	●	Medium	200kg
LSS390710	Polyester – Loop	LS390718	Mesh – Loop	●	Large	300kg
LSS390715	Polyester – Loop	LS390719	Mesh – Loop	●	X Large	300kg

Deluxe General Purpose Sling is registered in the ARTG, AUST R ARTG 224718



Clip Option

Aspire Deluxe General Purpose Sling with Head Support

Code	Description	Code	Description	Edging	Size	SWL
LSS390740	Polyester – Loop	LSS390756	Mesh – Loop	●	Small	200kg
LSS390745	Polyester – Loop	LSS390757	Mesh – Loop	●	Medium	200kg
LSS390750	Polyester – Loop	LSS390758	Mesh – Loop	●	Large	300kg
LSS390755	Polyester – Loop	LSS390759	Mesh – Loop	●	X Large	300kg
LSS3907551	Polyester – Loop	LSS3907591	Mesh – Loop	●	XX Large	300kg
LSS3907552	Polyester – Loop	LSS3907592	Mesh – Loop	●	XXX Large	300kg
LSS390900	Polyester – Clip	LSS390916	Mesh – Clip	●	Small	200kg
LSS390905	Polyester – Clip	LSS390917	Mesh – Clip	●	Medium	200kg
LSS390910	Polyester – Clip	LSS390918	Mesh – Clip	●	Large	300kg
LSS390915	Polyester – Clip	LSS390919	Mesh – Clip	●	X Large	300kg

Deluxe General Purpose Sling is registered in the ARTG, AUST R ARTG 224718



Aspire Deluxe Hammock Sling with Head Support

Code	Description	Code	Description	Edging	Size	SWL
LSS390800	Supersoft – Loop	LSS390860	Mesh – Loop	●	Small	200kg
LSS390805	Supersoft – Loop	LSS390865	Mesh – Loop	●	Medium	200kg
LSS390810	Supersoft – Loop	LSS390870	Mesh – Loop	●	Large	300kg
LSS390815	Supersoft – Loop	LSS390875	Mesh – Loop	●	X Large	300kg

Deluxe Hammock Sling is registered in the ARTG, AUST R ARTG 224720



Clip Option

Aspire Deluxe Low Friction Sling

Code	Description	Code	Description	Edging	Size	SWL
LSS390880	Parasilk – Loop	LSS390960	Parasilk – Clip	●	Small	200kg
LSS390885	Parasilk – Loop	LSS390965	Parasilk – Clip	●	Medium	200kg
LSS390890	Parasilk – Loop	LSS390970	Parasilk – Clip	●	Large	300kg
LSS390895	Parasilk – Loop	LSS390975	Parasilk – Clip	●	X Large	300kg

Deluxe Low Friction Sling (Loop) is registered in the ARTG, AUST R ARTG 224720
Deluxe Low Friction Sling (Clip) is registered in the ARTG, AUST R ARTG 224718

Aspire Deluxe Hygiene Access Sling

Code	Description	Edging	Size	SWL
LSS390760	Polyester – Loop	●	Small	200kg
LSS390765	Polyester – Loop	●	Medium	200kg
LSS390770	Polyester – Loop	●	Large	300kg
LSS390775	Polyester – Loop	●	X Large	300kg
LSS390776	Polyester – Loop	●	XX Large	300kg
LSS390777	Polyester – Loop	●	XXX Large	300kg

Deluxe Hygiene Access Sling is registered in the ARTG, AUST R ARTG 224719



Aspire Deluxe Hygiene Access Sling with Head Support

Code	Description	Code	Description	Edging	Size	SWL
LSS390780	Polyester – Loop	LSS390940	Polyester – Clip	●	Small	200kg
LSS390785	Polyester – Loop	LSS390945	Polyester – Clip	●	Medium	200kg
LSS390790	Polyester – Loop	LSS390950	Polyester – Clip	●	Large	300kg
LSS390795	Polyester – Loop	LSS390955	Polyester – Clip	●	X Large	300kg


Deluxe Hygiene Access Sling is registered in the ARTG, AUST R ARTG 224719



Aspire Deluxe Standing Sling

Code	Description	Code	Description	Edging	Size	SWL
LSS390720	Polyester – Loop	LSS390920	Polyester – Clip	●	Small	200kg
LSS390725	Polyester – Loop	LSS390925	Polyester – Clip	●	Medium	200kg
LSS390730	Polyester – Loop	LSS390930	Polyester – Clip	●	Large	300kg
LSS390735	Polyester – Loop	LSS390935	Polyester – Clip	●	X Large	300kg
LSS390736	Polyester – Loop			●	XX Large	300kg
LSS390737	Polyester – Loop			●	XXX Large	300kg

Deluxe Standing Sling is registered in the ARTG, AUST R ARTG 224698

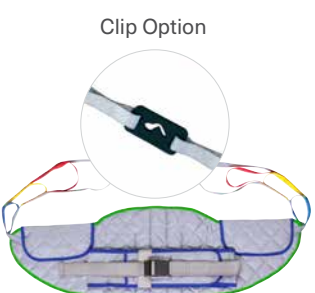

Clip Option

Aspire Hybrid Patient Sling

Code	Description	Edging	Size	SWL
LSS390830	Polyester – Loop	●	Small	200kg
LSS390835	Polyester – Loop	●	Medium	200kg
LSS390840	Polyester – Loop	●	Large	300kg
LSS390845	Polyester – Loop	●	X Large	300kg

Hybrid Patient Sling is registered in the ARTG, AUST R ARTG 306027



Aspire Single Use Patient Sling

Code	Description	Edging	Size	SWL
LSS391100	Non Woven Fabric	●	Small	200kg
LSS391110	Non Woven Fabric	●	Medium	200kg
LSS391120	Non Woven Fabric	●	Large	300kg
LSS391130	Non Woven Fabric	●	X Large	300kg

Single Use Patient Sling is registered in the ARTG, AUST R ARTG 306027



