

# NOVA STANDING LIFTER

The Nova Standing Lifter is a compact solution that weighs just 39.6 kg, yet offers a leading safe working load of 250 kg. Its modern design is ideal for facilities seeking a safer, more efficient way to perform standing lifts, even in tight spaces, and supports users of varying sizes. Designed for individuals with partial weight-bearing ability, the Nova Standing Lifter provides the necessary support to assist users in transitioning safely from a seated to a standing position.



## 0% TiMOTION®

### Bonus Second TiMotion Battery & External Charging Cradle

Allows off-unit charging to keep the lifter operational at all times

Product Code: MBT401050



**5** YEAR WARRANTY  
Frame

**3** YEAR WARRANTY  
Electronics

**1** YEAR WARRANTY  
Battery, Handset & Cables






### 39.6kg Aluminium Frame with Compact Base Design

Supports a 250kg safe working load, reducing carer effort in facilities, particularly where space is limited



### Motorbike-Style Handles

Provide enhanced user confidence and stability during standing transfers

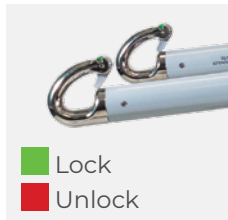

### Colour-coded Sling Locks

Patented Design enables quick visual checks for safer, more confident transfers



### Emergency Lowering & Stop Functions

Built-in features ensure safety in the event of power or system failure



### Over One Metre Lifting Range

Accommodates a broad range of user mobility levels, supporting safe and efficient clinical transfers



### Durable Leg Shrouds

Replaceable covers protect key components when operating under beds and furniture, preserving appearance and function



### Adjustable Knee Pads & Leg Strap

Designed for user comfort, with surfaces that support effective infection control



### Smart Handset

Displays battery level and service status indicators to support proactive maintenance and efficient servicing



### Powered Leg Opening & Low-Mounted Footplate


Facilitates smooth positioning around wheelchairs and chairs, improving access and safety



### 100mm Twin Front Castors & 125mm Twin Rear Castors

Smaller front castors enhance underbed clearance and manoeuvrability, with rear locks for added safety and control

## ACCESSORIES

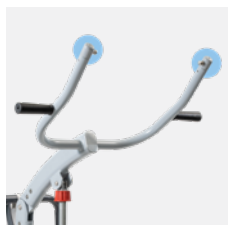

### Interchangeable Clip Boom

Optional accessory that enables compatibility with both loop and clip slings, offering greater flexibility in lift handling

LSS390660



### Broad Sling Compatibility

Compatible with all Aspire Lifting Slings



Explore full lifter range and compatible slings

**AS/NZS ISO 10535:2024 COMPLIANT**

Code	Lifting Height	Under Bed Clearance	Turning Diameter	Base Length	External Base Width	Internal Base Width	Overall Height	Front & Rear Castors	Lifter Weight	SWL
LSS390652	692 - 1770mm	136mm	1245mm	1100mm	682 - 1138mm	519 - 986mm	1095 - 1880mm	100 & 125mm	39.6kg	250kg