

07/10/21

ASPIRE COVE CHAIR – NOW AVAILABLE!

The new Aspire Cove Chair (the Cove Chair) helps to improve postural, pressure care and comfort outcomes for users that require medium to high level support.

It is available in two sizes and offers adjustable seating to accurately fit a wider range of users.


The Cove Chair replaces the incumbent Shoalhaven Chair (CHP198950) and is far more competitive in features and price.

KEY FEATURES

- **Anatomical Legrest Movement:** As the legrest raises, it extends outwards tracking knee joint range of motion and maximises lower limb support. Great for taller users without raising overall seat height.
- **Performance and Aesthetics:** main surface is DARTEX maximising skin microclimate and pressure redistribution. External surfaces are medical grade coated fabrics from Warwick Australia. These hardy yet stylish patterns offer an appealing aesthetic and blend seamlessly in with modern Aged Care interior designs.
- **Adjustable Pressure Seat Cushion:** The integrated seat cushion combines layered memory foam with a segmented hybrid air-foam cell that offers convenient pressure adjustment to carers and maintenance staff. Re-inflation is as simple as opening and closing a valve.
- **Ease of Movement:** High performance TENTE Levina castors have been chosen to reduce push/pull forces during transportation. Comes with two directional locking (front) and two total locking castors.
- **Ergonomic Carer Push Handles** are accessible throughout tilt range reducing WH&S risks.
- **Transitional Support Armrests** trace the user's arm movement throughout tilt range to maximise comfort and confidence when being repositioned.
- **Stowable Utility Table** (optional accessory) that remains horizontal during tilt.
- **Weight-Bearing Footplate** (optional accessory) redistributes pressure more evenly in upright sitting and improves transfer efficiency by enabling users to be positioned further back into the seat upon entry.
- **Custom Seating Options** include the Advanced Immersion Backrest cushion and compatibility with after-market high care seating surfaces such as ROHO cushions.



New Product Launch

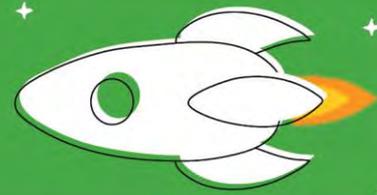

PRODUCT CODES

Product Code	Description
Chairs	
CHP198960	Chair Pressure Relief - Aspire COVE Chair - Grey/Andreas Pumice - Medium - 160kg
CHP198965	Chair Pressure Relief - Aspire COVE Chair - Navy/Bronx Steel - Medium - 160kg
CHP198970	Chair Pressure Relief - Aspire COVE Chair - Grey/Andreas Pumice - Large - 160kg
CHP198975	Chair Pressure Relief - Aspire COVE Chair - Navy/Bronx Steel - Large - 160kg
Accessories	
CHP198980	Stowable Utility Table - Aspire COVE Chair
CHP198985	Weight Bearing Footrest - Aspire COVE Chair
CHP198990	Advanced Immersion Backrest Cushion - Medium - Grey - Aspire COVE Chair
CHP198991	Advanced Immersion Backrest Cushion - Medium - Navy - Aspire COVE Chair
CHP198992	Advanced Immersion Backrest Cushion - Large - Grey - Aspire COVE Chair
CHP198993	Advanced Immersion Backrest Cushion - Large - Navy - Aspire COVE Chair

Pictured (left to right):



Weight-Bearing Footrest



Stowable Utility Table



Advanced Immersion Backrest Cushion

