

27/07/2021

## MALVESTIO DELTA 4 WARD BED – NOW AVAILABLE!

We are pleased to announce the arrival of our new Ward Beds from Malvestio.


Made in Italy with our production partners Malvestio, the Delta 4 Ward bed is a high-end bed solution that allows us to be very competitive within the Australian market, particularly in the Acute space.

Available in two versions, the Delta 4 is sure to tick all requirements for modern day healthcare:

- **Delta 4 Ward Bed – Standard**
- **Delta 4 Ward Bed W/DeltaDrive (electric drive wheel) & Weigh Scale**

### PRODUCT FEATURES

- o **INCENTRO System:** the combined back and pelvic translation nullifies horizontal forces, prevents slipping and keeps the body at a 30° angle, optimal for ventilation.
- o **Safe Exit Lights:** 2 x LED lights signal whether the bed is at a safe exit height, promoting safe egress for the patient (orange is unsafe and green is safe).
- o **Auto Regression:** in backrest automatically shifts backwards thus increasing the space on the pelvic area by creating an extra space of 13cm.
- o **Biocote:** metal components are epoxy powder coated and finished with an antibacterial additive for improved infection control.
- o **DeltaDrive:** electric centre 5<sup>th</sup> wheel with integrated controls in headboard provides increased manoeuvrability for hospital staff (*available on BEB050150 only*).
- o **Weigh Scale:** integrated system allows for precision weight measurements that can be set to automatically weigh a patient every 24 hours (*available on BEB050150 only*).

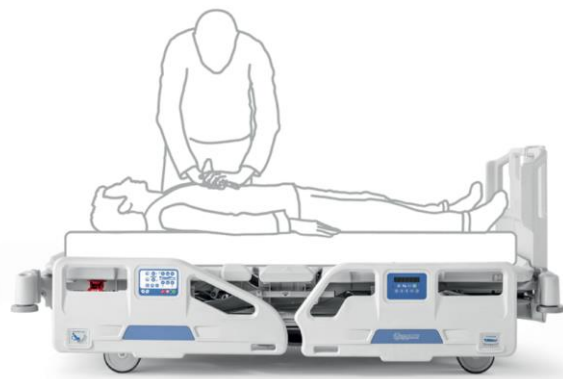

**New Product Launch**

# New Product Launch

## CONFIGURATIONS

---

- o Delta 4 Ward Beds are available in two configurations:
  - o BEB050140 – Delta 4 Ward Bed
  - o BEB050150 – Delta 4 Ward Bed W/DeltaDrive (electric drive wheel) & Weigh Scale



**BEB050140**  
**Delta 4 Ward Bed**

**BEB050150**  
**Delta 4 Ward Bed  
W/DeltaDrive &  
Weigh Scale**



***New Product Launch***