



# Technical Change Notice



09/11/2020

## BREEZY EVERY DAY WALKING FRAMES UPDATE!

A bit of a hidden feature – the Breezy Every Day Folding Walking Frames can actually be fixed in the open position using a simple C-clip which is included with each unit! When fitted correctly, this effectively converts the Folding model into a Non-Folding model, and is a common preference in some care environments.

Note: only one C-clip is needed per unit.

In addition, all Breezy Every Day Folding Walking Frames now come fully assembled from the factory, saving us and our distributors much assembly time. This was not possible for the Non-Folding Frames due to shipping space constraints.

We have therefore discontinued the Non-Folding Frames versions (S, M, L) and associated kits. See next page for code-specific details.

Please refer to the Assembly/Fitting Instructions at the bottom of this PAB for more info on fixing the frame with the C-clip. Or alternatively see pg. 5-6 or the Breezy Everyday User Guide attached.

### DESCRIPTION

---


- o The retaining C clip is shown below



Retaining C-Clip




Location of retaining clip

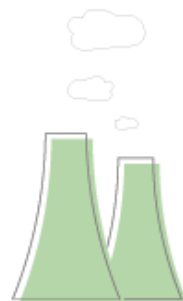

## Technical Change Notice

- o The following changes were made in GP:

Discontinued		Replaced by	
WAF704900	Walking Frame - Breezy Every Day - Alum Non Folding - Small	WAF706700	Walking Frame Assembled - Breezy Every Day - Alum Folding/Non Folding - Small
WAF704900K	Walking Frame KIT - w/ Front Castors and Rear Skis - BE - Alum Non Folding - Small	WAF706700K	Walking Frame KIT - w/ Front Castors and Rear Skis - BE - Alum Folding/Non Folding - Small
WAF704800	Walking Frame - Breezy Every Day - Alum Non Folding - Medium	WAF706800	Walking Frame Assembled - Breezy Every Day - Alum Folding/Non Folding - Medium
WAF704800K	Walking Frame KIT - w/ Front Castors and Rear Skis - BE - Alum Non Folding - Medium	WAF706800K	Walking Frame KIT - w/ Front Castors and Rear Skis - BE - Alum Folding/Non Folding - Medium
WAF704700	Walking Frame - Breezy Every Day - Alum Non Folding - Large	WAF706900	Walking Frame Assembled - Breezy Every Day - Alum Folding/Non Folding - Large
WAF704700K	Walking Frame KIT - w/ Front Castors and Rear Skis - BE - Alum Non Folding - Large	WAF706900K	Walking Frame KIT - w/ Front Castors and Rear Skis - BE - Alum Folding/Non Folding - Large

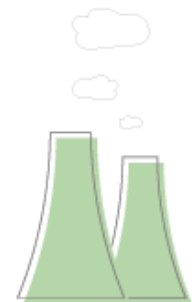

- o Extract from Breezy Every Day Aluminium Folding Walking Frames:

### ASSEMBLY / FITTING INSTRUCTIONS

#### TO OPEN


1. Holding the frame top in one hand and rail front in other, pull one away from the other.
2. Ensure the frame is in the fully open position before use (Fig 3).
3. Fit safety clip (Fig 6).

#### TO CLOSE

1. Remove safety clip.
2. Grasp the frame top in one hand and pull the wooden ball in the other hand (Fig 4).
3. When closed, the frame will fold into position (Fig 5).

As a safety feature to prevent the frame from closing during use, a safety clip must be in place.


The safety clip must be located into position when the frame is in the open position (Fig 6).

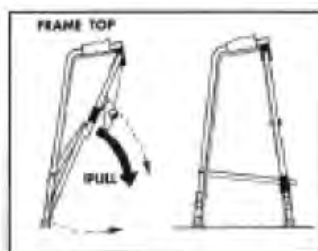

**FIG 3**



**FIG 4**

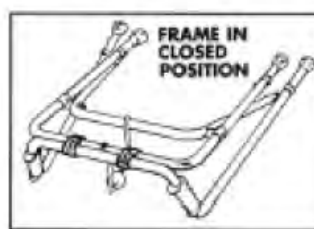

**FIG 5**



**FIG 6**

