

PREVENTATIVE MAINTENANCE

Aidacare's Preventative Maintenance and Asset Management Programs are tailored to meet our customer's requirements.

OUR COMMITMENT IS YOUR SATISFACTION

Aidacare is dedicated to the repair and on-going maintenance of healthcare equipment.

We offer a commitment to our customers to ensure we deliver the best service at competitive pricing.

Aidacare's maintenance includes advantages of:

FULLY TRAINED TECHNICIANS

All of our technicians are qualified and trained in our own manufacturing facilities, and by our suppliers. They are qualified to ensure your healthcare equipment is serviced according to the Australian standards, so it is compliant and safe to use.

ALL WORK GUARANTEED

Parts and labour are guaranteed. Using only quality parts and the latest diagnostic techniques, all repairs will be performed to the manufacturer's standards. We also offer the latest in product updates and refurbishments.

COMPETITIVE RATES

Our access to all major manufacturers' parts gives us the ability to offer the best pricing ensuring that our charges are competitive in every market every time.

MOBILE WORKSHOPS

For inspection, service and repair - we can come to you. Our service vans are fully equipped mobile workshops carrying a large range of spare parts. Our service technicians have quick access to service records and pricing details to minimise the downtime of your equipment.

PREVENTATIVE MAINTENANCE

Our service Agreements can cover a vast range of healthcare equipment and incorporate inspection and maintenance to ensure:

- · Risk Management
- Compliancy
- Maximising product lifetime
- Reassurance that the equipment is safe for patient care

Benefits of our Service Agreements:

- · Provides service history of the product
- · Compliancy and accreditation preparedness
- Asset Condition Profiling (Good / Poor / Replace etc.)
- Improved risk management via consistency of supply and control of service levels
- Improved financial planning for R&M budgets
- Improved CAPEX planning and insight into the condition of your assets via the Service Asset Register

Recommended Equipment Service Frequency

6 months

- · Mobile Lifting Hoist
- · Ceiling Hoist
- · Alternating Pressure Mattress
- Sanitisers
- Macerators
- Patient Slings

1 year

- · Hospital Bed (electric)
- Aged Care Bed (electric)
- · Medication Cart/Commercial Trolley
- · Pressure/Positioning Chair
- · Lift Recliner Chair
- · Wheelchair (manual)
- · Wheelchair (electric)
- · Wheelchair (tilt in space)
- Scooters (electric)
- · Shower Commode
- · Shower Commode
- Shower Trolley (manual)Shower Trolley (electric)
- Walking Frame
- Weigh Scales (not calibration)
- · Electrical Test and Tag (optional)

Maintenance Checklists



Lift Recline Chair

- √ Raising, lowering, heat and massage functions of chair
- ✓ Function, casing and cable of hand control integrity
- ✓ Integrity, condition and moving parts of frame structure
- ✓ Integrity, condition and tears in upholstery
- ✓ Movement and noise of actuators
- ✓ Electrical leads and connections



Pressure Position Chair

- ✓ Cable damage and connections of chair and tilt function
- ✓ For defects in cover such as wear and tear
- ✓ Movement and function of castors and brakes
- √ Frame for defects or damage
- ✓ Resistance from gas strut operation
- ✓ For leaks in bladders



Wheelchair

- ✓ Function, movement and folding frame of wheelchair operation
- ✓ Integrity, function and moving parts of wheelchair structure
- Movement, wear and pressure of tyres and castors
- ✓ Integrity, condition and movement of footplates
- ✓ Integrity and condition of upholstery
- ✓ Movement and adjust brakes



Electric Scooter

- ✓ Inflation and wear of tyres and anti-tip
- ✓ Function and wear of throttle and speed control
- ✓ Padding, movement and function of seat and seatbelt
- ✓ Function, bulb and wiring of lights
- ✓ Integrity and wear of frame, steering shroud and suspension
- ✓ Charge capacity, leads and connection of batteries and charger
- ✓ Function of electronics
- ✓ Function and/or adjust brakes



Shower Commodes

- ✓ Structural, corrosion, wear and tear of frame
- ✓ Movement, brakes and castors
- ✓ Integrity, condition and operation of armrests
- ✓ Padding and secured seat integrity
- ✓ Operation and integrity of footrests
- ✓ Operation, integrity, gas struts, cables of tilt and recline system



Shower Trolley

- ✓ Functions, movement, control system components of the shower trollev
- ✓ Structural, corrosion and movement of integrity frame
- ✓ Movement, wear, condition and clean castors
- ✓ Attachments of hoses and drainage
- ✓ Integrity, condition and tears of upholstery
- ✓ Movement and adjusting locking brakes

Maintenance Checklists



Hospital and Aged Care Beds

- ✓ Raising and lowering of bed operation
- ✓ Function, casing and cable
 of hand control integrity
- ✓ Movement, brakes and castors
- ✓ Integrity and condition of frame structure
- ✓ Movement and noise of actuators
- ✓ Electrical leads and connections
- ✓ Integrity and condition of mattress structure
- Structural and movement of bed ends




Mobile Lifter and Ceiling Hoist

- ✓ Raising and lowering of hoist operation functions
- ✓ Function, casing and cable of hand control integrity
- ✓ Movement, brakes and castors
- ✓ Integrity, condition and moving parts of frame structure
- ✓ Movement and load (weight) test of actuators
- ✓ Electrical leads and connections
- ✓ Test leads and connections of batteries and charger



Sanitiser Macerator

- ✓ Test run, controller and functions
- ✓ Solenoids, steam and pumps of control system
- ✓ Noise, function and bearings of wash pumps
- ✓ AO 600 steam verification test for thermal integrity
- ✓ Operation and wall charts in place
- ✓ Flush, wash, steam and cool of cycle times
- Auto, open and close of lid functions
- ✓ Leads, connections and switches of electrical items
- ✓ Pump, flow and close of detergent system

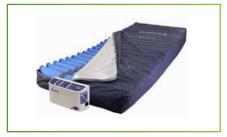

Alternating Pressure Mattress

- ✓ Airflow and alternating cycle of operation functions
- ✓ Function, casing and cable of alternating pump integrity
- ✓ Bellows, hoses and connections of mattress
- ✓ Filters, pump and moving parts of alternating pump
- ✓ Electrical leads, connections and electrical test and tag
- ✓ Material and zippers on mattress cover



Medication Carts & Trolleys

- ✓ Movement and functions of commercial trolley operation
- ✓ Integrity, condition and moving parts of frame and draw structure
- ✓ Integrity, condition and movement of castors
- ✓ Movement and adjusting locking brakes

Find your local service centre

With nationwide coverage, Aidacare is always available to provide quality service, ensuring your equipment is safe, compliant and maintained to the highest standard for best results and peace of mind.

Contact us today to book with your local service centre