

# ECO ASSIST

EMPOWER YOUR MOBILITY



REHASENSE




**2** YEAR  
WARRANTY

Exclusively distributed by  
**aidacare™**  
HEALTHCARE EQUIPMENT


# ECO ASSIST

## LIGHTWEIGHT POWER ADD-ON

ECO ASSIST is a powered add-on that transforms a manual wheelchair into a motorised mobility aid. Made from aluminium and steel, its lightweight and durable design is easy for users to dock to their wheelchair in seconds. Featuring a robust motor powered by an efficient lithium-ion battery, it's ideal for navigating city streets with ease. With adjustable controls and reliable breaking, ECO ASSIST allows users to enjoy everyday journeys with peace of mind.

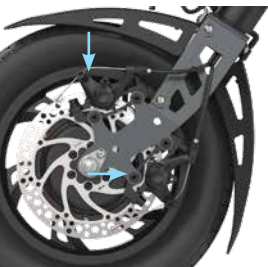

**Dual Brake System**  
Mechanical disc and electric braking system

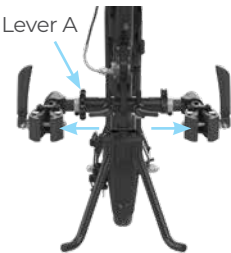

**Width Adjustable**  
Quick adjustment via manual clamps



**Leg Supports**  
Height adjustable between 260-420mm



**Quick Disconnect**  
Supports easier transportation and storage with a single point disconnection

Easy manual docking system allows for quick and easy mounting  
Compatible with both rigid and folding wheelchairs, with fixed front frames or removable legrests  
Universal clamp system allows fitting to most wheelchair types\*  
160mm disc brake for strong, reliable braking

Fits wheelchairs with a width of 210 – 480mm (between footplates)  
Clamp height adjustable from 260 – 420mm  
Walk-assist mode  
Traction control for improved grip  
Reverse gear for easier manoeuvring

### Accessories



**Additional Footrest**  
WCA872620



**Bag**  
WCA872268



Eco Assist	Wheel Size	Engine Power	Max. Engine Torque	Speed Modes	Max Speed	Driving Range	Overall Weight	Weight of Heaviest Part	SWL
WCA872000	12"	350W	39Nm	3	10km/h	Up to 30km	15.5kg	10kg	110kg

\*Not to be used with Carbon Fibre wheelchair frames; warranty will be voided.