

BT AND BT PLUS MEDICAL SIDE TABLES.

Exclusively Distributed by:








MEDICAL SIDE TABLE BT AND BT PLUS FOR HOSPITAL USE.

They roll in softly: the new GREINER BT and BT plus medical side tables. Developed and designed as top-quality additions to the MULTILINE NEXT DC and IT treatment chairs, these side tables ensure the focus is fi rmly on patient comfort. The BT with fixed tray and the BT plus with swivelling tray are flexible medical side tables that allow you to ensure your patients are as comfortable as possible when undergoing dialysis, chemotherapy or other treatments. They provide a unique combination of stability and lightness and can be adjusted with one hand. Everything you need to care for your patient is right at your fingertips.

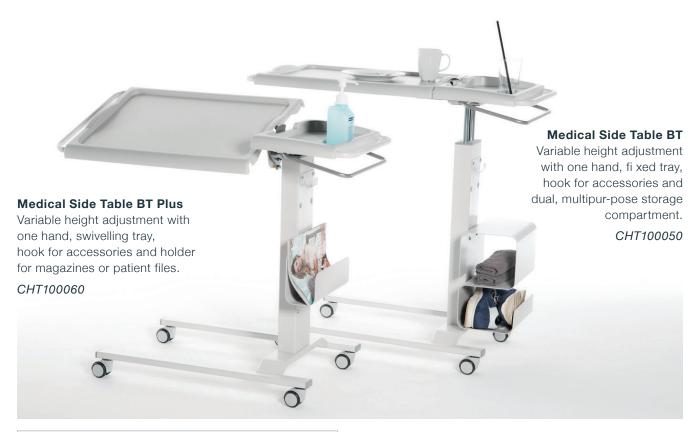



Perfect for combining with the MULTILINE NEXT DC and IT (IT shown here)



Height adjustment with dual control



Design and quality made by GREINER



Swivelling tray (BT plus)



Multipurpose accessories, see

The epitome of functionality, ergonomics, stability and hygiene, the new overchair tables have been meticulously designed to ensure they meet the most stringent requirements of patients and nurses alike. The sturdy frame in typical GREINER quality is combined with smooth-running castors to ensure the table is easy to move while still safely holding loads of up to 30 kilos. Powerful pneumatic springs make it easy to gently adjust the height of the tray with one hand – and provides added safety if the patient is in the emergency Trendelenburg position. On the BT plus the tray can be tilted at the touch of a button so that it is ergonomi-cally positioned for both patient and nurse. The removable trays have recessed grips and are the perfect size for hold-ing meal trays, patient items, clinical accessories and disin-fectants. The tough, resistant surfaces are quick and easy to clean and disinfect. A range of accessories are also avail-able to improve the patient experience still further.

□ Medical Side Table (LxW)	970 x 460 mm
□ Tray Surface (LxW)	560 x 390 mm
□ Side Tray (LxW)	370 x 210 mm
□ Variable Height Adjustment	250 mm (790 - 1040 mm)
□ Diameter of Castors	50 mm
□ Tilt Angle of Tray (BT Plus)	200 / 400 on each side
Maximum Load	30 kg
□ Weight	Approx. 22.5 kg
Colour	RAL 7047
Accessories	
Hook	For bags, etc.
Magazine Holder	For magazines, patient files, etc
Dual, multipurpose storage compartment	For shoes, wipes, books, etc,

GERMAN DESIGN AND QUALITY

Made by GREINER means made in Germany. GREINER has been supplying products for the medical field since the 1970s. A consistently high standard is applied in the development and production of GREINER couches. The same applies to our range of services. We achieve this by translating our flair for innovation into fully automated processes and handmade quality. High standards in design, materials, engineering and ergonomics combined with client-oriented flexibility make GREINER a partner to quality-conscious customers at the highest level.



we design quality.

Agreeably health-promoting seating has motivated us since 1922. For decades we have been making high-quality, technically sophisticated equipment within the medbest range. Our innovative chairs for examination, dialysis, oncology, blood sampling, care and transport are all made in Germany. Our quality commitment is verified by an independently audited quality assurance system: GREINER is certified under DIN EN ISO 9001 and 13485.



Test mark for quality management system

Exclusively Distributed by:





Aidacare P/L

3A, 1 Moorebank Avenue Moorebank NSW 2170 T – 1300 133 120

W - www.aidacare.com.au

