

MULTILINE NEXT AC. MOBILE PATIENT COUCH FOR OUTPATIENT SURGERY.

Exclusively Distributed by:







GREINER MEDBEST: MEDICAL DEVICES MADE IN GERMANY.

In response to the wishes of patients and those treating them, GREINER has taken up the challenge of developing inno-vative solutions that address clinical and economic needs in equal part. Made by GREINER means made in Germany. A consistently high standard is applied in the development and production of GREINER couches. The same applies to our range of services. We achieve this by translating our flair for innovation into fully automated processes and handmade quality. High standards in design, materials, engineering and ergonomics combined with client-oriented flexibility make GREINER a partner to quality-conscious customers at the highest level.

NEW COUCH CONCEPT FOR OUTPATIENT SURGERY.

Outpatient surgery centres depend on flexibility, precision of timing and handling, mobility and safety; the efficient use of space and cost-effectiveness are also critical. In short, patients should be fully involved in processes before and after an operation. In the modern clinical environment, cumbersome hospital beds make life needlessly difficult for healthcare teams and patients alike. Hospitals need couches that can admit prepared patients for operations with maximum flexibility and sta-bility, transport them safely and comfortably to the pre-op waiting room, enable flawless repositioning and safely and comfortably move patients to the recovery room after an operation. This is precisely what the new MULTILINE NEXT AC delivers.



MULTILINE NEXT AC MOBILE COUCH.

This new couch meets the patient positioning and mobilisation challenges posed by the growing number of outpatient operations. It is ideal for surgical preparation, transporting patients to and from the operating theatre, positioning after surgery and deployment in the recovery room after anaesthesia. Ergonomics, functionality, mobility and patient comfort form a unique combination of chair and couch. Complete ease of access and simple disinfection ensure maximum hygiene. The MULTILINE NEXT AC can also be installed in small spaces where conventional beds will not fit.

GREINER for Outpatient Surgery:


- Areas of application: Surgical preparation, post-operation rehabilitation and recovery couch, examination, treatment, mobilisation, transport
- Totally new concept in terms of design, technology and features
- Electric motor adjustment of backrest, seat and leg section as well as height (adjustment up to 900mm)
- Rechargeable battery, handles, castors and frame for mobile operation
- Emergency foot switch for Trendelenburg position, ensuring staff keep their hands free
- Independently certified product safety and superb reliability



MULTILINE NEXT AC MOBILE PATIENT COUCH.

Ergonomics, stability, mobility and safety combine to create the perfect solution for outpatient surgery. Access to the chair coupled with comfortable cushioning ensure patients and carers remain calm and composed at all times, even in stressful and critical situations. The powerful battery maintains mobility and safety for many hours, while adjustment to the appropriate table height and patient positioning to operation height ensure proper ergonomics. The two fully retractable armrests, which can easily be operated with one hand, ensure safe resting and transporting as well as ideal access when repositioning.



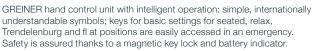

Comfort, functionality and ergonomics.

With stepless motorised adjustment and intelligent and intuitive operation via the selection and control keys of the GREINER control element, the MULTILINE NEXT AC mobile patient couch delivers outstanding seated and reclining comfort in all situations. From seated and relax positions to the fl at bed position and emer-gency Trendelenburg position, the couch quickly and smoothly moves to any required or preferred position at the touch of a button.

The GREINER design concept makes it easy for carers and patients alike to use the four motors for the backrest, seat and leg sections and height adjustment (up to 900mm) efficiently. The couch is mobile, stable and secure.



Seated position



Relax position



Bed/flat position



Trendelenburg and shock position



MULTILINE NEXT AC. NEW LEVEL OF COMFORT FOR OUTPATIENT SURGERY.

GREINER gives you the opportunity to individualise treatment couches according to your precise needs. For example, you can choose your preferred upholstery colour from the broad colour chart. Covering materials, which have a biocompatible surface, are resistant to urine and blood, scratch-proof, heat resistant and impervious to disinfectants and soiling. Equipment and accessories are designed with attention to detail.



Patient-friendly access at heights of 580-900 mm



Ease of access to powerful rechargeable battery



Handles for seated/reclined position



Central locking/directional stability



125 mm castors, perfect mobility



Pivoted infusion pole



Fully retractable armrests



Emergency foot switch



Standard rail



Wall charging point



HygieneCompliance with stringent regulations, ease of access, clean surface finish, tested surfaces, sparing use of seams

Ergonomics

Flush upholstery plus fixed and retractable armrests for simple repositioning before and after surgery

Mobility 125 mm castors for comfortable transporting, central braking system

straight line

MULTILINE NEXT AC: Equipment and accessories		
Product Codes (treatment chair & inclusions - listed below)		
CHT100350 - Multiline NEXT AC - Navy Blue Vinyl		
Equipment		
Head rest, sliding and removable		
Central braking system, one wheel with directional stability		
Foot switch/emergency switch		
Battery pack, includes storage tray and charging lead		
Head cushion		
Emergency foot switch		
Accessories		
Standard rails	Replacement battery	
Handle, rubber-coated	LED reading light	
Infusion pole	Infusion pole holder, pivoted, one-handed operation	
Wall charging point		

Dimensions and technical data	
Chair upright (floor space)	1070 x 770 mm
Recline position (floor space)	2000 x 770 mm
Total width	770 mm
Surface in reclined position (L x W)	2000 x 600 mm
Thickness of padding	90 mm
Armrest (L x W)	730 x 70 mm
Access height	580 mm
Height adjustment	320 mm (580 - 900 mm)
Castor diameter	125 mm (1 directional)
Backrest	- 5° to 78°
Seat section	0° to 22°
Leg section	0° to 90°
Maximum patient weight	220 kg
Weight of couches	128 kg
Voltage	100-240 V/50-60 Hz
Battery peformance	24 V DC, 1 x 2.9 Ah (up to 40 movements)





GERMAN DESIGN AND QUALITY

Made by GREINER means made in Germany. GREINER has been supplying products for the medical field since the 1970s. A consistently high standard is applied in the development and production of GREINER couches. The same applies to our range of services. We achieve this by translating our flair for innovation into fully automated processes and handmade quality. High standards in design, materials, engineering and ergonomics combined with client-oriented flexibility make GREINER a partner to quality-conscious customers at the highest level.



we design quality.

Agreeably health-promoting seating has motivated us since 1922. For decades we have been making high-quality, technically sophisticated equipment within the medbest range. Our innovative chairs for examination, dialysis, oncology, blood sampling, care and transport are all made in Germany. Our quality commitment is verified by an independently audited quality assurance system: GREINER is certified under DIN EN ISO 9001 and 13485.



Test mark for quality management system



Test mark for product certification

Exclusively Distributed by:





Aidacare P/L

3A, 1 Moorebank Avenue Moorebank NSW 2170 T – 1300 133 120 W – www.aidacare.com.au

