

MULTILINE NEXT IT

THE MULTIFUNCTION TREATMENT CHAIR FOR WARD, THEATRE AND CONSULTING SUITES



220kg

Safe working load



Exclusively Distributed by:



MULTILINE NEXT IT



MULTILINE NEXT IT.

This compact couch meets the specific needs of oncology day hospitals, infusion medicine and transfusion medicine. It combines space-saving design with the comfort and security of fully equipped models. Patients and those caring for them enjoy the full range of amenities and features, with the couch quickly conveying patients into the most suitable position for therapy. With a seat width of 60 cm and the perpendicular angle of the leg section, it can be used as a chair. The comfortable, extra soft upholstery, patient-friendly armrests that support therapy and wide range of fittings make the MULTILINE NEXT IT a next-generation couch.



COMPACT THERAPY COUCH FOR INFUSION AND TRANSFUSION.

Chemotherapy in outpatient oncology departments is a major issue of our times. With much of the treatment being administered while patients are seated, there is a demand for couches that offer maximum comfort in this position too. With the available space often confined, couches are needed that can demonstrate quality even in very enclosed areas. Upholstery and armrests are critical to patient comfort. Depending on age and health status, patient access in the seated position also plays a key role in acceptance by patients and those caring for them. Generally speaking, the demands are high – and the new MULTILINE NEXT IT meets them in full thanks to the chair-like access it offers.

GREINER for oncology:

- ◆ Totally new concept in terms of design, technology and features; smooth-running mobility with central locking if required
- ◆ Version with three motors for the stepless adjustment of backrest, seat and leg section into required positions, including the emergency Trendelenburg position
- ◆ Comfortable cushioning in a wide range of colours, high levels of seated and reclining comfort, maximum hygiene with a minimum of seams
- ◆ Unique design with 90° angle of leg section – perfect for chair-like access and exit, minimal space requirement
- ◆ Flexible design variants and configurations adapt to every practical purpose
- ◆ Independently certified product safety and superb reliability



MULTILINE NEXT IT. ONCOLOGY THERAPY COUCH.

The new therapy couch developed by GREINER packs the performance of technically sophisticated couches into a compact form for restricted spaces. Designed for patients who may be very weak or afflicted by extreme symptoms, the couch helps them through the difficult times. It features ease of access, smooth rollers, clear functions, stepless motorised movement to the required position, comfortable armrests and cosy upholstery in soothing colours. The emphasis is on patient comfort, ergonomics and efficiency to ensure the best treatment conditions for patients over hours of chemotherapy.



GREINER hand control unit with intelligent operation: simple, internationally understandable symbols; keys for basic settings for seated, relax, flat bed and Trendelenburg positions are easily accessed in an emergency. Safety is assured thanks to a magnetic key lock.

Comfort, functionality and ergonomics.

With stepless motorised adjustment and intelligent and intuitive operation via the selection and control keys of the GREINER control element, the MULTILINE NEXT IT for infusion and transfusion therapy promises high levels of seated and reclining comfort in all situations. The three-motor couch quickly and gently moves to any required or preferred position – from seated, relax and bed positions to the emergency Trendelenburg position – at the touch of a button. Access and exit is considerably eased by the vertical leg section and access height of 54 cm.

The GREINER design concept makes it easy for carers and patients alike to use the powerful motors for the backrest, seat and leg sections efficiently. The couch is mobile, stable and secure.



Seated position



Relax position



Bed/flat position



Trendelenburg and shock position



MULTILINE NEXT IT. NEW LEVEL OF COMFORT FOR INFUSION AND TRANSFUSION.

GREINER treatment couches are supplied with many inclusions as standard, however can be tailored to your precise needs. Premium covering materials and upholstery offer extra-soft, wear resistant surfaces; a dynamically flexible suspension system within the seating surface supports freedom of movement for the spine, pelvis and hips, thus relieving pressure to the benefit of patients. Innovative MULTILINE five-position armrests, the first to move in five directions, are setting new standards.



Low-level access from the front and side



Unique armrest technology with patient-friendly positioning



Head rest as standard



Fold-up armrests



LED reading light (available on request)




Handle with paper roll holder



Available in Navy Vinyl.
A range of custom colours are available upon request. (Lead times may vary)



Emergency foot switch

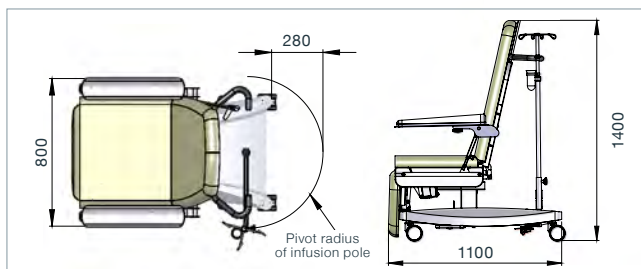

Flawless technology even in very enclosed spaces (measurements in mm)



Product safety
GREINER treatment
couches are independently
tested and certified

Hygiene
Design compliant with stringent regulations,
ease of access, clean surface finish, tested
surfaces, minimum of seams

Flexibility
Infusion pole holder can be pivoted
left and right with just one hand

Stability of movement
Fixing of backrest in any position
makes couch easier to push

Ergonomics
Unique armrest technology with
patient-friendly positioning - ideal
for infusion therapy

Comfort
Comfortable upholstery
with wide 60cm hand made seat



Access
Ease of access from front
or side with folded-up armrest,
access height of 54 cm

Seated position
Genuine chair position with
leg section at 90° angle

Mobility and braking
10 cm twin-wheel castors with
central locking if required

MULTILINE NEXT IT: Equipment and accessories

Product Codes (treatment chair & inclusions - listed below)

- CHT100110 - NEXT IT Treatment Chair - Navy Blue Vinyl

Inclusions

- Battery backup
- Handle with paper roll holder
- Padded armrests
- IV pole holder
- Head cushion
- Emergency foot switch

Additions

- Premium vinyl
- LED reading light
- Seat with pressure relieving springs
- Protective legrest cover

Dimensions and technical data

| | |
|--|--------------------------------|
| • Chair upright (floor space) | 110 x 78 cm |
| • Recline position (floor space) | 196 x 78 cm |
| • Total width | 78 cm |
| • Surface in reclined position (L x W) | 196 x 60 cm |
| • Thickness of padding | 7.5 cm |
| • Armrest (L x W) | 56 x 12 cm |
| • Access height | 54 cm |
| • Castor diameter | 10 cm |
| • Backrest | - 5° to 78° |
| • Seat section | 0° to 22° |
| • Leg section | 0° to 90° |
| • Maximum patient weight | 220 kg |
| • Voltage | 230 V/50-60 Hz, 110 V/50-60 Hz |
| • Weight of couches | 78 kg |



GERMAN DESIGN AND QUALITY

Made by GREINER means made in Germany. GREINER has been supplying products for the medical field since the 1970s. A consistently high standard is applied in the development and production of GREINER couches. The same applies to our range of services. We achieve this by translating our flair for innovation into fully automated processes and handmade quality. High standards in design, materials, engineering and ergonomics combined with client-oriented flexibility make GREINER a partner to quality-conscious customers at the highest level.



we design quality.

Our quality commitment is verified by an independently audited quality assurance system: GREINER is certified under DIN EN ISO 9001 and 13485.



Test mark for
quality management system



Test mark for
product certification

Exclusively Distributed by:



Aidacare P/L

3A, 1 Moorebank Avenue
Moorebank NSW 2170
T – 1300 133 120
W – www.aidacare.com.au

