

catch them before they fall

CENTRAL PROCESSOR (optional)


The Central Processor collects and stores the following information:

- Ward
- Bed
- Patient
- Alarm Time
- Alarm Source Patient Or Nurse

The information compiled can be used as a management tool to form a patient activity profile of peak times when a fall may occur.



INVISA-BEAM® Dual Beam Bed Monitor

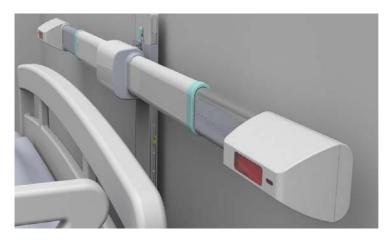

AIDS FOR FALLS PREVENTION

SUITABLE FOR

- Hospitals
- Aged Care Facilities
- Nursing Homes
- High Care

Wall Plate, CPU, Range Extender and powered remote also available.

Exclusively Distributed by:



: 1300 133 120

N: www.aidacare.com.au

INVISA-BEAM® products are made in Australia.

INVISA-BEAM® Dual Beam Bed Monitor

FEATURES

- Aids for falls prevention from beds.
- Alerts the nurse at the earliest stage, when the patient attempts to leave the bed
- Patient security—can help prevent unauthorised entry to the bed
- Alerts the staff of any unauthorised entry to the hed
- Monitors both sides of the bed with precise invisible beams
- Total wireless operation
- Bed free of equipment and cables
- No switches
- Hands-free—a bonus for infection control

- Monitor alarm cancels when the nurse enters the field of the bed
- Monitor alarm automatically resets when the nurse leaves the bed
- Activity at the bedside recorded— a management tool
- Can be used for lowered beds; also for bariatrics.

Code	
BEA007216	Dual Beam Monitor w/ Mono & Stereo Plug
BEA007217	Dual Beam Monitor w/ Remote Alarm
BEA007218	Dual Beam Monitor w/ Invisalert

FNTRY PASS

The Entry Pass is like a master key. It allows the nurse free access to the patient without worrying about switches.

When the nurse attends to the patient in response to an alarm:

- The Entry Pass cancels the alarm
- The bed monitor alarm automatically resets when the nurse leaves the bed



The Entry Pass may be worn with retractable clip supplied.

BEA007211	Entry pass
BEA007212	5 Card Charger
BEA007213	10 Card Charger

INVISALERT

FEATURES

- Displays the location of the alarm source
 Bed & Room
- Vibrating alarm
- Audible alarm

INVISALERT allows the nurse free access to the patient without worrying about switches.



When the nurse attends to the patient in response to an alarm:

- INVISALERT cancels the
- The bed monitor automatically resets when the nurse leaves the bed.

BEA007214	Invisalert Pager
BEA007255	3.7V Rechargable Battery Pack
BEA007256	Power Supply Charger

Exclusively Distributed by:

