

Botero

BARIATRIC BED



Exclusively Distributed by:








The **TRANSLATION** of the backrest **ENHANCES BREATHING** by reducing pressure on the plexus area.

Backrest TRANSLATION also helps to redistribute pressure more evenly for improved comfort and skin integrity.

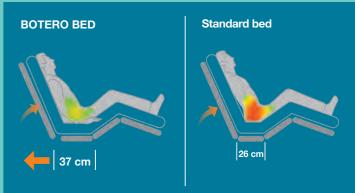

Scope for Changing Needs



Adjustable Mattress Platform

The mattress platform offers tool free adjustment from 102cm to 122cm in width.





Bed Extension

The mattress platform can be extended by up to 18cm with a padded mattress bolster available as an optional accessory.



Features

Low height of 38cm

4 x section electrically operated bariatric ward bed

Adjustable mattress platform with tool free adjustment from 102-122cm

Bed ends can be adjusted in width to suit mattress platform adjustment

Legrest adjustment is ergonomically synchronised with the knee break

Integrated bed length extension up to 18cm

IV & Self-Help Pole mounting brackets

Electric height adjustment with 2 x heavy duty actuators for a Safe Working Load (SWL) up to 400kg

Attendant control panel with intuitive button placement for bed control

Bed is equipped with battery backup as standard, allowing for power supply when transporting around wards or in the event of a power failure/outage

Translating backrest (Auto-Regression) shifts backwards by 11cm to compensate for patient movement towards the foot end of the bed as the backrest is raised.

Folding side rails as standard

Base frame is equipped with 150mm twin wheel castors featuring a centralised locking system that can be operated by foot pedals on either side of the bed end

Epoxy powder coating with antibacterial finish for improved infection control



Epoxy powder coating with infection control





226 - 244cm

102 - 122cm

OVERALL



38 - 85cm

BACKREST

KNEE BEND



UNDERBED HOIST







Technical data

Overall dimensions with side rails on	226 x 113 cm
Mattress platform dimensions at minimum extension	200 x 102 cm
Mattress platform dimensions at maximum extension	218 x 122 cm
Mattress platform height	min. 38 cm, max 85 cm
Trendelenburg/Reverse-trend.	+17° / -17°
Safe working load	400 kg
Degree of protection	IPX5

Medical device, in class I, in compliance with the EU Directive 93/42/EEC, modified by Directive 2007/47/EC

BEB050200	Botero Bariatric Ward Bed
BEM050200	Static Mattress w/Extension Bolster
BEA825100	Self-Help Pole
BEA825101	Telescopic IV Pole









Malvestio2018

The Company reserves the right to introduce all the technical and/or aesthetic modifications deemed necessary to their products without prior notice.











Exclusively Distributed by:



Aidacare P/L

4 Noonan Rd, Ingleburn NSW 2565 T - 1300 133 120

W - www.aidacare.com.au