

EDGE 3 shown in Blue Jean with iLevel® power seat elevation and Q-Logic 3 with Bluetooth®

## AVAILABLE WITH 30.5 CM (12") iLEVEL® POWER SEAT ELEVATION AT 7.24 KM/H

(4.5 MPH)

## Edge 3 Includes..

- USB charger
- Matching joystick shroud
- Fender lights
- Matching wheel inserts



quantumrehab.com

## Specifications

| •   |   |
|---|---|
| DRIVE WHEELS  | Standard: 35.6 cm (14") Solid<br>Optional: 35.6 cm (14") Pneumatic        |
| CASTOR WHEELS   | Front: 15.2 cm (6") Solid<br>Rear: 15.2 cm (6") Solid                     |
| ANTI-TIP WHEELS   | N/A   |
| MAXIMUM SPEED <sup>1</sup>  | Up to 9.6 km/h (6 mph)<br>Up to 7.24 km/h (4.5 mph) w/iLevel<br>activated |
| GROUND CLEARANCE  | 6.9 cm (2.7") (battery tray)  |
| TURNING RADIUS <sup>2</sup>   | 52.1 cm (20.5")   |
| OVERALL LENGTH <sup>2</sup>   | 82 cm (32.23") without front riggings                                     |
| BASE WIDTH  | 62.2 cm (24.4")   |
| SEATING SIZES   |   |
| Synergy®/TRU-Balance® 3 Static Seat                                       | Width: 25.4 -55.9 cm (10-22")<br>Depth: 25.4 -55.9 cm (10-22")            |
| TRU-Balance® 3 Power Tilt   | Width: 30.5 - 55.9 cm (12-22")<br>Depth: 30.5 - 55.9 cm (12-22")          |
| TRU-Balance® 3 Tilt and iLevel®   | Width: 30.5 - 55.9 cm (12-22")<br>Depth: 30.5 - 55.9 cm (12-22")          |
| SEAT-TO-FLOOR HEIGHTS <sup>10</sup>                                       |   |
| Synergy <sup>®</sup> /TRU-Balance <sup>®</sup> 3 Static Seat <sup>5</sup> | 41.9 - 47 cm (16.5-18.5")   |
| TRU-Balance® 3 Power Tilt <sup>5</sup>                                    | 44.5 - 49.5 cm (17.5-19.5")   |
| TRU-Balance® 3 Tilt & iLevel® <sup>5</sup>                                | 44.5 - 49.5 cm (17.5-19.5")   |
| MANUAL RECLINE <sup>6</sup>   | Yes   |
| BATTERY SIZE <sup>3</sup>   | NF-22 (2)   |
| BATTERY WEIGHT 7  | 17.24 kg (38 lbs.) each   |
| AVAILABLE ELECTRONICS 8   | 120A Q-Logic 3 EX   |
|   |   |
| BATTERY CHARGER   | 8A Off-board  |
| MOTOR PACKAGES  | 2-Pole, Mid-Wheel 6®  |
| WEIGHT CAPACITY <sup>4</sup>  | 136 kg (300 lbs.)   |
| BASE WEIGHT   | 63 kg (138.6 lbs.)  |
|   |   |


Features

Upgraded SRS (Smooth Ride Suspension)

- Gas-charged monotube shocks with springs mounted in a coil-over assembly for a smooth, comfortable ride and advanced stability
- The comfort of soft springs with the control of real dampening shocks
- All shocks function independently of each other
- Same springs and shocks as the 4Front™
- Close simulation to automative suspension
- 7.24 km/h (4.5 mph) at iLevel<sup>®</sup>
- Fender lights and USB charger
- Bold new colours
- Made in Australia of foreign and domestic parts



## 7 New Colours\* (12 in Total)

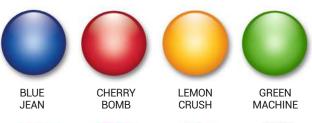



RASPBERRY BERET

**EVOLUTION** ORANGE

BACK IN BLACK



MIDNIGHT BLUE

NIGHT TRAIN

DFFP PURPLE

SUGAR

PLUM

GRFY STREET

3) Battery capacity may change when adding power positioning systems. 4) Weight capacity is dictated by the selection of seating.

Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure.
 Without front riggings

- 5) All measurements +/-3 mm (1/8") with 46 x 46 cm (18" x 18") seat, fixed towers, 136 kg (300 lbs.) weight capacity and solid tires (add 10 mm (3/8") with pneumatic tires). Widths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in the lowest position; does not include elevating system unless specified. 6) Maximum dimensions of the Manual Recline vary by power base and seating system. Call for details. 7) Battery weight may vary +/- 1.4 kg (3 lbs.) based on manufacturer.

- 8) All Quantum Rehab<sup>®</sup> units can be ordered with specialty controls.
  9) Call for availability
- 10) 41.9 cm (16.5") seat-to-floor height is not available with power tilt or iLevel®. 49.5 cm (19.5") seat-tofloor height is not available with static seat or power recline



Distributed by:

AFY146 240419

