

MANUAL WHEELCHAIR REFERENCE GUIDE

HALTI

AUSTRALIA



<mark>YOU DRIVE</mark> OUR DESIGN

QUICKIE wheelchairs are designed for you and by you. In development, we consider how you interact with your wheelchair to maintain maximum efficiency, performance, fit, and style.

ADJUSTABILITY AND VERSATILITY

Once you choose a QUICKIE, our goal is that it continues to be the right chair for you. We offer maximum adjustability in our designs, interchangeable parts, and the True Fit Growth Program* to help accommodate changes in your condition, environment, or lifestyle.

*Available on select QUICKIE wheelchairs

<mark>leading</mark> TECHNOLOGY


We are continually developing new products, improving our flagship products, and offering new options and solutions to enhance your lifestyle.

QUALITY AND DURABILITY

QUICKIE wheelchairs undergo rigorous testing that surpasses industry standards, so you can be confident that your product will hold up and perform throughout its lifetime.

BUILT-4-ME PERSONALISED

Our Built-4-Me Team is dedicated to creating product user-matched beyond what is listed on our product order forms. Experts in the engineering and manufacturing of configured rehab products will listen to your specific requests and build a chair as individual as your fingerprint!

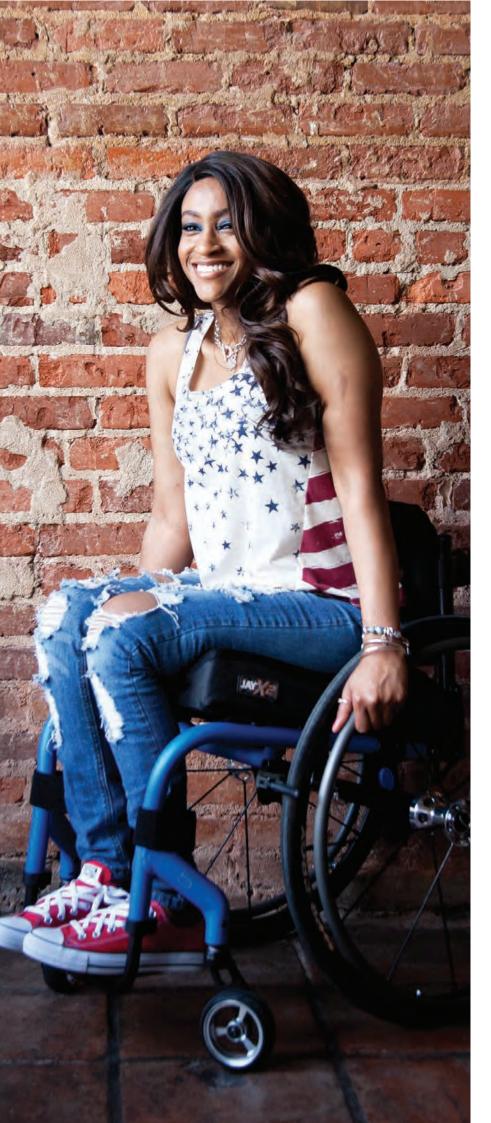



RIGID FRAME

Krypton R	5
7 ^R Series	6
5 ^R	7
GP Series	8
QRi	9
RIGID Comparison Guide	10

FOLDING FRAME

Krypton F	13
Xenon ²	14
Quickie 2	15
ΩXi	16
M6	17
FOLDING Comparison Guide	18
TILT IN SPACE	
Iris	20
SR45	21
TILT IN SPACE Comparison Guide	22

JAY [®] Seating	Selection	24





4

KRYPTON R CARBON FIBRE FRAME

REVOLUTIONARY CARBON FRAME WITH STAGGERING STRENGTH-TO-WEIGHT RATIO.

RIGID FRAME Ultra-lightweight from 3.65kgs in transit



CARBON BRAID TECHNOLOGY



This QUICKIE Krypton Range features a patented and fully automated BRAID process that makes it possible to build the perfect frame – strong, light and manufactured with a precise and reliable automated process.

RIGIDITY & RAZOR-SHARP AGILITY



It's not only the exquisite material or the famously seductive styling, but also the athleticism and agility of a unique rigid wheelchair. Infinitely adjustable centre of gravity, fine-tuning of the front and rear seat heights and backrest angle adjustment ensures the Krypton R is designed for something truly special.



7^R SERIES RIGID ACTIVE FRAME

STRENGTH IN MOTION DESIGNED TO MOVE WITH YOU!

Available in: 7^R - Adjustable 7^{RS} - Fixed



FREESTYLE BACKREST



The 7 Series' minimal, ultra-modern design takes backstage so YOU are in the limelight.

Designed to break barriers to independence by providing a full range of upper body motion and lowering the overall weight of your chair.

TOP OF THE LINE TECHNOLOGY



The genius behind the 7R Series' ultra lightweight frame is the 7000 Series aerospace aluminum that has been heat treated using ShapeLoc Technology.

The Result: thin, light tubing with strength & stiffness with a transport weight starting at 4.5kgs*

*Transport weight for the 7RS



Live without limits.

6

QUICKIE 5^R HIGHLY ADJUSTABLE FRAME

HIGHLY ADJUSTABLE TO ELEVATE YOUR MOBILITY

Frame Styles: CLASSIC

Traditional frame length for maximum stability **PERFORMANCE** Shortened length

for maximum responsiveness



AXLE PLATE DESIGN



The 5R's swept-back axle plate design isn't just for looks.

It is also incredibly quick to adjust and is available in two lengths to ensure you are getting the adjustments you need without adding extra weight.

XTR SUSPENSION



The XTR Suspension Frame helps you fearlessly conquer curbs, bumps, and tough terrain. Add Mountain wheels, Knobby tires, wide casters, and Frog Legs suspension caster forks for the ultimate off-road wheelchair experience.



7

GP SERIES RIGID WHEELCHAIR

A CLASSIC QUICKIE RIGID WHEELCHAIR

Available in:

GP - Standard Style

GPV - Tapered for increased maneuverability

GPSA - Swing Away for easy transfers



THE ORIGINAL RIGID FRAME



The GP's classic box frame and one piece construction exceed our superior standard for quality and deliver an unmatched rigid feel.

SPORTS OPTION



The tapered GPV front frame, optional Quickie Performance wheels, and up to 12° of wheel camber combine for exceptional maneuverability on and off court.



RIGID MANUAL FRAME

NO-NONSENSE UTILITY AND UNBEATABLE VALUE



COMMUTER FRIENDLY DESIGN



- Ultra lightweight and durable aluminum frame, WC-19 Transit approved up to 120kgs
- Open frame and fluid caster housing design for a stylish look
- Small footprint for maneuverability and easy transport

COMFORTABLE FIT



The QRi comes standard with 0" or 1" of frame inset, a highly adjustable axle plate for perfecting your Center-of-Gravity setting, and multiple front frame angle options.



9







Krypton R	7 ^R	5 ⁸
32 - 46 cm (12-18")	12-20"	12-20"
34 - 48 cm (14-20")	12-20"	12-22"
43- 55 cm (17-22")	16-21"	16-21"
38 - 50 cm (15-20")	13-20"	13-20"
Carbon Fibre Open Frame	Open Series 7000 Aluminium	Open Series 5000 Aluminium
25 - 45 cm (10-18")	8-20"	10-20"
-31° - 15°	-27° - 8°	-17° - 8°
92° / 100°	75° / 80° / 85°	75° / 80° / 85°
0"/ 1"	0"/ 1"/2"	0"/ 1"/2"
0 °/ 3 °	0 °/ 3 °/ 6 °	0°, 3°, 6°
125 kg	120 kg	120 kg
3.65 kg transit / 6.2kg	8.5 kg	9.8 kg
No	Yes*	Yes*
5 Finishes	29 colours	28 colours
Limited Life time	Limited Life time	Limited Life time
100934	100167	100167

*Conforms to ANDI/RESNA WC-19 Standards. Please refer to the Quickie & Zippie WC-19 Transit Securement Supplement

SEAT WIDTH SEAT DEPTH TRONT SEAT HEIGHT REAR SEAT HEIGHT FRAME STYLE HEIGHT BACKREST ANGLE INSET

MAX. USER WEIGHT PRODUCT WEIGHT TRANSIT APPROVED FRAME COLOUR WARRANTY

Live without limits.

QUICKIE

RIGID FRAME SPECIFICATIONS









GP	GP SA	GPV	QRi
12-22"	12-20"	12-22"	12-20"
12-22"	12-22"	12-22"	12-20"
15.5-23"	16.5-23"	15.5-23"	15-20"
14-20.5"	16-20.5"	14-20.5"	14-20"
Traditional	Traditional	Traditional	Open Series 5000 Aluminium
8-20"	8-20"	8-20"	8-20"
-4° - 8°	-4° - 8°	-4° - 8°	-8° - 17°
60° / 70° / 75°/ 80°	60° / 70° / 90°	60° / 70° / 75°/ 80°	75° / 80°/ 85°
N/A	N/A	1" V taper	0, 1" V inset
(Adjustable) Set at 3°	(Adjustable) Set at 3°	(Adjustable) Set at 3°	0° / 3°
113 kg	113 kg	113 kg	120 kg
12 kg	11 kg	11 kg	10 kg
No	ISO Only - thru frame	ISO Only - thru frame	Yes*
28 colours	28 colours	28 colours	24 colours
Limited Life time	Limited Lifetime	Limited Life time	Limited Life time
100167	100167	100167	100167











12

KRYPTON F CARBON FIBRE FRAME

REVOLUTIONARY CARBON FRAME WITH STAGGERING STRENGTH-TO-WEIGHT RATIO

FOLDING FRAME

Ingeniously designed cross-brace with an open frame feel!



CARBON BRAID TECHNOLOGY



This QUICKIE Krypton Range features a patented and fully automated BRAID process that makes it possible to build the perfect frame – strong, light and manufactured with a precise and reliable automated process.

LIGHTEST FOLDING WHEELCHAIR



From just 5.8kgs in transit, the QUICKIE Krypton F is our lightest folding wheelchair yet! From every angle, the Krypton F is conceived to portray beauty. The heart of the Krypton F is its ingeniously designed cross-brace that fits so neatly under the seat it's barely noticeable giving a very minimalistic, open-frame feel.



XENON² FOLDING FRAME SERIES

ONE OF THE LIGHTEST FOLDING WHEELCHAIRS IN THE WORLD!

Available in: FF - Open fixed front frame

DUAL - Reinforced fixed front frame

SA - Swing Away for easy transfers



EASY TO FOLD MECHANISM



The aluminium latches on the cross-brace will also keep your wheels in parallel when folding and makes sure the chair stays folded when lifting.

ADJUSTABLE TO SUIT YOU



You want to change your chair driving characteristic from passive to more active or the opposite way? Just change the centre of gravity setting by simply repositioning the axle stem. And don't forget to adjust the back angle for your ergonomic posture.



QUICKIE 2 FOLDING WHEELCHAIR

VERSATILE & MODULAR. THE INDUSTRY'S FAVORITE FOR OVER 25 YEARS!

Available in:

Q2 - Standard Modular frame

Q2 Lite - Ultra lightweight, modular frame

Q2 HP- High Performance



MODULAR FRAME DESIGN



Featuring a modular design, the QUICKIE 2 offers the broadest range of possible configurations in the industry: over 100,000 different complex rehab options and counting. Designed to meet all of your unique needs without compromising fit or performance, only one wheelchair can live up to the name QUICKIE.

QUICKIE 2 LITE



Love the Quickie 2, but looking for something with an edge? The Quickie 2's minimalist counterpart, the Quickie 2 Lite, features an extremely lightweight, modular frame with a rounded rear for unique style and either fixed or swing-away front hangers.



GXI FOLDING WHEELCHAIR

LOW MAINTENANCE & SIMPLICITY

DURABLE FRAME DESIGN



- Ultra lightweight and durable folding frame with a 120kgs weight capacity
- One-piece side frames reduce moving parts
- Three arm crossbrace sits inside frame for rigid ride

SIMPLE ADJUSTMENT & SERVICING



Most adjustments can be made in less than one minute, using only one tool - Center-of-Gravity (axle plate), Seat-to-floor height, Camber, Axle spacing and Back height



OUICKIE

MG FOLDING BARIATRIC

IDEAL FOR BARIATRIC CLIENTS WEIGHING UP TO 295 KG & UPTO 75CM WIDE

MODULAR & HIGH-STRENGTH



The M6's modular frame, features high strength steel for increased durability. Its interchangeable, heavy duty components accommodate changing conditions with simple parts replacements.

LOCKING CROSSBRACE & STABILSER BAR



The M6 utilises a double locking crossbrace and stabiliser bar to minimize frame flex and increase durability. The frame will hold an active user upto 295kgs.







Krypton F	Xenon ² SA [*]	Xenon ² FF [*]
32-46 cm (12.8-18.4")	32-50 cm (12.8-20")	32-46 cm (12.8-18.4")
34-48 cm (13.6 - 19.2")	34-50 cm (13.6-20")	34-50 cm (13.6 - 20")
43-55 cm (17.2-22")	43-56 cm (17.2-22.4")	43-55 cm (17.2-22")
37-50 cm (14.8 - 20")	37-50 cm (14.8-20")	37-50 cm (14.8 - 20")
Carbon Fibre Open frame	Reinforced with swing-away front	Open frame with fixed front
25 -47.5 cm (10-19")	25 -47.5 cm (10-19")	25 -47.5 cm (10-19")
75° to 103°	75° to 103° in 4° increments	75° to 103° in 4° increments
92° / 100°	70° / 80°	92° / 100°
0-1 cm	N/A	0-1 cm
0° / 2° / 4°	0° / 2° / 4°	0° / 2° / 4°
110 kg	125 kg	110 kg
from 5.8kgs transit / 8.3kgs	10.3 kg	from 8.8 kg
No	Yes	Yes
5 Finishes	32 colours	32 colours
Limited Lifetime	Limited Lifetime	Limited Lifetime
100934	100934	100934



Live without limits.

QUICKIE

.....

..... um i 1111 11111 1111 111111













Xenon ² Dual [®]	Quickie 2 Lite	Quickie 2	QXi	M6
32-50 cm (12.8-20")	14-20"	12-24"	14-20"	22-30"
34-50 cm (13.6-20")	14-20"	12-22"	14-20"	18-22"
43-55 cm (17.2-22")	14.5-21.5"	13.5 to 22"	12.5-19.5"	17-20"
37-50 cm (14.8 - 20")	13-19.5"	13 to 22"	12-19"	17-20"
Reinforced with fixed front	Traditional - Fixed or SA front	Traditional- Fixed or SA front	Traditional	Heavy Duty
25 -47.5 cm (10-19")	9-21"	9-21"	13-20"	17-20"
75° to 103° in 4° increments	-14° - +5° (optional)	-15° - +3° (Optional)	-15° - 3°	N/A
92° / 100°	70°/ 80°	60°/ 70°/ 80°/ 90°	60°/ 70° / 80°	N/A
0-1 cm	0 / 1" each side	-	-	N/A
0° / 2° / 4°	(Adjustable) Set at 3°	(Adjustable) Set at 3°	(Adjustable) Set at 3°	0°
140 kg	120 kg	136 kg Std / 158 kg HD	136kg	295 kg
9.5 kg	10.6 kg	11.5kg	13.25kg	24 kg
Yes	Yes	Yes	Yes	No
32 colours	24 colours	24 colours	24 colours	24 colours
Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime	Limited Lifetime
100934	100167	100167	100167	100167





IRIS TILT IN SPACE

INTELLIGENT ROTATION IN SPACE TECHNOLOGY

ife R

One of the best and lightest tilt-in-space wheelchairs on the market

INTELLIGENT ROTATION



The Quickie IRIS features up to 55° of Intelligent Rotation in Space Technology to achieve superior positioning and is designed to facilitate feeding and respiratory function, reduce pressures and improve visual alignment.

MONO BACKREST SYSTEM



The MONO Backrest System's single back post offers limitless seating and positioning capabilities. By mounting the backrest to a single central post, the backrest and chair widths may be selected independently to accommodate users with varying torso and hip widths.



IRIS WITH POWER TILT

HAVE MORE CONTROL OVER THE POSITION OF YOUR BODY IN SPACE

Available in:

- IRIS
- IRIS Heavy Duty
- improves comfort by assisting with sitting tolerance
- provides more independence to participate in daily activities



ATTENDANT CONTROL




Allows the carer to give control of the tilt function to the user. The battery gauge is displayed on the top of the control so you never have to guess when it's time to charge. Ability to conveniently hang the attendant control anywhere on the chair.

TOGGLE SWITCH



Can be mounted vertically or horizontally depending on needs. Charging port is conveniently located on the toggle switch for easy access.



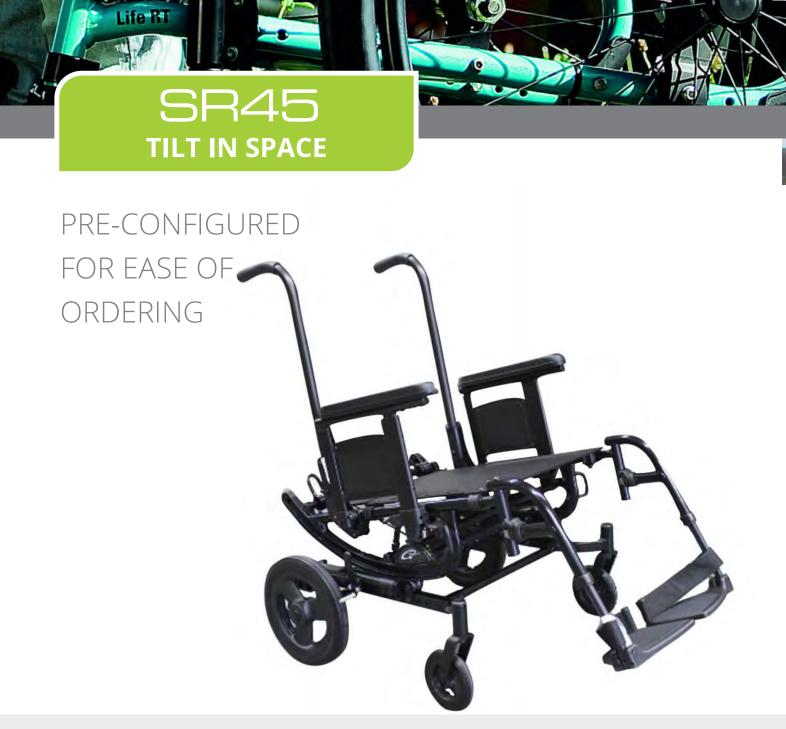

45° OF ROTATION-IN-SPACE



The Quickie SR45's seat frame rotates around the client's centre of gravity providing the shortest possible wheelbase without compromising stability. When the centre of gravity and centre of rotation are properly aligned, tilting is smooth, easy and controlled.

STANDARD SMART PACKAGE



SR45 comes standard with - Non folding back canes, Height adjustable back with adjustable stroller handles, Foot tilt, Transit tie downs, Dual post height adjustable flip back with full length arm pads, 70` Swing in/ out hangers, solid seat pan and rear anti tips and composite foot plates.



22

TILT IN SPACE SPECIFICATIONS





	IRIS	SR45
SEAT WIDTH	14-22"	16x16"
SEAT DEPTH	14-22"	18x18"
FRONT SEAT HEIGHT	12.5-19"	20x18"
REAR SEAT HEIGHT	12.5-19"	17"
FRAME STYLE	Rotation In Space	Rotation In Space
BACK HEIGHT	15-24"	21"
BACKREST ANGLE	-5° to +30°	-5° to +30°
LEG RESTS	60°/ 70°/ 80°/ 90°/ Elevating	70°
INSET	N/A	N/A
CAMBER	0°	0°
MAX. USER WEIGHT	136 kg / HD 181kgs	120 kg
PRODUCT WEIGHT	17.6 kg	17.6 kg
TRANSIT APPROVED	Yes	Yes
FRAME COLOUR	25 colours	Black
WARRANTY	Limited Lifetime	Limited Lifetime
ARTG	100167	100167



23

31

MANUAL REFERENCE GUIDE



Seating Selection



To learn about personalising JAY seating products, visit the JAY Your Way product page on www.SunriseMedical.com.au



Sunrise Medical Australia Pty Ltd 11 Daniel Street, Wetherill Park, NSW 2164 Phone: (02) 9678 6600 Email: enquiries@sunrisemedical.com.au www.sunrisemedical.com.au