








# THE POWER OF QUICKIE

A complete power range from basic mobility to complex rehab. Whatever you need - without compromise.

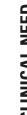



**BASES** 



SEDEO LITE CAPTAIN SEAT RWD



**Q400** 

BASES SEDEO LITE MWD CAPTAIN SEAT FWD RWD



Q500 / Q500 Mini

BASES

MWD

FWD

SEDEO PRO SEDEO LITE CAPTAIN SEAT RWD



**Q700** 

**SEATING** SEDEO ERGO SEDEO PRO CAPTAIN SEAT RWD

BASES

MWD

**FWD** 



Q700-UP

SEDEO ERGO UP MWD

**BASES** 

**FWD** 

RWD: Rear-Wheel Drive MWD: Mid-Wheel Drive FWD: Front-Wheel Drive

**BASE PERFORMANCE** 



# THE NEW QUICKIE Q-SERIES **POWER RANGE: A COMPLETE POWER PRODUCT PORTFOLIO UNDER ONE ROOF.**

A wide range of bases, drive tractions, seating systems, electronics, controls and accessories, that can be combined to create the perfect powerchair for each individual need.



# **BASES**

**SERIES** 

Q700 Q500 0.400

Q300

Front-Wheel Drive Mid-Wheel Drive Rear-Wheel Drive



**SEATING SYSTEMS** SEDEO ERGO UP SEDEO ERGO SEDEO PRO SEDEO LITE



### **CONTROLS**

**CAPTAIN SEAT** 

VR2

R-Net

**SWITCH-IT** Controls

Link-It Mounting System

Remote Seating App



# CHOICE WITHOUT COMPROMISE

The new QUICKIE Q-Series offers full flexibility to create the perfect powerchair for each individual user. Highly flexible configurations, interchangeable parts, quick and easy adjustments.











#### **BEST IN CLASS SEATING**

Q-Series offers a range of seating systems to cover any clinical need. A choice of five seating systems: SEDEO ERGO UP, SEDEO ERGO, SEDEO PRO, SEDEO LITE and CAPTAIN SEAT. The SEDEO range is fully compatible with JAY cushions and backs and WHITMYER headrests to assure the best seating and positioning.

### **BEST DRIVE BASES CHOICE IN THE MARKET**

From cost effective compact bases to high performance outdoor bases, all with the most extensive choice of drive tractions. An impressive range of base levels, and drive tractions that can accommodate all sort of user preferences, lifestyle needs and pricing/reimbursement requirements.

## WITH MOST ADVANCED TECHNOLOGIES

Advanced engineering and innovative technologies applied to improve people's lives: SpiderTrac<sup>™</sup> 2.0 and Full Track Suspension, Memory Seating, advanced anti-shear backrests, biometrical repositioning, Remote Seating app (pressure management), SWITCH-IT controls and Link-It mounting systems and advanced R-Net Controls.

### **EASIER TO LEARN**

The common design and the modularity create a simplified user experience, that makes it very easy to learn, manage and prescribe.

# **EASIER TO SERVICE**

The modularity and flexibility of the bases, seatings and controls allows an easy adaptation to different users, making it perfect for refurbishments and allowing a great fleet management.

# **SERIES OVERVIEW**

	Q700-UP	Ω700	Q500	Q500 Mini	Q400	Q300
Outdoor performance	••••	••••	• • • • •	• • • •	••••	• • • • •
Indoor performance	• • • •	• • • •	• • • • •	• • • • •	• • • • •	••••
Suspension	All Wheels	All Wheels	Drive wheels, Advanced	Drive wheels, Advanced	Drive wheels, Regular	Drive wheels, Basic
Motors	4 pole	4 pole	4 pole	2 pole	4 pole	2 pole
Speed	8,10 kph	8,10 kph	10 kph	10 kph	10 kph	10 kph
Batteries	60 Ah / 72 Ah	60 Ah / 72 Ah	60 Ah	34 Ah	60 Ah	34 Ah
Max. user weight	120 kg	160 kg	136 kg	136 kg	136 kg	136 kg
Drive bases	MWD / FWD	MWD / FWD / RWD	MWD / FWD / RWD	MWD	MWD / FWD / RWD	MWD / RWD
Seating	SEDEO ERGO UP	SEDEO ERGO / SEDEO PRO / CAPTAIN SEAT (MWD 136 kg only)	SEDEO PRO / SEDEO LITE / CAPTAIN SEAT (MWD only)	SEDEO PRO / CAPTAIN SEAT (MWD only)	SEDEO LITE / CAPTAIN SEAT (MWD only)	SEDEO LITE / CAPTAIN SEAT (MWD only)
Electronics	R-Net	R-Net/VR2 CAPTAIN SEAT	VR2/R-Net	R-Net	VR2	VR2
Crash tested	Yes	Yes	Yes	Yes	Yes	Yes

# CREATED INDIVIDUALLY FOR YOUR CLINICAL AND LIFESTYLE NEEDS.



# QUICKIE Q-SERIES OVERVIEW

### 0700 SERIES: THE ULTIMATE EXPERIENCE

Chairs designed to offer a superb outdoor experience while allowing a good indoor use. A choice of three high performance drive bases with market leading suspension technologies, powerful motors and long lasting batteries to provide the ultimate outdoor performance. Combined with a choice of complex rehab seating systems (SEDEO ERGO UP STANDING WHEELCHAIR, SEDEO ERGO & SEDEO PRO) to match users with specific or high clinical needs, and the most advanced electronics to access their world. Including also Captain Seat option.

# **0500 SERIES: MAKE YOUR OWN ADVENTURE**

Robust outdoor powerchairs with a choice of three drive bases with advanced suspension technologies and powerful motors, for a superior outdoor experience. Combined with a choice of rehab seating systems (SEDEO PRO, SEDEO LITE & CAPTAIN SEAT) to provide great support and positioning to users from mild to medium clinical needs, and the most advanced electronics to access their world.

# Q500 MINI: ULTRA COMPACT PERFORMANCE SUITABLE FOR COMPLEX REHAB NEEDS

Ultra compact wheelchair with highly adjustable and customisable seating and electronics, a combination to match users with high clinical needs, but tight indoor environments.

# **Q400 SERIES: ADVANCED OUTDOOR / INDOOR PERFORMANCE**

Robust outdoor powerchairs with a choice of three drive bases, with regular motors and suspension for a great outdoor experience with advanced electronics. Available with SEDEO LITE seating systems and Captain Seat.

# **Q300 SERIES: ULTRA COMPACT EVERY DAY PERFORMANCE**

Entry level chairs with ultra-compact mid or rear wheel drive bases for great indoor agility with the stability and traction for outdoor navigation and community access.





# **PERFORMANCE BASES**

# **DRIVE BASES OVERVIEW**

The new QUICKIE powerchairs range offers the most complete choice of drive bases on the market, available in different series and with many seating systems.

The choice of a drive base affects how the chair performs in different situations and environments. The right choice for each user will depend on personal preferences, experience, clinical, lifestyle and environmental needs.

	MWD	FWD	RWD
Manoeuvrability	••••	••••	• • • • •
Stability	• • • •	••••	••••
Kerb climbing	• • • •	••••	••••
Traction	••••	••••	• • • • •
Drive comfort	••••	••••	• • • • •
Ease of driving	••••	••••	••••
	MWD	FWD	RWD
Q700 UP	<b>~</b>	<b>~</b>	_
Ω700	<b>~</b>	<b>~</b>	<b>✓</b>
Ω500	<b>~</b>	<b>~</b>	<b>✓</b>
Q500 Mini	<b>~</b>	_	_
Q400	<b>✓</b>	<b>✓</b>	<b>✓</b>
Ω300		_	



# **MID-WHEEL DRIVE BASES** (MWD)

The mid-wheel drive bases have the best indoor performace. Their smallest turning circle provides great manoeuvrability. They are very intuitive and easy to drive, ideal for users with reduced hand control, driving with switches and users with poor spatial awareness.



# **FRONT- WHEEL DRIVE BASES** (FWD)

The front-wheel drive bases have great outdoor performance with good indoor manoeuvrability. Knee angle can be positioned at 90°. Maintains control on steep declines.



# **REAR-WHEEL DRIVE BASES** (RWD)


The rear-wheel drive bases are stable and have good kerb climbing ability. Intuitive drive that provides good control with non-proportional inputs.





	SEDEO ERGO UP	SEDEO ERGO	SEDEO PRO	SEDEO LITE	CAPTAIN SEAT
User clinical need	High / Standing	High	Medium to High	Medium to Low	Low
Seat width	400 - 500 mm	400 - 560 mm	380 - 560 mm	385 - 510 mm	16 - 20 Inches
Seat depth	400 - 540 mm	400 - 560 mm	400 - 560 mm	410 - 560 mm	16 - 22 Inches
Backrest height	560 - 610 mm	560 - 700 mm	440 - 620 mm	460 - 510 mm	25 inches
Configurability	Low	Low	High	Medium	_
Standing function	<b>✓</b>	_	_	_	_
Lay down position	<b>✓</b>	<b>✓</b>	_	_	_
Programmable seating	<b>✓</b>	<b>✓</b>	_	_	_
Powered recline	90° - 175° biomechanical	90°-175° biomechanical	90°-130°/90°-150° biomechanical	85°-120°	_
Manual recline	_	_	90°-120°	90°-120°	90°-115°
Anti-shear back movement	180 mm	180 mm	100 mm (optional)	_	_
Pressure release reminders	<b>~</b>	<b>✓</b>	_	_	_
Tilt angle	22°	50°	50°	30°/50° *	30°/50° *
Lift	200 mm	300 mm	300 mm	300 mm	300 mm
JAY & WHITMYER compatibility	<b>✓</b>	<b>✓</b>	<b>~</b>	<b>~</b>	JAY Cushion
Remote Seating app	<b>✓</b>	<b>✓</b>	_	_	_
Max. user weight	120 kg	160 kg	160 kg	136 kg	136 kg
Crash Tested	Yes, ISO	Yes, ISO	Yes, ISO	Yes, ISO	Yes, ISO

<sup>\* 50°</sup> tilt available on Q500 and Q700 Only

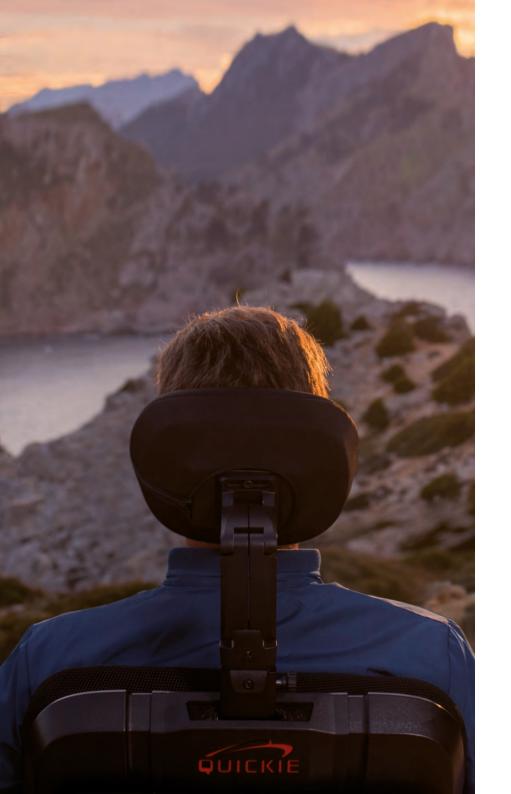



# **SEDEO ERGO UP**

High-end biomechanical seating, with many programmable seating positions and standing function. All the benefits of SEDEO ERGO seating with additional standing function. SEDEO ERGO UP provides a standing angle up to 85° with a highly rigid and stiff standing module.



# **SEDEO LITE**

Multi-adjustable seat frame with the option of either SEDEO LITE backrest cushions, Jay backrests or sling back upholstery. SEDEO LITE is designed for users with medium to low clinical needs that require adjustable seating system. The backrest, armrests and legrests can be adjusted to fit the user. Includes seat frame with seat rail for additional mounting solutions.



# **SEDEO ERGO**

High-end biomechanical seating, with many programmable seating positions. SEDEO ERGO covers the most demanding clinical needs with 180 mm anti-shear back movement, biometric repositioning function, infinitely adjustable seat and back, easy entry system, 90° knee angle in seated position and pressure relief reminders.



# **CAPTAIN SEAT**

Provides a relaxed sitting position and features a semi-recline, fold-down backrest, adjustable seat size, and convertible seat pan that allows you to use your preferred cushion. Available with power lift or 30°/50° tilt and power centermount legrest.



# **SEDEO PRO**

High-end fully adjustable and configurable seating frame with a wide range of comfort cushions and positioning options. Aimed at meeting moderate to complex clinical needs. With a wide range of armrests, armpads, legrests, seat and back cushions, headrests and positioning accessories, SEDEO PRO offers great flexibility to adapt to the individual user needs.



# WHITMYER & JAY SEATING

Compatible with all SEDEO seating systems, the JAY family of cushions and backs offers superior clinical seating for pressure management, positioning and comfort. WHITMYER headrests are also compatible with all SEDEO seating systems offering the best head support and positioning.

# **CONTROLS & CONNECTIVITY**

# VR2 & R-NET ADVANCED JOYSTICKS

The QUICKIE Q-Series wheelchairs come with the reliable VR2 control or the superior R-Net advanced joystick. A large, easy-to-read HD display and simple to control paddle switches allow for intuitive operation. (Depends on model).

## **RELIABLE**

Industry standard VR2 and R-Net systems offer easy to use and highly durable drive controls.

#### **PROGRAMMING**

VR2 and R-Net both provide advanced programming to customise the chair to your requirements and environments.

#### CONNECT

Seamlessly take control of your environment using the integrated Infrared functionality. Control your home entertainment system, TVs, DVDs, docking stations and other home entertainment systems simply through your joystick (R-Net only).

#### **BLUETOOTH & IDRIVE**

Control your phone, computer, tablet and any other Bluetooth and i-devices with this useful optional extra. Connect and control up to four devices (R-Net only).

# **QUICKIE'S PATENTED ASSIGNABLE BUTTONS**

Change the function on any button on your joystick to any command you desire, thanks to QUICKIE's patented assignable buttons. You can also assign a second (long press) function to buttons, giving you quicker, simpler and smoother access to functions you commonly use without having to use your joystick (R-Net only).



If driving with a traditional joystick is not the best option, alternative drive controls allow you to fully experience your chair's driving capabilities and benefits.


# **SWITCH-IT**INNOVATIVE CONTROLS

Conveniently access your 5 programmable positions with the Ctrl+5-a perfect combination with Memory Seating (Memory Seating with SEDEO ERGO and SEDEO ERGO UP only). It allows you to go straight to the function you need at the touch of a button (R-Net only).

# **LINK-IT**MOUNTING SYSTEM

Need to use additional switches or special controls? No problem, easily connect them using the universal Link-It mounting system.



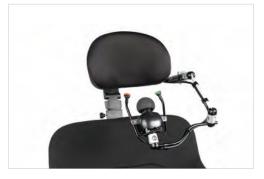















# DISTINGUISHED STYLE

With the Q-Series, we've taken style to the next level. With a Matte Black chassis finish, hidden wiring and components combined with sharp-tailored upholstery,

it's a sleek and a contemporary design with a variety of colours to choose from.

# Fancy a change?

A choice of contemporary colour inlays to perfectly complement its striking Matte Black bodywork.

Simply replace the colour inlays. Styling your powered wheelchair has never been so simple.



\*Q500 Mini:

\*Q300:



<sup>\*</sup> The number of colours available may differ depending on the seat and base combination.



















Sunrise Medical Australia Pty Ltd. 6 Healey Ciruit Huntingwood NSW 2148 Phone: (02) 9678 6600 Fax: (02) 9423 6969

Email: enquiries@sunrisemedical.com.au www.SunriseMedical.com.au

© 2019 Q-Series Range Brochure AU Rev B

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult **Sunrise Medical** with any queries you may have.

Pictures and illustrations may be overseas models and some seating and accessories not available in Australia.

If you are visually impaired, this document can be viewed in PDF format at SunriseMedical.com.au.

