



Product Information

# SOFIHUB home



SOFIHUB home Care has been developed for a single person home and provides greetings, reminders and alerts.

SOFIHUB home Care is packaged with a hub and 8 sensors.


The unique intuitive technology assists in visualising behaviours and then raises an alert when behaviours are different.

Some alerts may relate to such things as,

- Late to wake up
- Late to bed
- Long bathroom duration

Carers of the resident are provided with a SOFIHUB home login and password which they can use to login to the SOFIHUB portal online interface and customise the reminders.

- INCLUDES:
- 1 SOFIHUB home + 8 Sensors
  - 4G Connectivity + WiFi
  - SOFIHUB portal Access
  - Alerts/Notifications

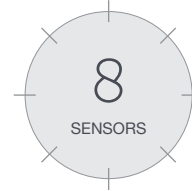

INTELLIGENT MESSAGES:

**Greeting**

*“Good morning...”*

Every morning SOFIHUB home lets you know she is there with a greeting and today’s weather.

**Medication**

*“It looks like you have forgotten to take your 12 o’clock medication.”*

SOFIHUB home medication sensor knows when you forget to take your medication and gives you a gentle reminder.

**Hydration**


*“Good afternoon, it’s very hot today, make sure you drink plenty of water.”*

**Doctor’s Appointment**

*“Good Afternoon! Don’t forget, you have an appointment with Dr. Jones at 3pm.”*

**Bathroom**

SOFIHUB home will send your carer/family member a notification if you are in the bathroom for longer than expected and that you might need some help.

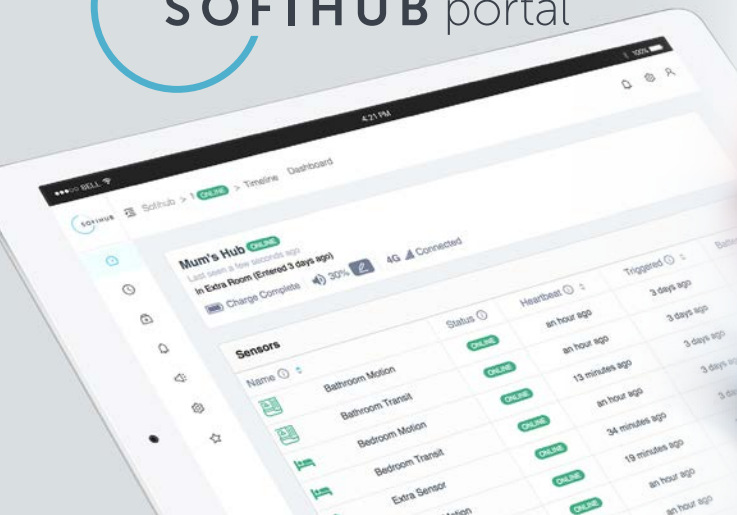

## What is Intelligent Message?

**Time based messages:** Play at the specified time.

**Medication reminders:** A reminder message will play if the medicine cabinet is not accessed at the scheduled time.

**Occupancy based messages:** If the room that the hub is located in is occupied between the start and the expiry times, the message will be played.

Features	Benefits
Portal - online access for SOFIHUB home and SOFIHUB beacon	Realtime access for families and carers and notifications
Timeline reporting	View daily activities
Ability to enter reminders and add intelligent messages	Improvements to routine



# SOFIHUB home

## Specifications

### Hub

<b>Size</b>	95 x 95 x 230mm high
<b>Power</b>	240 Volt with battery with up to 4 hour backup
<b>Speaker</b>	105db maximum
<b>Connectivity</b>	4G Connectivity + WiFi

### Sensors

<b>Size</b>	45 x 45 x 40mm high
<b>Power</b>	Battery powered
<b>Mounting</b>	Self adhesive mounting system
<b>Battery</b>	CR123A



## Installation Information

SOFIHUB home is supplied fully packaged with the main hub and sensors. The sensors are placed in key locations to assist in visualising regular movement patterns.

SOFIHUB home will be fully operational out of the box only requiring the hub to be connected to mains power. The SOFIHUB home sensors are battery operated. The SOFIHUB portal will provide the status of the battery, of each sensor and will alert you when the batteries require replacement.

SOFIHUB home is fully connected via the mobile network to send alerts and receive software updates. In case of power cuts SOFIHUB home has a built in power back up that will last up to 4 hours.

Once installed SOFIHUB home can connect to WiFi if available and operate under mobile network as a back up.

Clear instructions are provided to assist you with the sensor placement and connecting the sensors to the hub. No drilling or wiring is required.

SOFIHUB can also provide a full installation service by our own qualified technicians.

## SOFIHUB portal

Carers of the resident are provided with a SOFIHUB home login and password which they can use to login to SOFIHUB home online to interface and customise the reminders. Carers are able to view history, view timeline and learn the resident's routine.

