VICAIR[®] ADJUSTER 02

THE OPTIMAL CUSHION FOR WHEELCHAIR USERS WITH BODY ASYMMETRY

Pre-ischial bar limits sliding forward | converts shear into normal forces Follows the body's contours for equal pressure distribution | support ITcoccyx area | managing pelvic asymmetry

_					
Thi	σh	I SU	Ini		rt.
	5.	30	ואי	50	

Individual solution for asymmetry

Minimizes sliding and shear forces

Optimal adjustment to the body's contours for equal pressure distribution

The Adjuster O2 is extremely suitable for use with asymmetric body contours and amputations. It automatically manages small obliquities and is individually adjustable to support more extreme asymmetry. When you're looking for a high comfort level, optimal skin protection and stability: the Adjuster O2 is always a safe choice. The design allows total immersion to equally redistribute pressure and limits sliding forward. 100% machine washable and breathable, for optimal hygiene and perfect temperature and moisture regulation. As a result the Adjuster O2 offers you optimal protection against skin breakdown.

PRODUCT PERFORMANCE

Skin protection / Anti decubitus	
Pressure redistribution	
Shear force reduction	
Micro climate control	
Impact & vibration damping	

Positioning

Positioning	
Lateral stability	
Frontal stability	
Sliding reduction	
Allows manual adjustments	
Allows easy transfers	

Maintenance & Safety	
Durability	
Reliability	
Functionality when misused	
Low maintenance level	
Level of hygiene	
Functionality & Application	

runedonancy o Application	
Comfort	
Low cushion weight	
Available size range	
Maximum user weight	
Low seat to floor height	
Auto adjustments	

Mid height zip closure for easy adjustment when client is seated | cushion keeps its shape, enhancing the user's skin protection



USED WITH

Hemiplegia (CVA), amputations, foot propulsion, kyphosis, pelvic obliquity/ tilt/rotation, Spina Bifida, spasms, Fibrodysplasia Ossificans Progressiva (FOP)

TECHNICAL SPECIFICATIONS		
Height	6 and 10 cm	
Average weight	725 grams	
Max. user weight	250 kg	
Machine washable	60°C max.	

CUSHION SIZES

LL	CUSHION SIZES											
	Standard (Sizes >53 cm height 10 only)											
	On	req	ues	st								
												width
Inch		14	15	16	17	18	19	20	21	22	23	24
	СМ	35	38	40	43	45	48	50	53	55	58	60
14	35											
15	38											
16	40											
17	43											
18	45											
19	48											
20	50											
21	53											
22	55											
23	58											
24	60											

7

THE MOST STABLE CUSHION FOR WHEELCHAIR USERS

VECTOR 02

Ischial well - Two compartments for superior pressure reduction in the IT-coccyx area & extra positioning ... possibilities

VICAIR

Trochanter support firmer density to provide lateral stability | redistributes pressure to help create an ischial well

Front compartments provide leg alignment and stability | easily adjustable for individual positioning needs **Mid height zip closure** for easy adjustment when client is seated | cushion keeps its shape, enhancing the user's skin protection

Nine adjustable compartments for optimal stability

Expanding the user's reach

Extra positioning of the thighs

Sometimes it is necessary to reduce pressure in the vulnerable area around the IT's (sitting bones) and coccyx (tailbone). The Vector O2 does just that. Its unique compartment layout further provides extra pelvic stability, excellent leg guidance, optimal skin protection and comfort.

100% machine washable and breathable, for optimal hygiene and perfect temperature and moisture regulation. The Vector O2 offers support where needed, resulting in improved functionality of the upper body and expanding reach.

PRODUCT PERFORMANCE Skin protection / Anti decubitus Pressure redistribution Shear force reduction

Positioning

i obitioning	
Lateral stability	
Frontal stability	
Sliding reduction	
Allows manual adjustments	
Allows easy transfers	

Maintenance & Safety			
Durability			
Reliability			
Functionality when misused			
Low maintenance level			
Level of hygiene			

Functionality & Application			
Comfort			
Low cushion weight			
Available size range			
Maximum user weight			
Low seat to floor height			
Auto adjustments			




USED WITH

Spinal cord injuries, neuromuscular diseases, hypermobility syndrome, Ehlers Danlos Syndrome (EDS)

TECHNICAL SPECIFICATIONS			
Height	6 and 10 cm		
Average weight	725 grams		
Max. user weight	250 kg		
Machine washable	60°C max.		

CUSHION SIZES

Standard (Sizes >53 cm height 10 only)
 On request

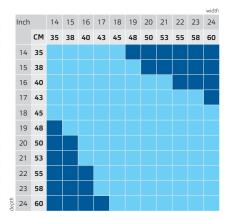

FOR THE ACTIVE WHEELCHAIR USER

ACTIVE **0**2

Ischial well - Two compartments for superior pressure reduction in the IT-coccyx area & extra positioning possibilities

VICAIR®

Trochanter support firmer density to provide lateral stability | redistributes pressure to help create an ischial well

Ergonomically-shaped foam for easy transfers



Perfect for users who transfer frequently

Excellent lateral and frontal stability

Features a nifty storage pouch for personal items

The Active O2 offers the unique combination of ergonomically shaped foam and Vicair SmartCells. The four rear SmartCell compartments provide a high level of comfort and skin protection, and reduce pressure on the sitting bones and tailbone area. The foam front makes transfers easy and feels steady and improves wheelchair handling. The foam inlay is protected by an easy-toclean fluid-resistant liner. The rest of the cushion is 100% machine washable and breathable for optimal hygiene and perfect temperature & moisture regulation. The Active O2 has an easy accessible storage pouch for precious belongings (keys, phone, wallet). In short: everything you need to maintain your active lifestyle!

PRODUCT PERFORMANCE Skin protection / Anti decubitus

Shear force reduction	
Micro climate control	
Impact & vibration damping	

Positioning

i ositioning	
Lateral stability	
Frontal stability	
Sliding reduction	
Allows manual adjustments	
Allows easy transfers	

Maintenance & Safety	
Durability	
Reliability	
Functionality when misused	
Low maintenance level	
Level of hygiene	
Functionality & Application	
Comfort	

Connort	
Low cushion weight	
Available size range	
Maximum user weight	
Low seat to floor height	
Auto adjustments	

SmartCell – Foam barrier limits sliding forward



USED WITH

Spinal cord injuries, active users with a high transfer frequency

TECHNICAL SPECIFICATIONS			
6 and 9 cm			
1250 grams			
200 kg			
60°C (without foam)			

CUSHION SIZES	
Standard	(Sizes >53 cm height 9 only)
On request	

_	_				_			_		_		width
Inch		14	15	16	17	18	19	20	21	22	23	24
	см	35	38	40	43	45	48	50	53	55	58	60
14	35											
15	38											
16	40											
17	43											
18	45											
19	48											
20	50											
21	53											
22	55											
23	58											
24	60											

VICAIR ALL

ALLROUNDER 02

EXPLORE YOUR POSSIBILITIES





Easily adjustable straps to provide comfort and to follow the body

Stays with you during activities

Skin protection

Extremely light to support optimal mobility

Active as you are, not all of your activities require or allow the use of a wheelchair. In these situations the Vicair AllRounder O2 keeps you comfortable and helps to protect your skin. As it is easily buckled onto your body, it follows you around wherever you go.

Feel protected and secure during your sports and outdoor activities, when you sit in a normal chair or on any hard surface, or when you play on the ground. Afterwards you can simply put the entire AllRounder O2 in the washing machine then hang it up to dry and you are ready for your next adventure. Experience the freedom and explore your possibilities!!

PRODUCT PERFORMANCE


Skin protection / Anti decubitus	
Pressure redistribution	
Shear force reduction	
Micro climate control	
Impact & vibration damping	

Positioning

FUSICIONING	
Lateral stability	
Frontal stability	
Sliding reduction	
Allows manual adjustments	
Allows easy transfers	

Maintenance & Safety	
Durability	
Reliability	
Functionality when misused	
Low maintenance level	
Level of hygiene	

Functionality & Application	
Comfort	
Low cushion weight	
Available size range	
Maximum user weight	
Low seat to floor height	
Auto adjustments	

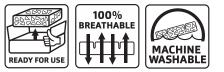

USED WITH

Activities outside of the wheelchair | When seated on hard surfaces for long periods in time.

TECHNICAL SPECIFICATIONS			
Height	n.a.		
Average weight	800 grams		
Max. user weight	100 kg		
Machine washable	60°C max.		

CUSHION SIZE

COSTINUTION	
Small	< 38 cm
Medium	38 – 42,5 cm
Large	43 – 48 cm

Determining the correct hip width Measure the width of your hips to check which size matches the size of your body. When seated, hold a tape measure horizontally above the widest point of your upper legs, closely against the front of your hips. Do not follow your leg contours, keep the tape measure straight. Measure from A to B (or B to A). Choose the larger size if the measured size is in between two sizes.

