



aidacare

WARD BED

Integrated Bed Extension

One Button Chair Positioning

Infection Control Mattress Platform

Linak Electronics & Tente Castors

Safe Working Load of 250kg



www.aidacare.com.au

1300 133 120



aidacare



The Aidacare Ward bed specifically caters to the demands of busy ward environments. The bed frame is Australian designed and offers high quality components including Linak OpenBus electronics and Tente castors.


The bed is feature rich and versatile, offering solutions for infection control, emergency resuscitation, lockout functionality and power outage.

- ✓ IEC 60601-2
- ✓ Ward Bed is registered in the ARTG, AUST R ARTG 239179

QUALITY
 Designed to IEC 60601-2 International Standards.

AUTOMATIC BATTERY BACK-UP
 Ensures bed can function in the event of a power outage.

Bed Functions

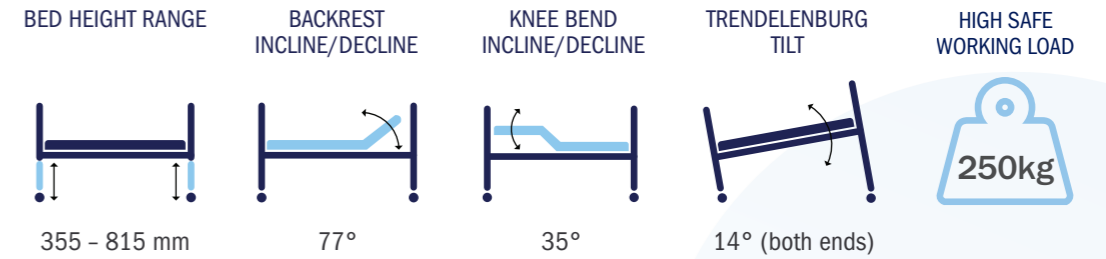

REMOVABLE MATTRESS PLATFORM
 Facilitates easy cleaning and infection control compliance.



INTEGRATED BED EXTENSION
 Provides an additional 195mm of length for changing needs.



WALL PROTECTION BUFFERS
 Rotating discs on all corners to prevent damage and injury.

WARRANTY

FRAME	10 YEARS
ELECTRONICS	3 YEARS
HANDSET	1 YEAR



ATTENDANT CONTROL PANEL (ACP)
 Facilitates patient lockout for added safety.



CENTRAL LOCKING CASTORS

German made TENTE 150mm twin castors with central locking and steering assist.



ONE BUTTON CHAIR POSITIONING FEATURE

Offers multiple clinical benefits, such as the prevention of respiratory issues.



BUILD QUALITY

Australian frame design, Linak electronics, TENTE castors and gas assisted side rails (when fitted).



IN-RAIL PATIENT EGRESS HANDLES

Assists with patient independence and bed mobility (option).



CPR QUICK RELEASE BACKREST

Emergency CPR release with assisted lowering on both sides of the bed for rapid response.



Accessories

- DELUXE SOFT DROP SIDE RAILS
- 3/4 LENGTH FOLDING UNIVERSAL SIDE RAILS
- SELF HELP POLE - STAINLESS STEEL
- IV POLE - STAINLESS STEEL
- BEECH HEAD AND FOOT BOARDS
- OXYGEN BOTTLE HOLDER - C SIZE



Dimensions & Adjustment Range

Code	BEB046000
Height Range	355 - 815mm
Overall Length	2195mm
Overall Length - Extended	2390mm
Mattress Surface Length	2030mm
Mattress Surface Length - Extended	2225mm
Overall Width	1020mm
Backrest Incline	0 - 77°
Knee Bend	0 - 35°
Trendelenburg	0 - 14°
Reverse Trendelenburg	0 - 14°
Product Weight (no side rails)	123kg
Safe Working Load	250kg
Recommended Mattress Size (L x W x H)	2000 x 900 x 170mm (max H)
Warranty	Frame - 10 years Electronics - 3 years Handset - 1 year

Accessories

Code	Description
BEA029810	Deluxe Soft Drop Side Rails
BEA029700	3/4 Length Folding Universal Side Rails
BEA029705	Self Help Pole – Stainless Steel
BEA029710	IV Pole – Stainless Steel
BEA029815	Beech Head and Foot Boards
BEA029820	Oxygen Bottle Holder – C Size



Laundry Trolleys



Patient Weigh Scales




Emergency, Theatre & Medicine Trolleys



Lifters & Slings



Instrument & Dressing Trolleys

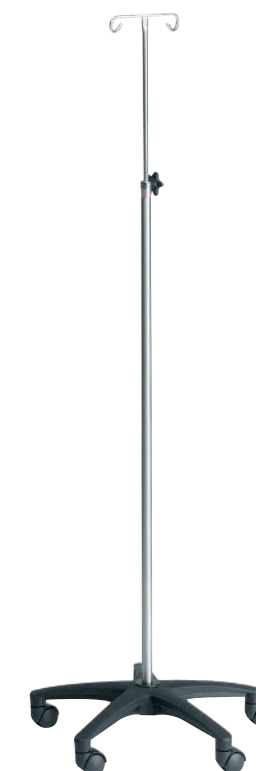

IV Stands



Overbed Tables


MATTRESS RANGE



LIFECOMFORT

Dual layer support surface incorporating a full size memory foam overlay and durable high density foam base. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for maximum infection control.


Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046410	Lifecomfort - AC Series Single	2000mm	900mm	130mm	5 Years Foam 2 Years Cover	120kg
BEM046425	Lifecomfort - Small	1920mm	820mm	130mm	5 Years Foam 2 Years Cover	120kg
BEM046460	Lifecomfort - King Single	2030mm	1050mm	130mm	5 Years Foam 2 Years Cover	120kg

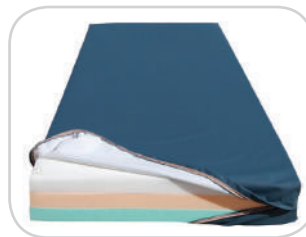

LIFECOMFORT DELUXE

Three layer PUPPS compliant support surface incorporating a full size memory foam overlay and two layers of high quality foam for maximum comfort and support. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for optimum infection control.


Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046510	Lifecomfort Deluxe - AC Series Single	2000mm	900mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046525	Lifecomfort Deluxe - Small	1920mm	820mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046455	Lifecomfort Deluxe - King Single	2030mm	1050mm	150mm	5 Years Foam 2 Years Cover	150kg

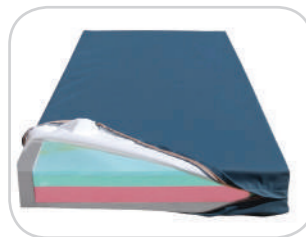

LIFECOMFORT BARIATRIC

Four layer foam, designed for larger users with a maximum user weight of 300kg. Full size memory foam overlay and side walls for added stability and ease of transfers. European fabric cover with fully welded seams, water resistant zippers and reinforced waterfall flaps for optimum infection control.

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046700	Lifecomfort Bariatric - King Single	2030mm	1050mm	200mm	5 Years Foam 2 Years Cover	300kg

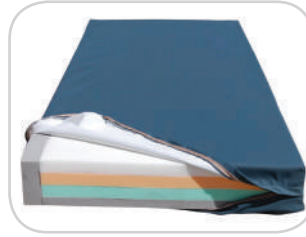

LIFECOMFORT PRESSURE REDUCING

Top of the range PUPPS2 compliant mattress featuring a quad layer support surface with additional side wall supports. European fabric cover has fully welded seams and water resistant zippers with reinforced waterfall flaps for optimum infection control.

Code	Description	Length	Width	Height	Warranty	S.W.L
BEM046610	Lifecomfort Press Red - AC Series Single	2000mm	900mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046625	Lifecomfort Press Red - Small	1920mm	820mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046660	Lifecomfort Press Red - King Single	2030mm	1050mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046675	Lifecomfort Press Red - Double	2030mm	1370mm	150mm	5 Years Foam 2 Years Cover	150kg
BEM046680	Lifecomfort Press Red - Queen	2030mm	1520mm	150mm	5 Years Foam 2 Years Cover	150kg

ASSET MANAGEMENT

Covers have individual serial numbers heat transferred into base material. This allows accurate age tracking on fixed asset registers.



INFECTION CONTROL

1. Fully welded seams
2. Reinforced waterfall flaps
3. Water resistant zippers



Benefits

Australian made covers using the highest quality European two way stretch fabric.

- Australian made foam
- Fire retardant
- Unique serial number

Removable Cover:

1. Easy inspection
2. Facilitates better hygiene
3. Fully Washable

Mattresses do not require rotation, minimising OH&S handling risks.



ALTERNATING PRESSURE SYSTEMS

FEATURES

- Compact, powerful and ultra-quiet control units
- Two way stretch fully welded top cover – high temperature washable
- Waterfall flaps over all zippers for protection and infection control
- Medical Grade power cords run inside covers to avoid hazards
- CPR emergency release function
- Static head sections to eliminate dizziness

Air-Lo and Theraflow Alternating Pressure Systems offer the latest in pressure prevention and treatment. The unique patented technology is designed to provide maximum support, release and reperfusion to prevent and heal pressure ulcers in patients up to 350kg.



AIR-LO 5 OVERLAY

The Air-Lo 5 Alternating Overlay leads the field in the prevention of pressure areas and skin breakdowns. Air-Lo 5 supports patient loads of up to 150kg and should be used in conjunction with a foam mattress or underlay. 18 pure TPU transverse air cells alternately allow pressure to be released across the skin surface. When only the very best prevention and support is required, choose Air-Lo 5 to safely care for patients assessed as at high risk of developing pressure areas.

Code	Description	Length	Width	Height	Warranty	S.W.L
PTM636630	Air-Lo 5	2000mm	900mm	130mm	1 Year Cover & Softgoods, 2 Years Pump	150kg
PTM636631	Air-Lo 5 King Single	2000mm	1070mm	130mm	1 Year Cover & Softgoods, 2 Years Pump	150kg



AIR-LO 8 MATTRESS REPLACEMENT

The Air-Lo 8 Alternating Mattress Replacement system leads the field in the critical prevention and treatment of pressure areas and skin breakdowns. Air-Lo 8 will support patient loads of up to 200kg with 20 pure TPU transverse air cells which alternately allow pressure to be released across the skin surface. The 8 inch depth allows the Air-Lo 8 to be used without an underlying foam mattress or underlay. When only the very best prevention, therapy and support is required, choose an Air-Lo 8 Mattress Replacement to safely care for patients assessed as high risk or with existing pressure ulcers. Also available in King Single size.

Code	Description	Length	Width	Height	Warranty	S.W.L
PTM636635	Air-Lo 8	2000mm	900mm	210mm	1 Year Cover & Softgoods, 2 Years Pump	200kg
PTM636636	Air-Lo 8 King Single	2000mm	1070mm	210mm	1 Year Cover & Softgoods, 2 Years Pump	200kg



THERAFLOW ALTERNATING MATTRESS REPLACEMENT

The Theraflow Alternating Mattress Replacement system utilises patented cell design technology to rapidly accelerate the deflation of cells. This alternately removes almost all interface pressure under nearly 80% of the patient's body resulting in more pressure free time. The at risk heel area is completely elevated for 30% of the alternation cycle. Theraflow supports patient loads up to 350kg with 28 pure TPU paired air cells which alternately allow pressure to be released in 14 minute cycles. Four static head cells reduce motion sensation while all outer surfaces are fully welded for infection control.

Code	Description	Length	Width	Height	Warranty	S.W.L
PTM636640	Theraflow	2000mm	900mm	250mm	1 Year Cover & Softgoods, 2 Years Pump	200kg
PTM636645	Theraflow King Single	2000mm	1070mm	250mm	1 Year Cover & Softgoods, 2 Years Pump	200kg
PTM636646	Theraflow Bariatric	2000mm	1070mm	250mm	1 Year Cover & Softgoods, 2 Years Pump	350kg

ASSOCIATED PRODUCTS



LIFTERS & SLINGS



TRANSFER AIDS



DAILY LIVING AIDS



THERAPY & EXERCISE



CRUTCHES,
FRAMES & STICKS



SCOOTERS



OVERBED TABLES



BED STICKS & POLES



FOAM MATTRESSES



AIR MATTRESSES



COMMODOES
URINALS & BED PANS



OVER TOILET AIDS



SHOWER STOOLS & CHAIRS



RAISED TOILET SEATS



SHOWER COMMODOES



BATH BOARDS & BENCHES



SEAT WALKERS



WHEELCHAIRS



STD POWER WHEELCHAIRS



PRESSURE CUSHIONS



DAY CHAIRS



LIFT & RECLINE CHAIRS

Distributed by:



ABBR03.20 110117

